



User Guide for Pulsar Operating on NUCLEUS and ED.X Devices

Version 1.0

2023.08.01

TABLE OF CONTENTS

Chapter 1.....	6
Name.....	6
Calibration.....	6
Harness Application and Usage.....	10
Home Use.....	11
PULSAR demographics part 1.....	14
Starting the Program.....	14
Connecting the Device to your Computer.....	16
Demographics.....	17
SOC questionnaire.....	20
Birth Information.....	22
Chapter 2.....	24
Testing.....	24
Subfiles.....	26
Menu bar.....	27
Auto Focus Zap.....	28
Blue info grid.....	29
Chain Reaction Function.....	34
Columns in the main matrix.....	37
Main Matrix.....	38
Manfly.....	39
VARHOPE Score Interpretation.....	41
Alarm reaction.....	44
Chapter 3.....	45
Risks Profile – Mineral Panel.....	45
Risks Profile – Virtual Doctor Test.....	50
Risks Profile.....	51
Short Sarcodes Feedback Rx.....	53
Chapter 4.....	70

TABLE OF CONTENTS

Bodyviewer Training Program Activation.....	70
Iridology.....	75
Chapter 5.....	78
Disease Probability and Dictionary.....	78
Chapter 6.....	81
Close Out of Programs.....	81
Shut Down Your Device.....	82
Check Current Rectifications.....	82
Close Out of ED.X and NUCLEUS SOFTWARE.....	83
Closing a session with your client.....	83
Chapter 7.....	88
Modifying the SOC Index of a Previous Client.....	88
Loading a Previous Client.....	89
Modifying Birth Date.....	91
Chapter 8.....	91
Auto Frequency Therapy / Auto VARHOPE.....	91
Auto Meridian / Hyper-Reactants.....	93
Auto Scalar.....	95
Auto Trivector.....	96
Spiritual Resonance.....	97
Trivector and Hololinguistic Therapy.....	98
Two Suggested Training Programs.....	99
Auto Colour Therapy.....	100
Chapter 9.....	101
Superlearning.....	101
Timed Therapies.....	102
Chapter 10.....	103
ECG Fast Fourier Transform.....	103
Superconscious Reduction Panel.....	104
Autonomic Nerval System.....	105

TABLE OF CONTENTS

Chapter 11.....	108
Bones, Nerves, Circulation, Lymph, Digestion.....	108
Spinal Muscular Re-Education, Sarcode Feedback.....	109
Timed Biofeedback Therapies.....	111
Chapter 12.....	112
Hormonal Training Program.....	112
Internal Enzymes.....	113
Nutrition and Homeopathy Feedback.....	114
Stress and Cortisol.....	115
Chapter 13.....	115
Stimulate Emotional Growth Stress Reduction.....	115
Chapter 14.....	116
Other Biofeedback Choices.....	116
Reduce Stress in Nerves.....	117
Chapter 15.....	118
Body Scan and Face Therapy.....	118
Appendix A.....	119
EDUCATIONAL LANGUAGE.....	119
TERMINOLOGY for PRACTITIONERS.....	122
Appendix B.....	122
Electrode Data Check Box Does Not Appear.....	123
BASIC TROUBLESHOOTING.....	123
Problems With Graphics or Multimedia.....	124
Biofeedback Systemic Relaxation.....	124
Disease Probability And Dictionary.....	125
Hints And Warnings.....	125
Setting The Com Port.....	126
Special ED.X and NUCLEUS Programs Or AFE.....	127
Training Programs Prone To Lock Up.....	128
Tricks To Prevent And Deal With Freeze-Up.....	128

TABLE OF CONTENTS

Using Ctrl+Alt+Delete Tricks..... 128

What to do when computer screen is frozen..... 129

Updates and chnages..... 129

Brain wave pattern therapy..... 129

Changes made with Update no.43 & 44..... 132

Changes made with Update no.45 & 46..... 138

Changes made with Update no. 46..... 145

Chapter 1

Name

Step 1: NAME

Type in the name of the *Client*.

If doing a Practice Session:

for Client Name, type the words *Test Client 1*.

(for later sessions, you may type in Test Client 2, etc.)

If doing a Client Session:

for Client Name, type the full name (first, last and middle name if desired)

Click *Male* or *Female* or *Both*


Calibration

Now that we have our client set up in demographics, we are here in the main panel where the calibration option is now bold.



The Calibration process is our first hand-shake or connection with the device. This is the process whereby VARHOPE numbers are calculated.

Click **Calibration** (from the *Main Menu*)

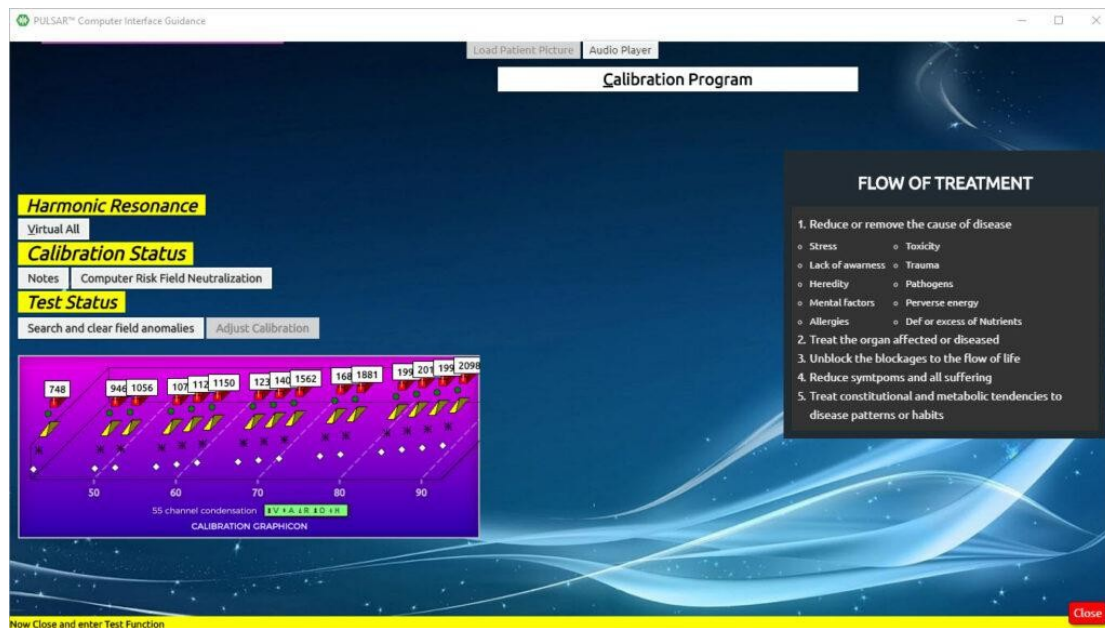


On this panel, you also have the option to choose **VIRTUAL ALL**. This disables the harness from the patient, especially in the case you work on a person with a pacemaker.


In the option **Computer Risk Field Neutralization**, you can check your room for geopathic stress, and clear it if you wish. This is a good idea especially when you start or finish your day at your office, or even between clients

Other options when you can click **Computer Risk Field Neutralization**:

- Neutralize and diminish computer stress
- Check and clear area/room for Geopathic Stress – it is suggested to start your day or end your day with this, or at least to do it every week, or when working in a new area
- The Schuman Field and Subspace protection field can be activated and deactivated here

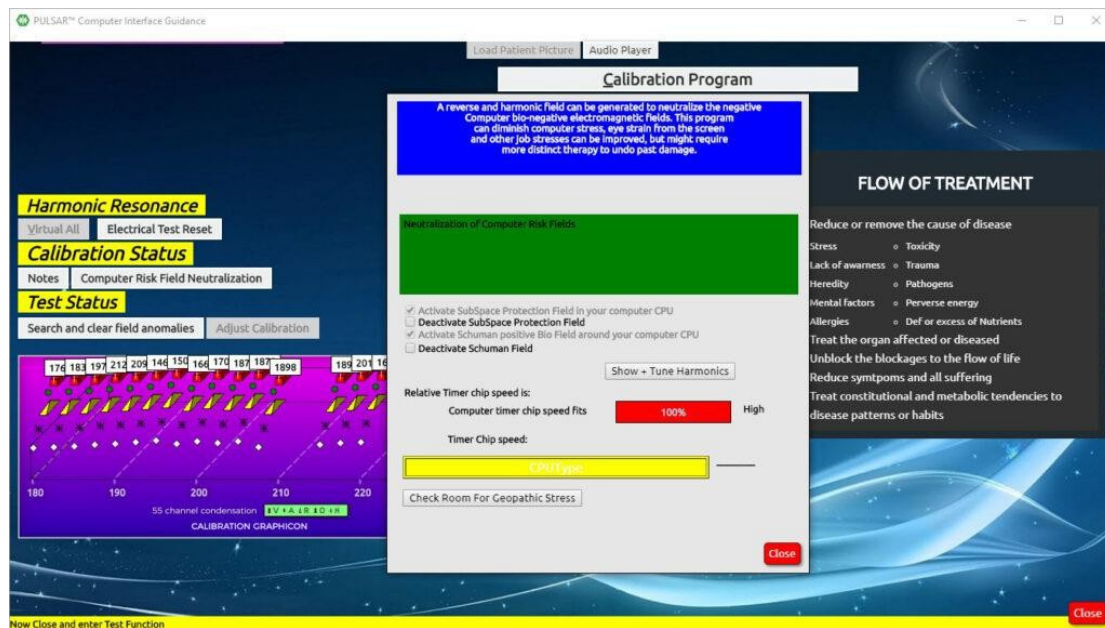

Click close on the **Computer Risk Field Neutralization**

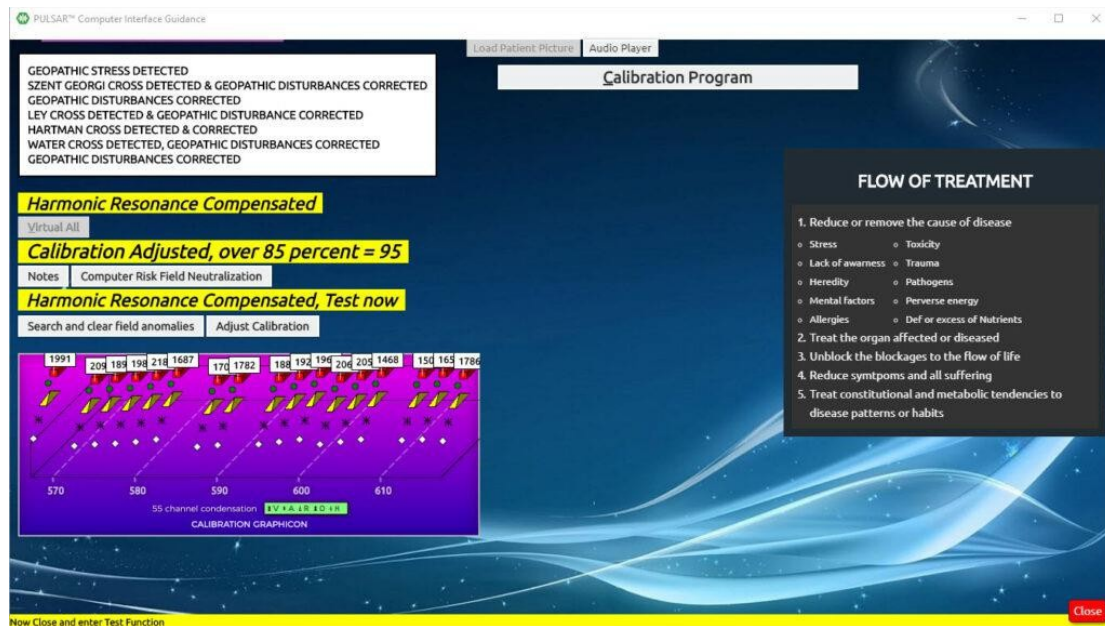


Now you may either proceed with the **Calibration Program** directly, or you may Load the **Client's Picture** first if you did not do it in **Demographics**.

We also have the option of loading in a picture (**Load Patient Picture**)

The PULSAR software loads the picture of the client much easier than previous devices, and there is no need to alter the picture in any way. All you need is a saved .jpeg photo.

Click on **Load Patient Picture**


This will automatically load the patient's picture in the background of the calibration interface

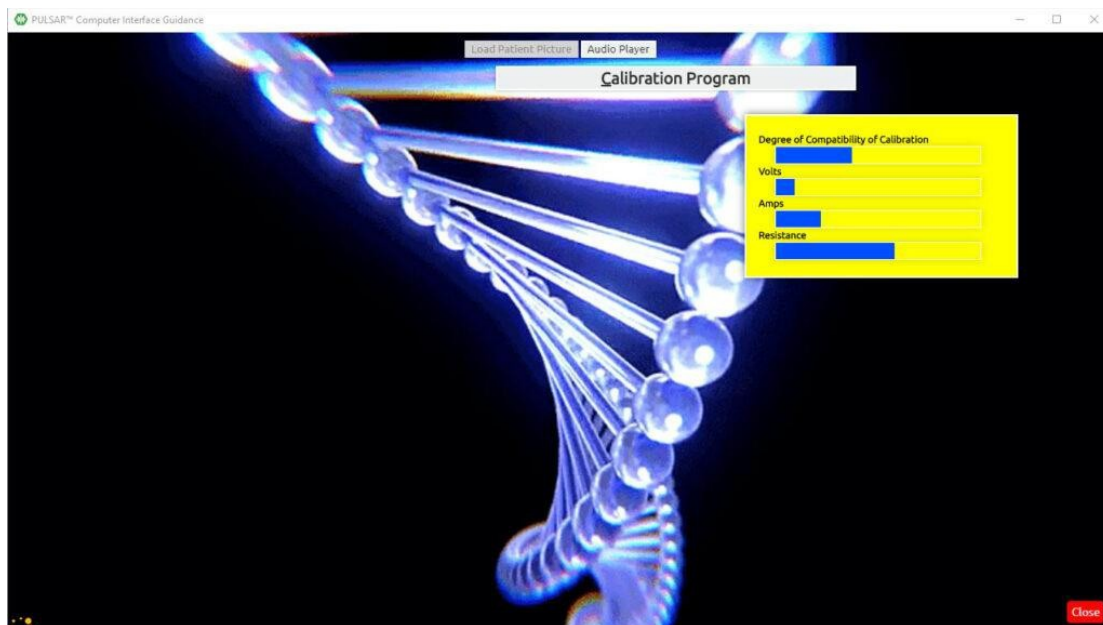

Click (Calibration Program)

Calibration begins

If a message comes up you just read it and convey it to your client and Click OK

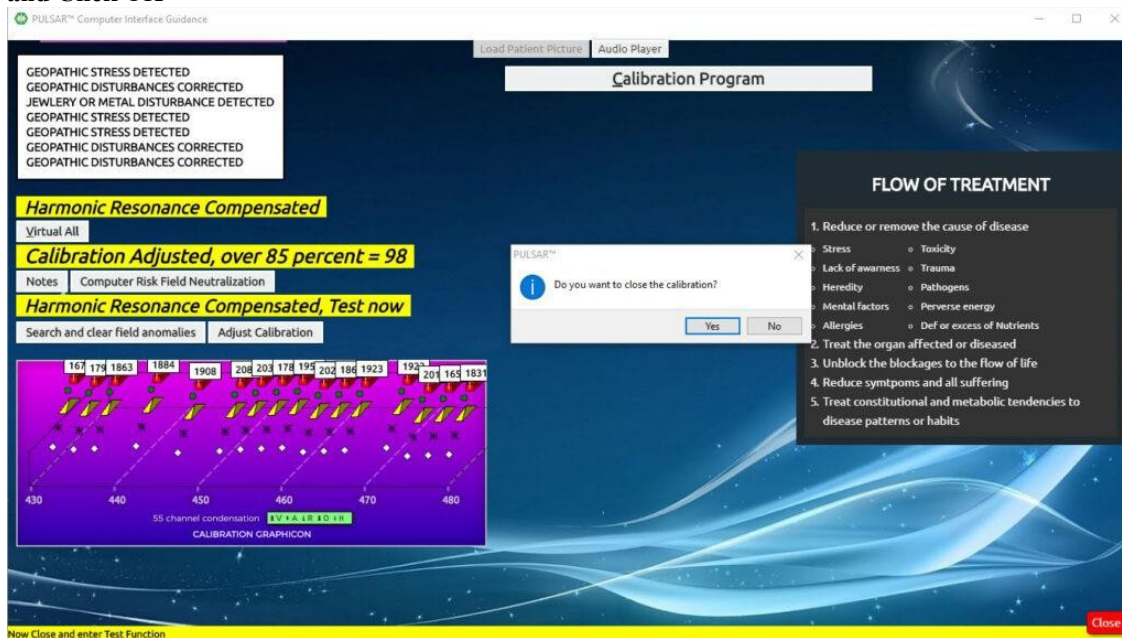

Click Yes (Do you want to close the calibration?)



And you are returned to the Main Panel where you can continue.

Harness Application and Usage

Have the client remove all jewelry not worn at night, loose change, keys, and cell phone.

Attach the tray, head harness, and limb straps to the device accordingly in the back of the device matching the colors of each cable.

Headband – Wrap the head harness around the client’s head, with the White cord to the Right side of the head and the eight electrodes on the forehead. If you wish, you may use a wet single-ply tissue between the forehead and the headband. It needs to be comfortable but not too tight.

Limb Straps – Attach limb straps around wrists and ankles, in a comfortably snug fit. The wires should be connected in the following manner:

Wrists – **RED** on Right – **YELLOW** on Left

Ankles – **BLUE** on Right – **BLACK** on Left

- **Do not apply over broken skin**
- **Do not** use on clients with pacemakers
- **Do not** use on clients with extreme electrical allergy
- **Be cautious** with clients who have had electric shock therapy

- **Wipe and clean** the harness after each use (consider a microfiber cloth as it is naturally disinfecting without chemicals). You can also use **diluted Everclear** or **diluted Vodka** as they leave no chemical residue. **Do not use Isopropyl Alcohol** as it dries out the harness and leaves a chemical residue.
- Surrogate testing can be done on infants with the harness connected to the mother and the child in the mother's arms
- If you suspect some difficulties with the client, do not connect the harness to the client, but have them place their hand on the tray or work in virtual.
- For clients who :
 - Are Pregnant
 - Have a history of epilepsy
 - Who has pacemaker

Please do not do therapy on them until you have acquired enough experience and have become a Biofeedback Specialist, in order to have a complete understanding of the type of protocols that can be followed.


- Note:

Attempting to open the box will interfere with the operation and void any warranty on your device.

Home Use

All buttons on the Main Menu panel are now bold and ready to open. At this point, you have the option to open any of the buttons, but more commonly, you would go to either **Home Use** or **Test**.

The **Home Use** is a program that is automatic and can use unconscious choices based on the calibration.

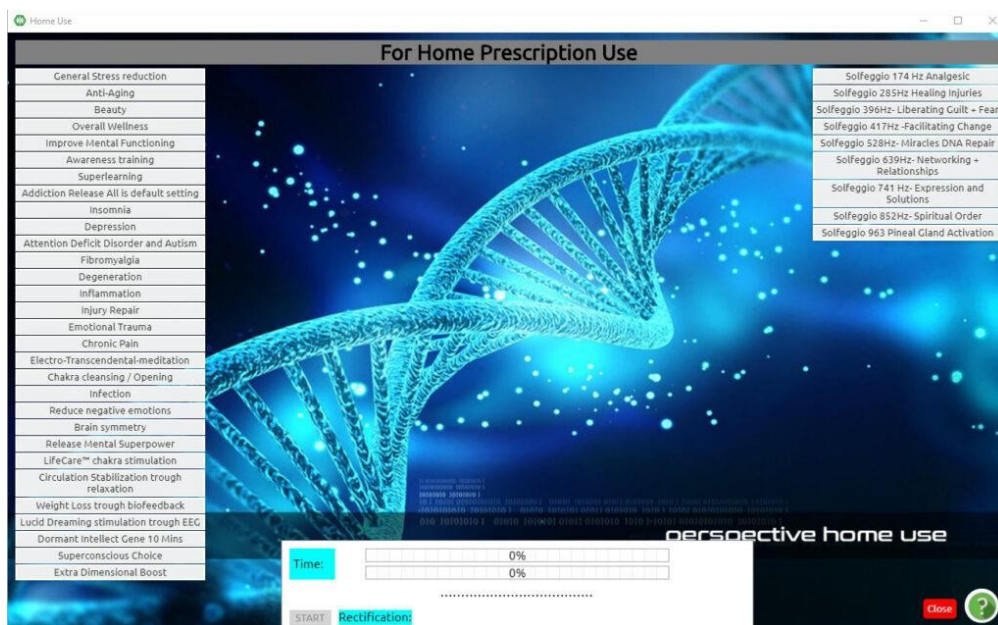

Click **HOME USE** ONCE

In order to select the program, you can:
click on **SuperConscious Choice** a few times (on the bottom of the list of programs), so that the PULSAR shows you the unconscious needs of the client.

Or you can select an individual therapy of your choice based on your professional experience.
The **Home Use** is a program that is automatic and can use unconscious choices based on the calibration.



By default the timer for each therapy is set at 1 minute, however, you may change the time according to what you think is best. The maximum time per individual therapy is 20 minutes.

The total time of the programmed therapies is shown in blue at the bottom part of the panel. The maximum of programmed therapies allowed is one and a half hours.

When you are ready with your selections and times, **Click on Start** on the white bottom part of the Home Use panel.



While the program is working, in the bottom white part of the panel you will be able to see which therapy is being done, how much percent and time are left to finish the therapy and the average percent of rectification of all programs.

Upon completion of each therapy, a percentage is shown on the side of the timer.

The **Home Use** is a program that is automatic and can use unconscious choices based on the calibration.



Upon completion of the selected therapies, you will see rectification results beside each selection and an average rectification result at the bottom.

If you wish to continue with more therapies you can click on other individual programs once to make them bold.

To de-select a therapy, all you need to do is click on the selection once and the items will be de-selected.

Click Close

Here you can either click OK to proceed to the EEG/ECG panel or Cancel to proceed to TEST.

PULSAR demographics part 1

From the MAIN MENU, press on DEMOGRAPHICS. (This will access the demographics panel).

Start by clicking once on the option NEW PATIENT. (That enables us to enter new details for a new patient).

Important: The client can be a person, an animal an inanimate object, a group of people, land, etc, depending on your scope of practice.

At the NAME, introduce the name of the PATIENT. For the gender, you shall select according to what you were working on, male or female, or if it is an inanimate object, a group of people with both male and female, or depending on physical or personal identity, or expression, you may choose **BOTH**.

Fill the rest of the fields with specific information: Date of birth, Place of Birth, and Country. The remaining fields are optional.

From the menu bar (on the top of the window) select: **BIRTH TIME AND PLACE**, then EDIT BIRTH DATA. (IF you are unsure of birth time, you can calculate. It's usually quite precise)

From the **GEOGRAPHY** tab, select EUROPE, and double-click on the country of birth (on the showing map). A window will pop up to tell you DATA RECEIVED, then press OK. This will send you back to the main demographics panel. From here, press **SAVE CURRENT CLIENT**. A window will pop up then press OK.


Starting the Program

Note: freeze ups may occur if clicks are too rapid, if you double click when only one click is needed, or if you click when and where no clicks are needed, with experience and practice you will avoid this from happening. If this happens wait and see if it clears, or you may have to restart. Please keep in mind that the software needs its time from program to program, do not rush.

DBL Click the PULSAR Icon on your desktop.

It displays as *ED.X* and *NUCLEUS* on your Task Bar.

This brings up the *Welcome Screen*.

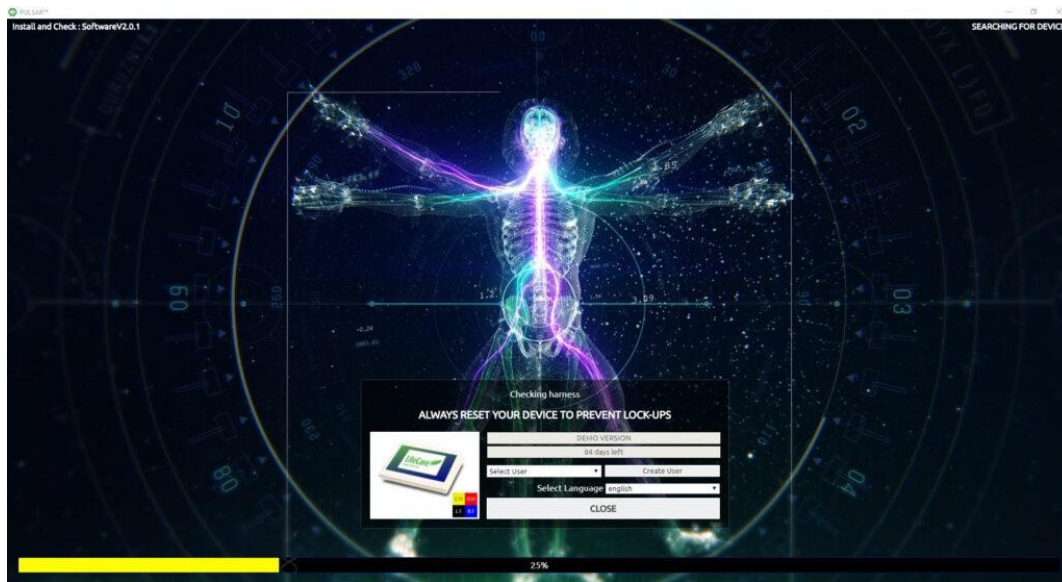

At this point, you can also choose your language preference at the bottom center of the Welcome Screen with a drop-down button.

At 25% of loading, you will see the Electro Data check appear on the screen.

If this panel does **not** appear then it indicates that the *ED.X* or *NUCLEUS* device is still searching.

If it continues searching, please make sure that the *ED.X* or *NUCLEUS* device is on, and that the USB cable is connected, otherwise please check **Appendix B Troubleshooting** for help.

Click **Close** on the Electro Data check panel and the program will continue to load.


The **Main Menu** screen now appears:



Click Password

Type 0210 into the password box

Click OK

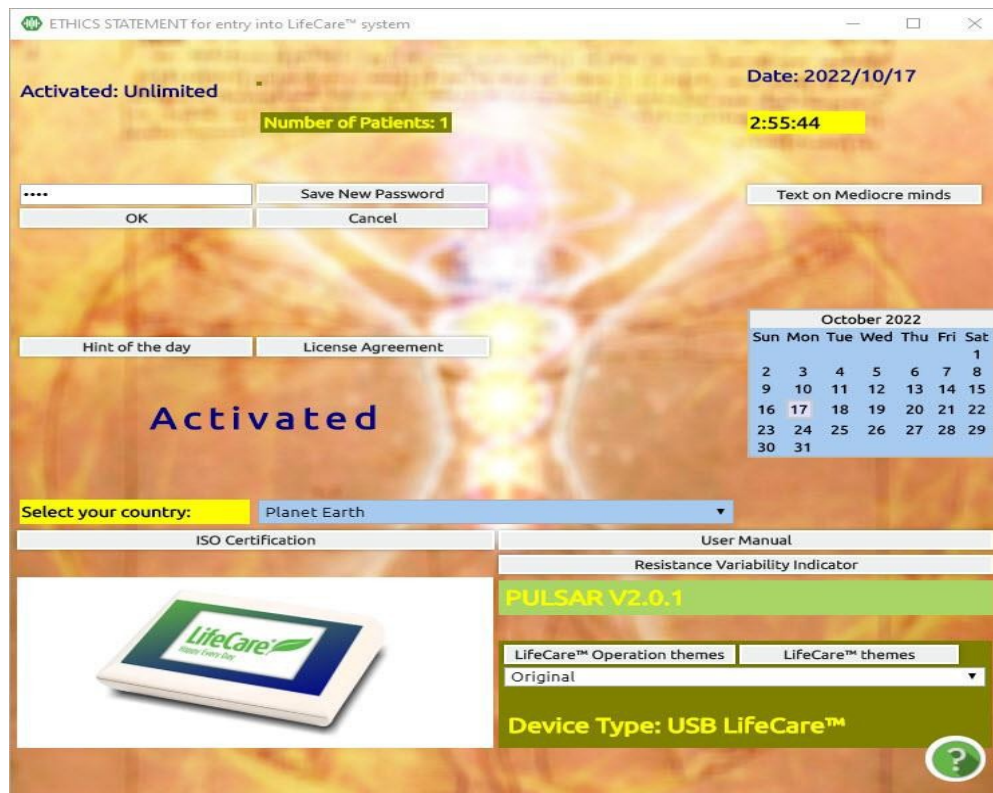

Click OK

Note: from here on in the software, you will find various text buttons that you can read for educational purposes. On this panel, you will see (as in the above picture) Text on the License Agreement, Text on Mediocre minds, the user manual, and the troubleshooting manual.

Connecting the Device to your Computer

1. You have received a USB3 type A male to type B male cable in your Starter Pack.
2. Insert one end of the cable to the back of the device.
3. Insert the other end into the designated USB3 port on your computer.
4. Note: It is a good idea that the laptop you will use for your new device is specifically for this usage so that you have better results and output. Your computer must be configured to communicate with the device.
5. Download the PULSAR from <https://quantum-medical.eu/index.php/downloads/> after you have connected your device to the computer, turn the device on. After that, you can run the downloaded installer. For more detailed instructions click here: <https://quantum-medical.eu/wp-content/uploads/2022/04/PULSAR-Installation-and-activation-instructions-1.pdf>
6. If you have issues with installation or activation, please contact the technical department of Quantum Medical SRL for the configuration of your laptop and activation of your software.
7. Turn your laptop on.
8. When the start-up of Windows on your laptop has been completed, push the **POWER button** on the back of your ED.X and NUCLEUS device to **ON**.
9. You will see graphics appear on the top panel of the device where it will also appear ED.X or NUCLEUS is on.

Demographics

Click **Demographics** (from the **Main Menu**)

While learning to navigate, it is important to practice on Test Clients (fictitious) in order to navigate freely. This enables us to explore the program without the need to complete each Training Protocol properly and without concern for imbalances.

Click **New Patient**

Cursor will go directly to the patient data so you can enter details.

Note:

For previous already entered patient

Click the name of the patient in the list (or you can search for a past client by pressing the button on the left to search for your client, whereby you can then click the name of the patient)

And the relative client shall be loaded

Click the small close button (shown above)

(please refer to Chapter 7 in this manual to load/modify previously existing clients)

Step 1: NAME

Type in the name of the **Client**.

If doing a Practice Session:

for Client Name, type the words **Test Client 1**.

(for later sessions, you may type in Test Client 2, etc.)

If doing a Client Session:

for Client Name, type the full name (first, last and middle name if desired)

Click Male or Female or Both

Step 2: Birth Information

Fill in the Birth information:

For a Practice Session (Test Client):

Date of Birth: you can put today's date

Place of Birth: Type in any City, Province, or State

Country: Type in any country (wherever you are is fine)

There is no need to fill in the Address

For a Client Session: Enter your own (if working with yourself) or your client's date of birth

Make sure to enter the date in the same format in which your computer is set:

Month/day/year 03/19/2008

Or-

Day/month/year 19/03/2008

Place of Birth: Type in the City, Province or State or other necessary identifiers Country: Enter country of birth

Note:

In order to achieve the correct calibration of the individual, it is important to enter the client's date of birth and place of birth with accuracy.

Patients

Today is: 2022/10/17 Visit # 1

Name: John Doe Use Gender specific history: General complaints First visit:

Date of Birth: 19/03/2008 ☐ Male ☐ Female

Place of Birth: Amsterdam Subscribed Disorders:

Country: Netherlands Amount of negativity in 100:

Address: City/Country/ST: Tel/Fax: List under this name: Search for next client:

Street: Zip Code: Insurance Company:

Buttons: New Patient, Modify Patient Data, Save Current Patient, Delete Patient, Cancel, Hide Patients

id	name	birthDate	birthPlace	country	sex	examCount	soc	birthTime	linked	creationDate
1	John Doe	19/03/2008	Amsterdam	Netherlands	Male	1	100	19/03/2008	1	19/03/2008

SOC INDEX: 100 John Doe to be tested

Device used:

What part of your body do you use this device on?

What positive things did you observe for?


How would your life change if your device was gone?

Buttons: Calculate SOC, List Patient Picture, Delete Patient Picture, Add more Data to Report, Report scores, Bio rhythms, Add Symptom Profile, Check Body Fat Percentage, Print interview help, Linkless inventory, Home testing, Close

- **Select Birth Time and Place in the top menu bar, and select Edit Birth Data**

If the time of birth is known, please enter it in the space provided where it says Enter Birth Time in Nautical Time Otherwise, if it is unknown, Select If unsure of Birth Time, calculate, and enter the time which will appear in the designated box under the latitude marker in the space provided for entering **Birth Time in Nautical Time**

- **Select the drop-down button for Geography, and select the continent in which the client was born**

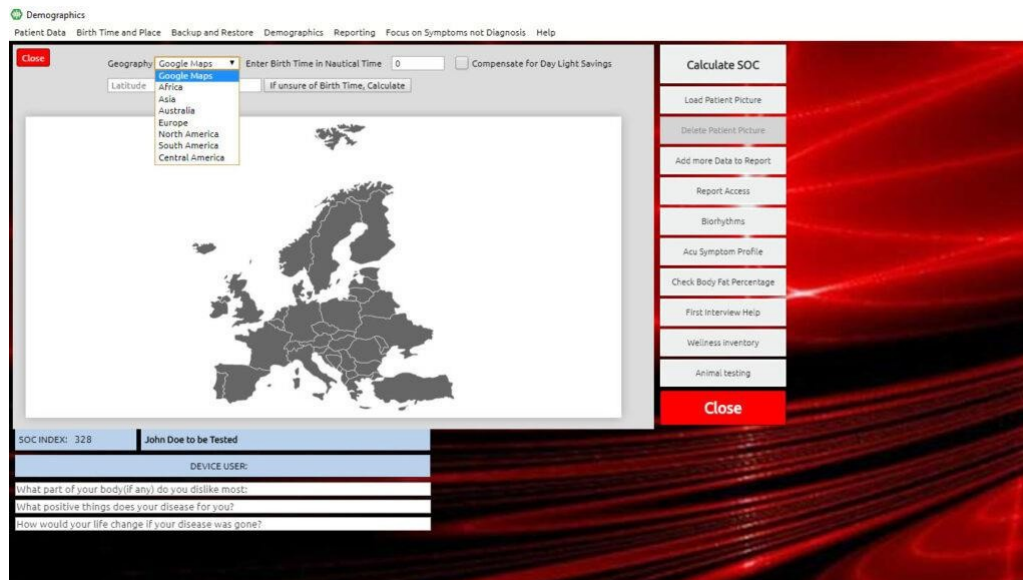

DBL Click on the area of the country in which the client was born

Click OK

The panel with the Demographics of the client details will appear. **Click the small Close**

Step 3: SOC INDEX QUESTIONNAIRE

The SOC questionnaire is done for educational and not medical purposes. It is used to raise client awareness and acquire an understanding of any areas of stress. When the SOC information is entered, the PULSAR program can now individuate to the specifics of the client. This ensures that the program will be more specific to the client as opposed to running generically. Certain answers may prompt other buttons to appear. Filling in the SOC will also provide an important opportunity to dialogue, educate and develop a rapport with the client. Here we are building a relationship with the client.

In further training, you will learn more about the details of filling in SOC details. For our purposes, fill in the SOC details as the client answers you accordingly.

Click the big Close button

If a pop-up box displays a *Biorhythm Concern for Today*
Click **No**

Note:

See Chapter 7 to learn how to return to a client you have *previously* entered in your Data List or how to *modify* a client's demographic information

Patient data is automatically saved when closing.

SOC questionnaire

The SOC questionnaire is done for educational and not medical purposes. It is used to raise client awareness and acquire an understanding of any areas of stress. When the SOC information is entered, the PULSAR program can now individuate to the specifics of the client. This ensures that the program will be more specific to the client as opposed to running generically. Certain answers may prompt other buttons to appear. Filling in the SOC will also provide an important opportunity to dialogue, educate and develop a rapport with the client. Here we are building a relationship with the client.

In further training, you will learn more about the details of filling in SOC

details. For our purposes, fill in the SOC details as the client answers you

accordingly.

Click the big **Close** button

If a pop-up box displays a *Biorhythm Concern for Today*

Click **No**

Note:

See Chapter 7 to learn how to return to a client you have previously entered in your Data List or how to modify a client's demographic information

Patient data is automatically saved when closing.

NOTE:

The importance of filling in the data entry of the SOC, is for the device to operate at its full capacity.

Entering SOC data is required by the system, as the entered values of each factor listed can affect results. High or low entries depending on the factor, which turns red, alert the system's power settings and safety features. The fields pictured in white do not affect the system.

The system calculates and takes note of the numeric entries entered to alert the system accordingly or to compensate for the factors. For example, the number of injuries a client or patient lists may require an adjustment to the level of resistance energy the system uses in order to ensure both safety and efficiency. Another example such as the number of glasses of water, if left blank, the system will artificially lower the general resonant frequency setting for the client for reasons of safety and optimal results will be reduced.

This panel is also an opportunity for the client or patient to be educated about lifestyle factors that are the client's responsibility to remain aware of and maintain for optimal results.

Hovering over the data entry fields reveals hints:

Number of organs removed	Don't forget tonsils, adenoids, appendix, uterus or ovaries, and one kidney. Only the majority of teeth removed would count as an organ removed.
Number of synthetic drugs	Count all taken, even OTC which stands for "over the counter" medications. Some medications can suppress symptoms and might reduce biofeedback reactivity strength to significant items, such as pain signatures. Users who are unlicensed in medicine seeking certification with a scope of stress management, however, are prohibited from advising on medications.
Amount of times you smoke	Use daily average over the last two months
Number of steroid-type drugs	Count again if in the second category

Number of amalgam fillings	Any metal fillings in teeth
Number of street drugs	Total number of different types use double # if excessive
Number of all known allergies	List of all major confirmed allergies
Number of unresolved mental factors	Anger, greed, desire, sadness, resentment, and fear are examples
Amount of fat in the diet yields	Most diets are 40%, and 20% is ideal. Note: this is one field that uses a scale of 1-10 to represent 10%-100%. Therefore if the client in diet intake is 30%, you enter 3.
Personal Stress 0-10	Max is the number of most tolerated, higher is possible
Number of Sugar Products per day	Men should not have over 3, women over 2
Number of cups of coffee per day	Caffeine is in chocolate and other foods too
Number of extreme toxic exposures in the past year	Count each chemotherapy and radiation treatment as well as any accidental exposures
Number of major injuries in the past year	Count all emotional, physical, or other trauma
Number of major infections	Count all major health-threatening infections
Number of glasses of water per day	Does not reveal a hint
Number of kilos overweight	Window to Patients self-view is designed for patients to answer

For more specific concerns you may click on More Input or Add More Data. Please note that the input in these panels is not saved from session to session.

Birth Information

Fill in the Birth information:

- For a Practice Session (Test Client):
- Date of Birth: you can put today's date
- Place of Birth: Type in any City, Province, or State
- Country: Type in any country (wherever you are is fine)
- There is no need to fill in the Address
- For a Client Session: Enter your own (if working with yourself) or your client's date of birth
- Place of Birth: Type in the City, Province or State or other necessary identifiers
- Country: Enter country of birth

Note:

In order to achieve the correct calibration of the individual, it is important to enter the client's date of birth and place of birth with accuracy.

Demographics

Patient Data Birth Time and Place Backup and Restore Demographics Reporting Focus on Symptoms not Diagnosis Help

Patients

Today is: 2022/10/17 Visit # - 0 Number of Patients:1

Name: Use Gender specific Matrix: ☐

Date of Birth: ☐ Male ☐ Female

Place of Birth: ☐ Pregnant

Country: Inherited Disorder: ☐

Amount of negativity 0-10: BMI:

General complaint first visit:

Address

City/Country/St. List order By: Search for past Client

Street: Tel/Fax:

Zip Code: Insurance Company:

- Select **Birth Time and Place** in the top menu bar, and select **Edit Birth Data**
- If the time of birth is known, please enter it in the space provided where it says **Enter Birth Time** in Nautical Time
- Otherwise, if it is unknown, **Select If unsure of Birth Time**, calculate, and enter the time which will appear in the designated box under the latitude marker in the space provided for entering **Birth Time in Nautical Time**

Demographics

Patient Data Birth Time and Place Backup and Restore Demographics Reporting Focus on Symptoms not Diagnosis Help

Geography: Enter Birth Time in Nautical Time: ☐ Compensate for Day Light Savings

Latitude: ☐ If unsure of Birth Time, Calculate

Africa
 Asia
 Australia
 Europe
 North America
 South America
 Central America



SOC INDEX: 328 | John Doe to be Tested

DEVICE USER:

What part of your body (if any) do you dislike most:

What positive things does your disease for you?

How would your life change if your disease was gone?

- Select the drop-down button for Geography, and select the continent in which the client was born



- DBL Click on the area of the country in which the client was born
- **Click OK**
- The panel with the **Demographics** of the client details will appear. Click the small Close

Chapter 2

Testing

This process will reveal the body's response to over 10,000 frequencies, showing potential weaknesses and strengths.

After the Calibration process mentioned above, we are ready to proceed with Testing.


Click *Test* (from the ***Main Menu***)



After clicking **TEST**, you enter the panel where you can prepare for Testing. In this panel, you will notice the VARHOPE panel in red, which are the first reactions that come up after Calibration. Please refer to the VARHOPE INTERPRETATION section for more information.

Click **Prepare** [top right and center of the screen (by the arrows)]

As soon as **Prepare** is pressed, a small electrical pulse is sent to the client to prepare them for testing. Once preparation is done, this button now gives you the prompt for **TEST**.

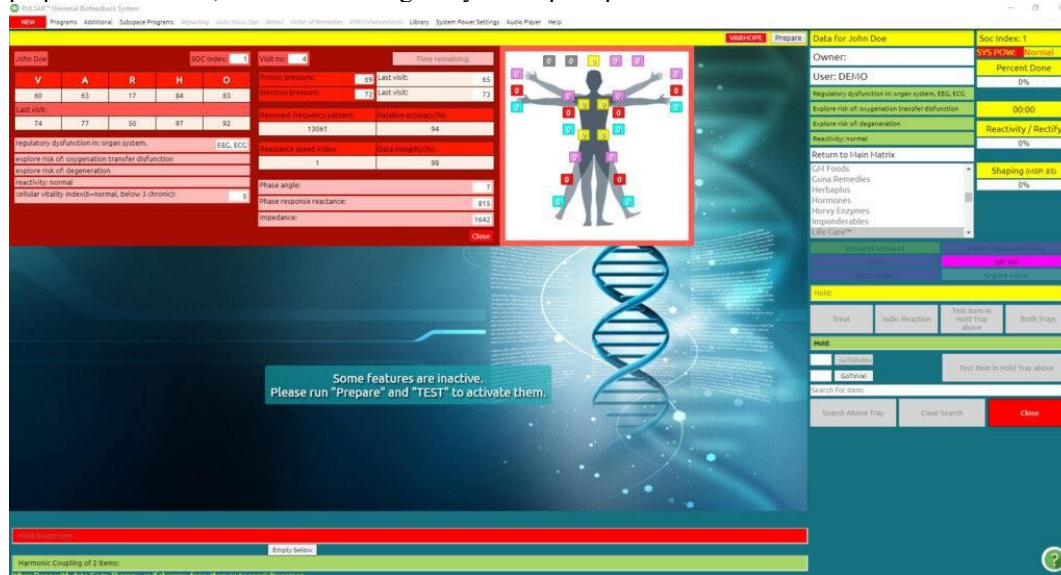


Note:

During the testing process, the client should remain relaxed, and still and refrain from speaking. It is also important to refrain from crossing arms, legs, or interlacing fingers during the entire session; this helps to keep the meridian flow open.

The word **TEST** now appears in red.


Click **TEST**

WAIT for the test to complete, it takes approximately 3 minutes. There is a yellow box in the upper right corner of the VARHOPE panel that tells you the amount of time left for the testing process.

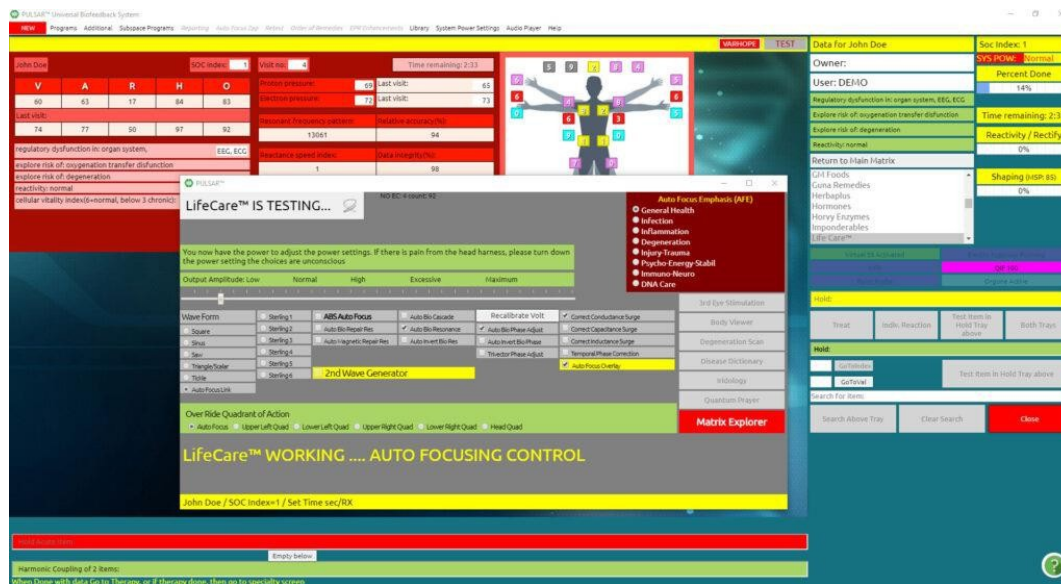

When Testing is over the screen is automatically transferred to the main matrix that shows all the results from the test, and the blue info grid (“Info Report”) and the **VARHOPE** panel are also automatically opened. You can close or reopen them at any time with the buttons above the Matrix called “Info Report” and “VARHOPE”.

You can place the “**Info Report**” with your mouse, anywhere you want on the screen. It is “draggable”. The VARHOPE Panel is not draggable, it will stay at the top of the screen.

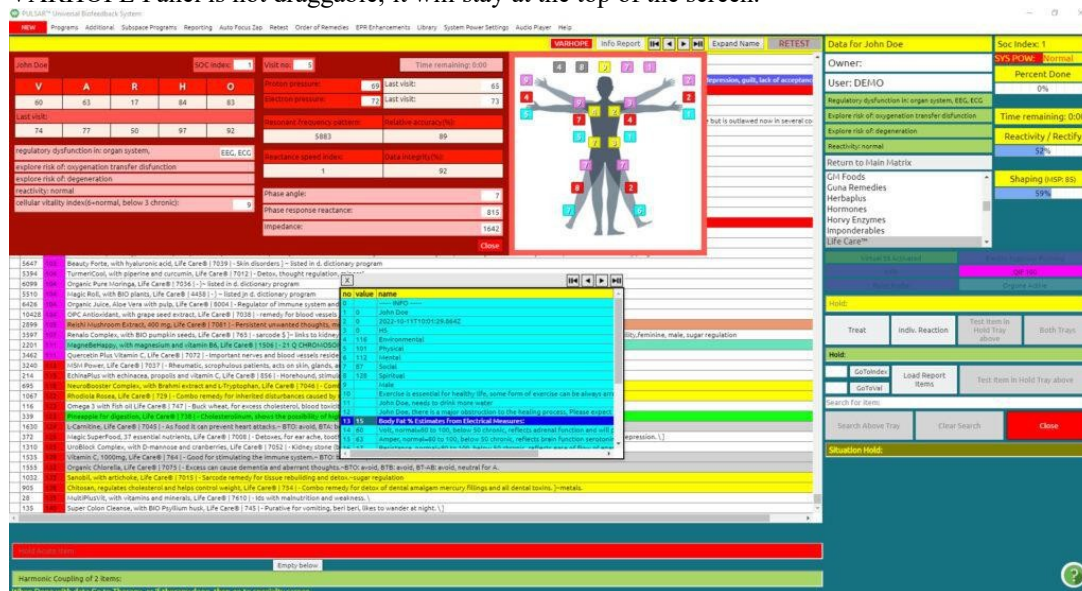

Subfiles

The *subfiles* are specifically *filtered* items from the *Main Matrix*. They are grouped into different *categories* in which all items are ranked by the degree of reactivity. Subfiles are located on the right side of the main Matrix page.

The **Main Matrix** menu bar contains various options for programs that have drop-down menus offering many choices thereafter. The main options are:

- NEW
- Program
- Additional
- Subspace Programs
- Reporting
- Auto Focus Zap
- Retest
- Order of remedies
- EPR Enhancements
- Library
- System Power Settings
- Audio Player
- Help

Auto Focus Zap

From the Main Matrix screen:

Click *Auto Focus Zap* (top menu bar)

Click on any *function* from the list

[illegible]

This will begin a Training Protocol for items having the top reactivity values for each particular category. Most of these functions last about 5 minutes, with the exception of ***Zap top Pathogens***, which will engage a 12-minute training protocol.

Note:

Remember to click **Return to Main** (first item under the Subfiles) before initiating any **Auto Focus Zap** Training: to return to the original full **Main Matrix** results.

After you have entered the **Risk Profile** (Chapter 3) and closed with the bottom close button **Zap Close**, and thereafter go to this drop-down menu, you will see an “x” mark in the most needed autofocus bio resonance zap training. These “X” marks give you a hint in general of also what is most needed for the client.

Blue info grid

The highest numbers indicate the general area of reactivity for that day.

The highest numbers indicate the general area of reactivity for today.

You can also retrieve this panel by simply clicking on the “Info Report” button above the Matrix, or:

Click **Reporting** (from the main matrix drop-down menu)

Click **Show INFO**

The blue info box can also be utilized as a “system tray function” for multiple Zapping, or a “Piggyback field” for entraining multiple items together.

Click **Reporting** (from the top bar in the **Main Matrix**) drop-down menu.

Click **Empty INFO to start a new one.**

The screenshot shows the PULSAR Universal Biofeedback System interface. The main window is a table of items, with a blue 'Show INFO' box open over item 10114, 'Life Impulse® ECO Pure Noni (8115)'. The right sidebar contains a 'Data for John Doe' section with various status indicators like 'SYS POW: Normal', 'Percent Done: 0%', and 'Time remaining: 0:00'. Below this is a 'Situation Hold' section with buttons for 'Treat', 'Indiv. Reaction', 'Test item in Hold Tray', and 'Both Trays'.

A prompt will come up:

“Are You Sure? Empty INFO Table?”

Click **OK**

The screenshot shows the PULSAR Universal Biofeedback System interface after the 'Empty INFO' prompt. The main window shows the 'Empty Info Table' dialog box. The right sidebar is the same as the previous screenshot, but the 'Data for John Doe' section now shows 'SYS POW: Normal' and 'Percent Done: 0%'. The 'Situation Hold' section is also visible.

The **INFO** box shall empty.

You can now **drag** and **drop** anything you wish in the INFO box in order to make a cocktail therapy. If you are working as a system tray function for multiple Zapping, it is a good idea to add up to a maximum of eight items in the blue INFO box.

You can repeat this process as many times as you wish.

On the other hand, when working with the blue INFO box as a “Piggyback field” for entering multiple items together, you can add as many items as you wish.

The screenshot shows the PULSAR Universal Biofeedback System interface. The main window is a table of remedies with columns for index, value, and name. The table is color-coded with red, yellow, and green rows. On the right side, there is a 'Data for John Doe' panel with various status indicators like 'SVS PCW: Normal', 'Percent Done: 0%', and 'Time remaining: 0:00'. Below this, there are buttons for 'Treat', 'Indiv. Reaction', 'Test item in Hold Tray', and 'Both Treat'. At the bottom, there is a 'Situation Held' section and a 'Harmonic Coupling of 2 items' section.

Option:

From the Reporting drop-down menu, you have the option to also load the top 200 remedies from the matrix to the INFO box, or the bottom 100 items.

The screenshot displays the PULSAR Universal Biofeedback System interface. The main window is titled 'PULSAR® Universal Biofeedback System'. It features a top menu bar with options like 'New', 'Programs', 'Additional', 'Substance Programs', 'Reporting', 'Auto Focus Zap', 'Reheat', 'Order of Remedies', 'BPR Enhancements', 'Library', 'System Power Settings', 'Audio Player', and 'Help'. The central area is divided into three main sections: a 'Library' on the left, an 'INFO Report' in the center, and a 'Data for John Doe' panel on the right. The 'Library' lists various items with their indices, values, and names. Item 1077, 'Life Impulse® B12 Pure (100%)', is selected and highlighted in blue. The 'INFO Report' for item 1077 is displayed in the center, showing detailed information about the item, including its description, ingredients, and usage instructions. The 'Data for John Doe' panel on the right contains fields for patient information, including name, age, sex, and test results. Below this panel are buttons for 'Treat', 'Test', and 'Test item in Hold Tray above'. At the bottom of the interface is a 'Harmonic Coupling of 2 Items' section with a text area and a 'Empty below' button.

Once you click on the “load report” button, you will notice on the purple tray “Report items loaded” and now the button shows “unload items”

This means that the items in the blue INFO box are now loaded for

- Treating,
- Zapping,
- Testing (individual reaction of the group of items),
- Test item of hold tray above (that is the report items loaded in the blue info box) or piggyback training

You may click from the **Auto Focus Zap drop-down menu**, Auto Focus Zap 3 min, or Auto Focus Zap 9 min. This will switch the Treat button on the side to either ZAP 3 min or ZAP 9 min.

By having the report loaded, you may either Treat (which is a few seconds), ZAP 3 Min, or ZAP 9 min., or just the fact that the items are loaded, they work as piggyback training during your entire session.

Additionally, you may choose to do an Individual Reaction on these loaded items, it will give you an average profile of the group of items tested together in the INFO box, Click Test Report.

Here you can interpret Results, you have the option to do a 30-sec zap, or you can do a 3-minute zap. We have added a “Custom time” for zapping where you can set the length of the gap between 1 and 15 minutes, on the slider that appears on the right side of your screen:

Hold:

ZAP 1 min	Indiv. Reaction	Test Item in Hold Tray above	Both Trays
-----------	-----------------	------------------------------	------------

Hold


Go

[illegible]

As mentioned in the purple tray, the button **Test Items in Hold Tray above** is also activated so as to test against filters. This will give you more background information concerning the items you are testing.


Chain Reaction Function

The Chain Reaction Function is one of the latest features that the ED.X and NUCLEUS present. It assists greatly in the therapy, making the therapy right from the start more focused on what you wish to work on, getting to the basic cause and issue of the person and the age at which the issue began, very quickly, precisely, and only within a few seconds after testing.

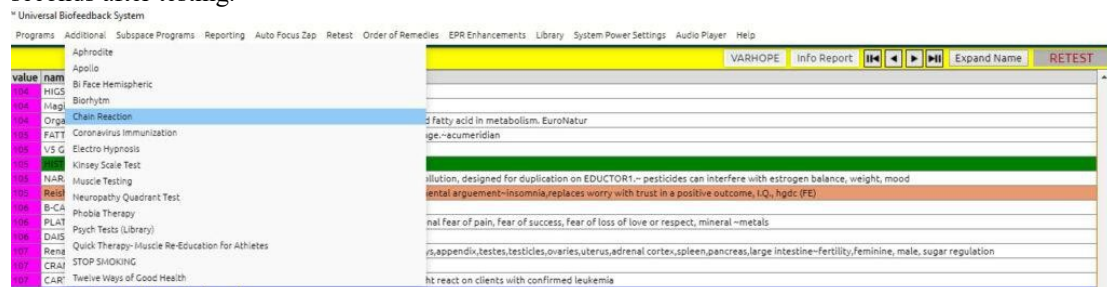

From your main program, you first want to search for the item of the main concern/issue that your client may have. Let's say in our example it is eczema.

Put the word eczema in the search for item field

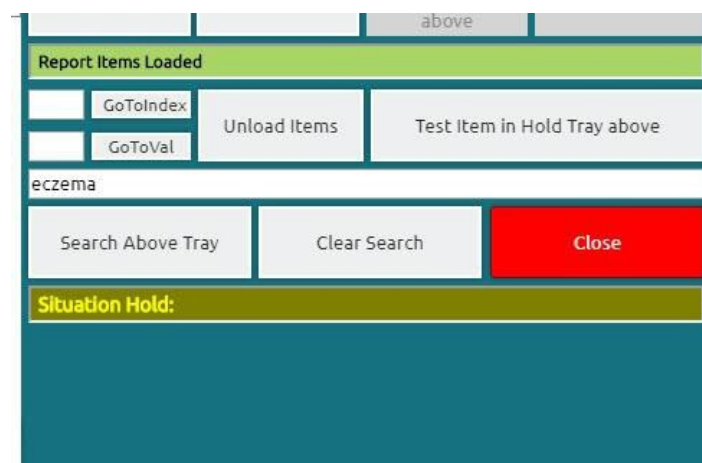

And
Search above the tray

The screenshot shows the PULSAR Universal Biofeedback System interface. On the left, a list of items is displayed with their reactivity values. The highest reactivity value is 99 for 'eczema'. On the right, a panel shows the search results for 'eczema' and a list of items with their reactivity values. The highest reactivity value is 99 for 'eczema'.

By searching the area of concern (in our example “eczema”), items reacting concerning our issue are filtered in the matrix. We notice the item that comes up with the highest reactivity value is 99 in our example. Please remember or keep this number in thought.

The next step you would want to do is
Return to Main. From the subfilters on the side

The screenshot shows the PULSAR Universal Biofeedback System interface. A red arrow points to the 'Return to Main Matrix' button in the subfilters on the side.

This will take you back to the main matrix panel.

The above panel appears on top of the blue info grid panel.

Put the highest reactivity of the item of concern in the third box where it mentions the Highest concern value. In our case, I will put 99 as per our above eczema example.

Automatically it will tell me that the main issue now for the area concern (eczema) began at the age of 12 and that we can find the causes of our issue within the reactivity values of 32. Which appear in the area below.

At this point, we can look at all the items with a reactivity value of 32 to see what our main issue was at that period in the person’s life. Additionally, we can discuss with the client possible difficulties that occurred at that specific age of 12 as per our example.

VARHOPE

Info Report

Expand Name

RETEST

Month old:

503

Close

Chain index:

3.59

Highest concern value:

99

Reactivity number:

27.58

Years, at chain reaction started:

18.24

Add to Report

index	value	name
7812	118	=ResonVal/React=27 / Complex of minerals and trace elem
2330	119	=ResonVal/React=27 / WAVE A2 Composite, dental isode
2940	121	=ResonVal/React=27 / GRANULOMA PURULENT Growth of
5830	122	=ResonVal/React=27 / MYOPIA ~ listed in d. dictionary prog
7106	123	=ResonVal/React=27 / CLOTOSOL Dental Materials Tempor
1423	123	=ResonVal/React=27 / HISTIDASE Enzyme involved in anti-in
4286	124	=ResonVal/React=27 / Empty
829	124	=ResonVal/React=27 / EASTERN OAK Allergy or sensitivity t
6654	125	=ResonVal/React=27 / DENTOT0X16 Post surgery. GUNA s.r.
3319	125	=ResonVal/React=27 / LM 21 NEUROLOGICAL N Sarcose S
4681	126	=ResonVal/React=27 / ISOPOGON (FE)4681 Able to learn fr
8759	127	=ResonVal/React=27 / Schisandra, Energizing, Life Care® 75
958	127	=ResonVal/React=27 / Complex with vitamin D3 (2000 IU), K2

Add to report

This will add all the items with a reactivity value of 26 into the blue info grid box to work as a piggyback field therapy throughout your session.

Load report items on the middle side of your screen so that you load your items in your blue info grid box to your piggyback therapies. : Report items loaded” will appear on your purple tray on the side.

VARHOPE

Info Report

Expand Name

RETEST

Month old:

503

Close

Chain index:

3.59

Highest concern value:

99

Reactivity number:

27.58

Years, at chain reaction started:

18.24

Add to Report

index	value	name
7812	118	=ResonVal/React=27 / Complex of minerals and trace elem
2330	119	=ResonVal/React=27 / WAVE A2 Composite, dental isode
2940	121	=ResonVal/React=27 / GRANULOMA PURULENT Growth of
5830	122	=ResonVal/React=27 / MYOPIA ~ listed in d. dictionary prog
7106	123	=ResonVal/React=27 / CLOTOSOL Dental Materials Tempor
1423	123	=ResonVal/React=27 / HISTIDASE Enzyme involved in anti-in
4286	124	=ResonVal/React=27 / Empty
829	124	=ResonVal/React=27 / EASTERN OAK Allergy or sensitivity t
6654	125	=ResonVal/React=27 / DENTOT0X16 Post surgery. GUNA s.r.
3319	125	=ResonVal/React=27 / LM 21 NEUROLOGICAL N Sarcose S
4681	126	=ResonVal/React=27 / ISOPOGON (FE)4681 Able to learn fr
8759	127	=ResonVal/React=27 / Schisandra, Energizing, Life Care® 75
958	127	=ResonVal/React=27 / Complex with vitamin D3 (2000 IU), K2

Owner:

SYS PCWME

Normal

User:

DEMO

Percent Done

0%

Time remaining

0:00

Reactivity/Rectify

52%

Shaping

59%

Report items loaded

Get Item

Unload Item

Test Item in Hold Tray

Search for Item

Search Above Tray

Clear Search

Close

Shedding Hold

So now you are ready!

This is an excellent tool to get focused on your therapy, right from the start, on the root cause of your issues. It helps you to start your discussions with your client concerning the age of when it all got started.

The **Value** column indicates the degree of reactivity to a specific item. These values are calculated during the **TEST** process and they change from session to session.

FUSAR™ Universal Biofeedback System

These values are calculated during the **TEST** and they are unique for each session and client, depending on their reactivity at that specific time of testing.

Keep in mind that apart from the electrical reaction, **SOC** results also influence the values and reactions to the items in the matrix.

The following is a general guideline concerning the values:

Any number above 100 is significant, current, or cents

Below 50 reveals low reactivity, potential exhaustion, and **possible chronic or long-term stress.**

Red: These values indicate the most reactive items to which the person's energy field showed the most significant

Pink: These values indicate slightly less reactive items. (The reactivity is more than **two** standard deviations from the norm).

Yellow: These values indicate slightly fewer reactive items than the Pink ones. (The reactivity is more than **one** standard deviation from the norm).

Blue: These values indicate the reactivity with little or no response. Items above 100 may have significance.

The values are reactivity responses of the body to either present, past, or future needs or stressors. Keep in mind that reactions are not fixed, and do vary.

Before going to Individual Reactions, which is a more in-depth method for electrical investigation, it is advisable to go to the **Risks profile**, in order to establish priorities and interpret reactions. The Risk Profile represents a correlation of all reactions into categories.

Please keep in mind that a strong reaction to a main matrix item signature named after a disease is not a diagnosis of that disease, but rather a clue to compare against the rest of the matrix reaction profile.

Repetitions of strong or high reactions of items within the Main Matrix, of the same nature, can be an indicator of the area of focus for your session.

Main Matrix

Working directly from the main matrix:

From the **Main Matrix** screen:

Click **Item of Choice** from the **Main Matrix** (this will highlight the item)

Click **Treat (Train)**

This will deliver a seven-second frequency of the selected item (the button is on the right side of the screen, below the sub-files, under the yellow Hold Tray).

To select a longer frequency Training Protocol:

Click **Auto Focus Bio-Resonance Zap** (from the top menu bar)

Click **Auto Focus Zap 3 min**

The **Treat (Train)** button is now transformed to **Zap 3 min**

The screenshot shows the 'Universal Biofeedback System' interface. The top menu bar includes 'File', 'Programs', 'Additional', 'Subprograms', 'Reporting', 'Auto Focus Zap', 'Reset', 'Order of Remedies', 'BPR Enhancements', 'Library', 'System Power Settings', 'Audio Player', and 'Help'. The main window is divided into three sections: a list of items on the left, a detailed description of the selected item in the center, and a control panel on the right. The control panel includes buttons for 'Zap 3 min', 'Indiv. Reaction', 'Test Item in Hold Tray', and 'Both Trays'. The status bar at the bottom indicates 'Harmonic Coupling of 2 Items' and 'When done with data go to Therapy, or if therapy done, then go to specialty screen'.

— Or —

Auto Focus Zap 9 min

The Treat (Train) key is now transformed to Zap 9 min

Note:
Remember to return the **Zap 3 min or 9 min** button to its regular **Treat (Train)** button, when complete with the 3 or 9 min. selection. To do this:
Click **Auto Focus Bio-Resonance Zap** (top menu bar)
Click **Unzap**

[illegible]

Programs/Universal Biofeedback Therapy on the Main Menu


The screenshot displays the 'Universal Biofeedback System' interface. On the left, a list of programs is shown, including 'Biofeedback Acupuncture Therapy', 'Biofeedback Sensory Reeducation', 'Blood Analysis Library', 'Body Scan And Face Therapy', 'Cervical Sacral Reflexology', 'Dark Field Analysis Library', 'Dental - Oral Muscular Re-Education', 'Electro-Physiological Oscillation', 'Feedback Neurophysiological Reeducation', 'Hypnosis', 'Muscular Re-Education for athletes', 'NLP Emotional Growth Stress Reduction', 'Nutritional and Herbs Therapy', 'Pain-Quantitative Biofeedback Treatment', 'Rife Profile', 'Short Service Feedback Rife', 'Spinal Muscular Re-Education, Sensory Feedback', 'Stress + DNA Reactivity', 'Test reactivity (test only)', 'Time-Quantitative Feedback Therapy, Music/Superlearning', and 'Universal Biofeedback Therapy'.

The center panel shows the details for the 'Universal Biofeedback Therapy' program, including its description, a list of ingredients, and a list of related programs. The right panel displays the 'Data for John Doe' section, which includes a 'Bio Index 1' table with columns for 'Owner', 'User', 'Percent Done', 'Time remaining', 'Reactivity', and 'Rectify'. The table shows data for 'John Doe' and 'User: DEMO', with a 'Percent Done' of 74% and a 'Time remaining' of 0:00. Below the table, there are sections for 'Return to Main Matrix', 'Grid Foods', 'Grid Remedies', 'Herbaplus', 'Hormones', 'Honey Enzymes', 'Imponderables', and 'Life Care'.

Multi-Media on the bottom Task Bar


The screenshot displays the 'Universal Biofeedback Therapy' interface. The top section shows 'percent successful 0%' and a 'Report' button. Below this, there are sections for 'Test Time Spent', 'Start', and 'Text Explaining Rife and Bioresonance'. The 'TriVector and Holinguistic Therapy' section includes a 'Primary Disease' dropdown, an 'Organ relation' dropdown, and an 'Emotional Blockage' dropdown. The 'Rife Generator' section includes 'Upper Freq' and 'Lower Freq' dropdowns. The 'Scale Generator' section includes 'Upper Freq' and 'Lower Freq' dropdowns. The 'Bioresonance Bicom + Mora Like' section includes 'Amplify in Db' (0=Off), 'Set Program #', and 'Amplify healthy Freq', 'Amplify Disease Freq', 'Automatic Interaction', 'Wobble On', and 'Intermittent' options. The bottom task bar includes a 'Power Norm' button, a 'Double' button, a 'Spiritual Healing Resonance Panel' button, and a 'Multi-Media' button. The 'Multi-Media' button is highlighted in yellow.

On the bottom left of the *Universal Biofeedback Therapy* screen, posted in the yellow bar, you will see *Multi-Media*. A Pop up box will appear.

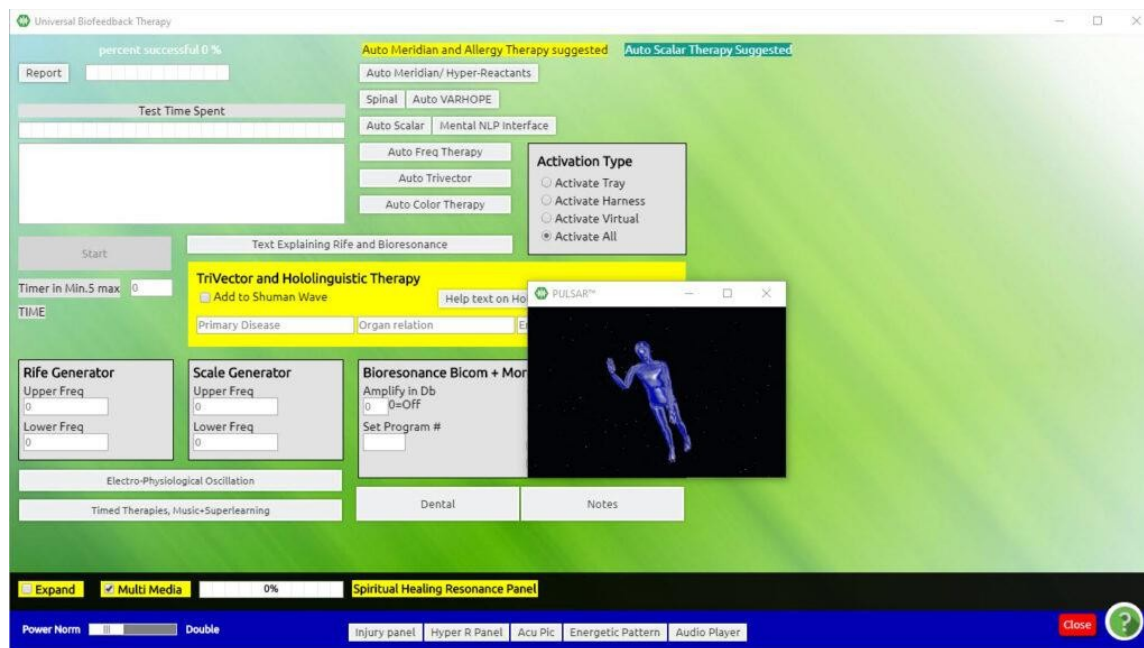

Now minimize the pop-up box by clicking anywhere on the therapy screen, so that it will run in the background.

VARHOPE Score Interpretation

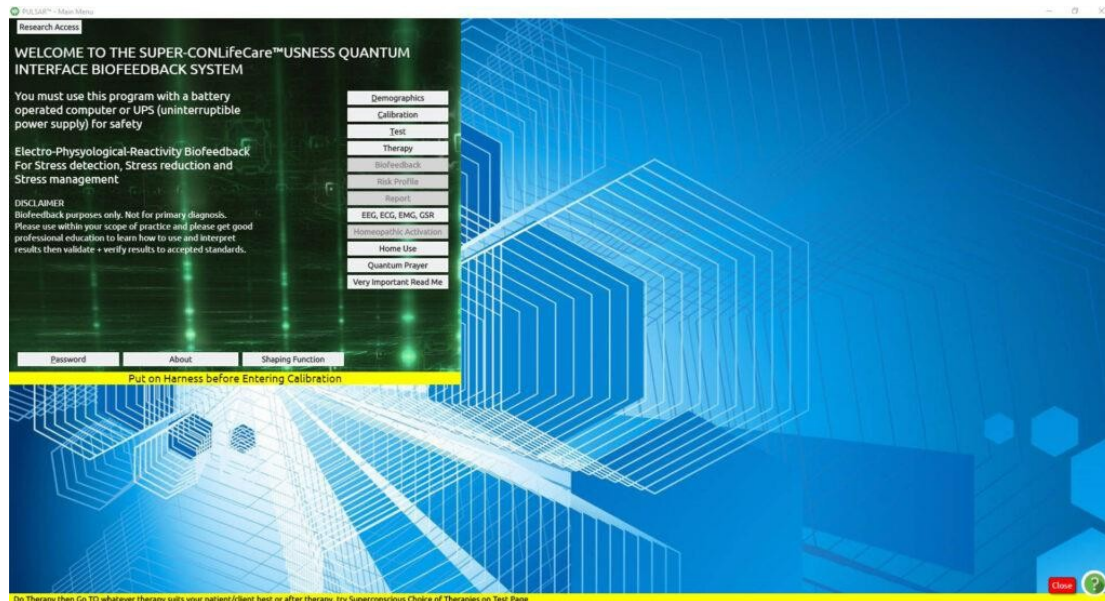

The VARHOPE panel first appears by clicking TEST in the main panel, after you have gone through the process of Calibration. They are the initial reactivity numbers after the first handshake of the Calibration.

The screenshot displays the VARHOPE software interface. On the left, a table lists various health metrics and their corresponding values. In the center, a human figure diagram shows energy flow points. On the right, a detailed report for 'John Doe' is visible, including personal information, system status, and a list of remedies.

V	A	R	H	O
60	45	17	84	93

Regulatory dysfunction in organ system: ECG, ECG
 Explore risk of: degeneration transfer dysfunction
 Explore risk of: degeneration
 reactivity: normal
 cellular vitality index(normal, below 3 chronic): 5

VARHOPE values (appearing in the red access panel) are for educational purposes. They are an indication of energetic parameters. Low VARHOPE values generally indicate an energetic weakness i.e. stress. These values may also give you information about the client's lifestyle habits.

You can also access the **VARHOPE** screen, from the **VARHOPE** button

The screenshot shows the VARHOPE screen with a list of remedies and their associated VARHOPE values. The table includes columns for index, value, and name.

index	value	name
1715	26	BLIND HBP 1 - Monastery of Herbs For high blood pressure/ circulation /
3274	25	BIO Lemon Essential Oil, Life Care® 4606 COLD AND FLU NKF
1922	26	GIARDIA (Flagellated parasite from bad water, produces dysentery (Stomach, colon, intestine) at first, affects liver later *) -digestion
1224	26	GIARDIA (Flagellated parasite from bad water, produces dysentery (Stomach, colon, intestine) at first, affects liver later *) -digestion
6471	26	MenstruEasy, with night light oils and vitamin E, Life Care® 1511 -Premenstrual Syndrome, Natureplant - Ortomolecular Vitamin J-feminine

The **VARHOPE** values (appearing in the red access panel) are for educational purposes. They are an indication of energetic parameters. Low **VARHOPE** values generally indicate an energetic weakness i.e. stress. These values may also give you information about the client's lifestyle habits.

Information such as client name, number of visits, and last visits VARHOPE values is also mentioned on this panel. There is also information on “regulatory dysfunction”, “explore the risk of” and “reactivity” fields that can give you hints and ideas on how you can plan your session.

The screenshot shows the VARHOPE screen with a list of remedies and their associated VARHOPE values. The table includes columns for index, value, and name.

index	value	name
1715	26	BLIND HBP 1 - Monastery of Herbs For high blood pressure/ circulation /
3274	25	BIO Lemon Essential Oil, Life Care® 4606 COLD AND FLU NKF
1922	26	GIARDIA (Flagellated parasite from bad water, produces dysentery (Stomach, colon, intestine) at first, affects liver later *) -digestion
1224	26	GIARDIA (Flagellated parasite from bad water, produces dysentery (Stomach, colon, intestine) at first, affects liver later *) -digestion
6471	26	MenstruEasy, with night light oils and vitamin E, Life Care® 1511 -Premenstrual Syndrome, Natureplant - Ortomolecular Vitamin J-feminine

Regulatory Dysfunction is a primary electrical concern that needs to be addressed that day. Common examples that may come up include spinal energy flow, EEG, ECG, and metabolic error dysfunction.

The Explore Risk items are secondary electrical concerns of that day that also need to be addressed.

Reactivity response ideally may be either normal, hyporeactive which is slow reactivity, under the normal, or hyper-reactivity indicating over the normal reactivity. Reactivity Dysfunction may indicate both hyper and hypo reactivity. In order to have a good understanding of the readings, it takes experience and good knowledge of each parameter and combination to look for.

Let us see however the Interpretation of the VARHOPE values as follows:

V = Voltage:

Energetically reflects stress related to adrenal function and willpower

This shows us possible chronic fatigue or fatigue in general.

The normal range is 80 – 100, Below 50 is chronic.

A = Amperage:

Energetically reflects stress related to brain function, serotonin, life force

The normal range is 80 – 100, Below 50 is chronic.

R = Resistance:

Energetically reflects overall stress related to ease of energy flow in the meridians of the body

The normal range is 80 – 100, Below 50 is chronic.

H = Hydration:

Energetically reflects stress related to the body's water absorption on a cellular level.

Not necessarily how much water one may drink but actually what they absorb.

The normal range is 80 – 100, Below 50 is chronic.

O = Oxygenation:

Energetically reflects stress related to the ease of oxygen flow in the body reflected in cellular vitality.

This can be interpreted into oxygen level, blood issues, hyperventilation, nutritional issues

The normal range is 80 – 100, Below 50 is chronic.

P = Proton Pressure:

Energetically reflects stress related to the body's energetic

pH. The normal range is 65 – 70.

The overall energetic pH is considered to be the most balanced when the numbers between Proton Pressure and Electron Pressure are within 10 points or less of each other.

E = Electron Pressure:

Energetically reflects stress related to the whole of the body's electrical functioning. This may be an energetic indication of the number of available electrons.

A good range is 70 – 75.

Cellular Vitality:

Associated with the present energy, as estimated electrically for the cell. Components of cell vitality are linked to oxygenation, protein digestion or blood issues/anemia, and essential fatty acid intake.

6 and above is normal, whereas below 3 is chronic.

Resonant Frequency Pattern:

Represents the reactivity of the nervous system, tri vector signatures, and stored stress such as from unprocessed emotions, aging, sugar intake, or a spiritual issue.

Below 1.000 is associated with a risk of degeneration, and above 15.000 is possible nervous energy or psychic pain. When the Resonant Frequency Pattern is out of the normal range, it is recommended to train via the ECG and the ANS (Automatic Nervous System training programs).

Reactant Speed:

Represents the speed at which the body resounds, such as to pathogen tri vector signatures.

Below 15 is ideal, below 40 is normal, and above 100 is stressed.

Phase Angle:

Represents the ability of the body to regenerate cells as a cellular aging index.

8 plus is ideal, 6 is normal and below 3 is chronic.

Phase Response:

Represents the ability of all the systems of the body to work together.

900 plus is ideal, 700 is normal and below 300 is of concern.

Impedance:

Represents the ability of energy to flow around and through the

body. 1600 is ideal, below 1200 is weak and below 800 is of concern.

At the end of each session, you can go re-test the VARHOPE improvement made during the therapy.

From the main menu go to **Retest**

Click on **Test VARHOPE improvement**

The screenshot shows the VITA LIFE Universal Biofeedback System interface. The main window displays a list of 100 items with columns for index, value, and name. The right sidebar shows 'Data for John Doe' with fields for Owner, User, and various status indicators like 'SYS POWER' and 'Percent Done'. The bottom section shows 'Harmonic Coupling of 2 Items' with a 'Retest' button.

This will do a 1.5 min Retest and give you new values on your VARHOPE along with **Today's Improvement Percent**.

This percentage represents how much the person has reacted and improved electrically from their session.

The screenshot shows a detailed report for John Doe. The report includes a table of V, A, R, H, O values, a table of test results (Proton Pres, Electron pres, Res Freq Pat, Reactance Speed Index, Phase Angle, Phase Response Reactance, Impedance), and a diagram of a human body with numbered points (1-10) indicating specific areas of interest.

Please keep in mind that a low or even negative percentage value of improvement does not mean that the session was not good. However, it may mean that the patient may need some time, a day or two, for the body to self-adjust electrically and absorb the therapy that they received.

Alarm reaction

Please be aware that if an **Alarm Reaction** occurs at any time during your session, this simply indicates that a stress response needs to be addressed.

There are 5 key reasons that alarms may occur:

The relaxation and/or re-education of routes of elimination were not adequately supported, either from Short Sarcodes (Chapter 3), or Iridology and Body Viewer (Chapter 4).

Too many relaxations and muscle re-education training were performed without proper preparation.

The client's body has had enough training in that area and is done for the day.

The client may need more support with minerals, amino acids, and other nutrients.

You have received reactions from a level the client is not ready to work with today, such as from emotional stress.

If a client alarms excessively, the session may have exceeded the tolerance level of the client, and you should proceed with caution.

Double Click **on the red *Alarm Reaction* indicator (wherever it may show in any of the Training Programs)** when complete it will either indicate *Alarm, No More Therapy*, or *Alarm Response Corrected*.

If the response is corrected, continue the training protocol you were doing.

If the response is not corrected, discontinue the training protocol you were doing and move on to another area of the program.

Or-

If this panel does not have the possibility of addressing the *Alarm Reaction* training protocol, you will need to rectify it by returning to the *Main Matrix* screen.

On the *Main Matrix* screen, you will find a flashing *Alarm Reaction* button on the upper right side of the screen.

Double Click **Alarm Reaction**

Click **OK** in the pop-up box

In the center of the red *VARHOPE* Panel, the message *Correcting Alarm Response* will appear.

After a short training time, it will change to either *Alarm Response Corrected* or *Alarm Response Not Corrected*.

If the response is corrected, continue the training protocol you were doing.

If the response is not corrected, discontinue the Training Protocol you were doing and move on to another area of the program.

Double Click the background of the red *VARHOPE* panel to remove

If a box suggests: *Please use Patient (Client) Superconscious Panel*:

Click **OK**

(Refer to the Table of Contents for instructions on using the *Patient (Client) Superconscious Reduction Panel*)

Chapter 3

Risks Profile – Mineral Panel

RISKS PROFILE – ORIENTAL HERBS, MINERALS, AMINO ACIDS, AROMA THERAPY, PERVERSE ENERGY

The following Training Protocol can be performed for any of these programs: Minerals, Oriental Herbs, Amino Acids, Perverse Energy, Aromatherapy

From the *Risks Profile* screen:

Click on either *Oriental Herbs, Minerals, Amino Acids, Aroma Therapy, or Perverse Energy* (top of the panel)

Risk Profile

Treat: HORMONAL, LIVER, FUNGUS

BLOOD

Oriental Herbs

58	AntiHelminthic, parasites
121	AntiAsthmatic, regulates lungs, chills fever
149	AntiPhlogistic, excess mucous, wind cold dampness
81	AntiPyretic for fever, expels toxic heat, wind, chills
162	AntiRheumatic, rheumatic pain, activates blood flow
159	Blood, replenishes blood, restores vitality
88	Carminative, relieves gas and stagnant food
73	Digestive Stimulator, invigorates stomach, spleen
88	Expectorant, clears lungs, dispels heat in lungs
68	Female, regulates menstrual, vitalizes blood
165	Male, regulates virility, vitalizes blood
47	Hepatic, removes stagnation of liver + chi
68	Hypertensive, reduces heat in liver, reduces yang
38	Laxative, relieves deficiency of spleen fluid
99	Mental, invigorates Kidney meridian, treats senility
25	Sedative, nourishes blood in liver, treats fatigue
89	Bone, for degenerative bone conditions
46	Cold and Flu, stimulates chi and liver
41	Esophagus, relieves spasms and treats esoph cancer
89	Kidney, treats kidney deficiency and degeneration
110	Large Intestine, treats meridian weakness and bowel
47	Liver, treats chi and liver weakness
103	Lung, treats lung degeneration adds energy
98	Lymph, treats degeneration in the Lymphatics
26	Preventative, use to prevent cancer
108	Stomach, treats deficiency of stomach energy
150	Uterine, treats problems of uterine cellular degeneration

Double click on an item for individual reactance and resonance
Single click to add to Hold tray, remove extra text and leave product name.

Organ: Test Organ For Involvement

Reactance = 0
Resonance = 0
For Oligo use, best Potency: 1 (Double click to treat)

Treat:

Treat + Test Hold Items Make Report of edit

Make Report of above

Close

Rectified = 0

RESPIRATORY

Any item with reactivity over 100 is significant and should be treated.

Double Click on each high reactivity item (one at a time) until it is rectified to 85 or above.

– Or –

Click on a high reaction item, it will be added to the bottom hold tray (up to five may be added).

Click **Treat+Test Hold Items** to initiate the training protocol for all the items in the tray. Rectify to 85 or above.

Risk Profile

Treat: HORMONAL, LIVER, FUNGUS

BLOOD

Oriental Herbs

58	AntiHelminthic, parasites
121	AntiAsthmatic, regulates lungs, chills fever
149	AntiPhlogistic, excess mucous, wind cold dampness
81	AntiPyretic for fever, expels toxic heat, wind, chills
162	AntiRheumatic, rheumatic pain, activates blood flow
159	Blood, replenishes blood, restores vitality
88	Carminative, relieves gas and stagnant food
73	Digestive Stimulator, invigorates stomach, spleen
88	Expectorant, clears lungs, dispels heat in lungs
68	Female, regulates menstrual, vitalizes blood
165	Male, regulates virility, vitalizes blood
47	Hepatic, removes stagnation of liver + chi
68	Hypertensive, reduces heat in liver, reduces yang
38	Laxative, relieves deficiency of spleen fluid
99	Mental, invigorates Kidney meridian, treats senility
25	Sedative, nourishes blood in liver, treats fatigue
89	Bone, for degenerative bone conditions
46	Cold and Flu, stimulates chi and liver
41	Esophagus, relieves spasms and treats esoph cancer
89	Kidney, treats kidney deficiency and degeneration
110	Large Intestine, treats meridian weakness and bowel
47	Liver, treats chi and liver weakness
103	Lung, treats lung degeneration adds energy
98	Lymph, treats degeneration in the Lymphatics
26	Preventative, use to prevent cancer
108	Stomach, treats deficiency of stomach energy
150	Uterine, treats problems of uterine cellular degeneration

Double click on an item for individual reactance and resonance
Single click to add to Hold tray, remove extra text and leave product name.

Organ: Test Organ For Involvement

Reactance = 87
Resonance = 5
For Oligo use, best Potency: 1 (Double click to treat)

AntiHelmintic,, AntiAsthmatic, AntiPyretic for Fever, AntiPhlogistic

Treat + Test Hold Items Make Report of edit

Make Report of above

Close

Rectified = 24

RESPIRATORY

Click Close

TIP:

Minerals play an important role in the biological response to stress. Therefore, being adequately mineralized enables the client to respond not only to stress itself but in general, including Biofeedback.

Training with the mineral tri-vector signatures can often have a positive relaxation effect and is highly recommended in each session. If all or most mineral values are high, suspect:

Poor diet

Low stomach acid

Malabsorption of the small intestine, or a combination of all three.

You can test the Mineral reactivity to a specific Organ or gland sarcode tri vector signature by entering the organ name in the white box (Organ:).

The screenshot displays the 'Risk Profile' software window. The 'Treat' field is set to 'HORMONAL, LIVER, FUNGUS'. The 'BLOOD' tab is selected, showing a value of 138. The 'Minerals' section is highlighted in yellow and lists 27 elements with their associated functions. Below the list, there are instructions on how to use the interface, a text input field for 'Organ:', and buttons for 'Test Organ For Involvement', 'Make Report of above', 'Close', 'Treat + Test Hold Items', and 'Make Report of edit'. The 'Rectified' value is shown as 97. The 'RESPIRATORY' tab is visible at the bottom with a value of 78.

#	Name	Treat: HORMONAL, LIVER, FUNGUS
20	BLOOD	138
35	Minerals	
10		
36	44	Calcium (weak bones,nerves,adrenals)
31	31	Potassium (fatigue,nerves,energy,heart)
29	55	Sodium (depression,nerves,fatigue,digestion)
20	20	Magnesium (adrenal regulation, oxygen energy)
147	147	Sulphur (detox,energy,emotions,unfocused)
145	145	Manganese (nervous and muscle disease)
94	94	Chromium (sugar regulation, muscle)
34	31	Zinc (immunity,oxygen,metabolic disease)
39	39	Selenium (detox,nerves,energy,skin)
146	146	Iodine (thyroid function,energy,metabolism)
93	93	Phosphorous (cellular function,energy, thought)
82	82	Boron (nervous regulation, dementia, depression)
30	30	Molybdenum (thyroid function,lactation,fatigue)
159	159	Silicon (bone,skin,nerves)
142	142	Cobalt (anemia, immunity)
50	50	Lithium (nerves, thought, energy, hypothalamus function)
113	113	Germanium (nerves,oxygenation,skin,energy)
40	40	Arsenic (energy,nerves)
38	38	Antimony (skin,detox,bowel function)
109	109	Tin (nerves,muscles)
143	143	Carbon (energy,life regulation)
121	121	Vanadium (liver function, heart muscle)
129	129	Aluminum (nerves, thought regulation)
156	156	Copper (nerves,energy,detox)
23	23	Nickel(liver function, heart, muscle)
159	159	Silver (immunity,energy,detox)
92	92	Gold (rejuvenation,skin,depression)
77	77	Iron (anemia,immunity,energy)

Double click on an item for individual reactance and resonance
Single click to add to Hold tray, remove extra text and leave product name

Organ: Test Organ For Involvement Make Report of above Close

Reactance = 98
Resonance = 92
For Oligo use, best Potency: 4 (Double click to treat)

Calcium , Manganese , Moly Tin , Lithium , Molybdenum , Phosphorous , Zinc , Manganese , Iron

Treat + Test Hold Items Make Report of edit

Rectified = 97

RESPIRATORY 78

Risk Profile

#

Name

Treat: HORMONAL, LIVER, FUNGUS

20

BLOOD

136

10

Aroma Therapy

36

82

Balsam oil, antiseptic and expectorant

29

88

Basil aroma oil, aids digestion, antispasmodic

108

108

Cajeput aroma oil, gastrointestinal, pulmonary antimicrobial

22

84

Caraway aroma oil, after meals prevents flatulence, stims digestion

95

95

Cedar aroma oil, skin irritations, dermatosis, dry skin eczema

5

79

Cinnamon aroma oil, antiseptic deodorant, nausea, constipation

74

74

Clove aroma oil, antiseptic, toothache, indigestion, aphrodisiac

34

66

Coriander, aids digestion, prevents fermentation, bloating, circulation

97

97

Cypress aroma oil, astringent, contracts veins, varicose, legs feet

32

99

Eucalyptus aroma oil, antiseptic for respiratory tract, antiviral

66

66

Geranium aroma oil, eczema, acne, cuts, dermatosis, insect repellent

18

110

Grapefruit aroma oil, heart tonic, oily skin, cleans digestive tract

115

115

Green Anise aroma oil, digestive, bloating, intestinal tract

12

76

Juniper aroma oil, rheumatic pain, detox, diuretic, stims kidney

78

78

Lavender aroma oil, calming, detox, skin and nerve relaxant

8

118

Lemon aroma oil, stims immunity, weightloss

96

96

Mandarin aroma oil, epilepsy, antispasmodic, palpitations, relaxant

13

122

Orange aroma oil, relaxant, sedative for children

62

62

Oregano aroma oil, protects immune, sudorific

14

96

Peppermint aroma oil, nervous digestion calmer, NOT if breast feeding

38

91

Thyme aroma oil, intestinal, pulmonary, muscular pain

95

95

Rosemary aroma oil, tonifier, revitalizes liver gallbladder

26

109

Sandalwood aroma oil, urinary antiseptic, cystitis, prostate

85

85

Tea Tree aroma oil, skin, detox, infection, sinus

28

68

Vervain aroma oil, relaxant, antiseptic, dissolve kidney stones

63

63

Ylang Ylang aroma oil, hypertension, aphrodisiac

17

0

—

4

0

—

27

Double click on an item for individual reactance and resonance

40

Single click to add to Hold tray, remove extra text and leave product name

3

Organ:

Test Organ For Involvement

15

Reactance = 89

Resonance = 98

7

For Oligo use, best Potency: 4 (Double click to treat)

31

Balsam oil, Cajeput aroma oil

2

Treat + Test Hold Items

Make Report of edit

24

Rectified = 116

37

RESPIRATORY

78

Oriental Herbs and Aroma Therapies are used for prevention. High reactivity values indicate strong reactions that would be good to consider treating.

Risk Profile

Name Treat: HORMONAL, LIVER, FUNGUS

BLOOD

Amino Acids

64	Phenylalanine, (pain control, nerves)
86	Alanine, (kidney and nerves)
101	Asparagine, (nutrasweet, neural destruction)
61	Cysteine, (kidney utilization)
75	Leucine, (mood control, emotions)
56	Isoleucine, (emotional control, nerves)
43	Serine, (carbohydrate energy conversion)
92	Tryptophan, (serotonin and relaxation)
80	Histidine, (antiinflammation, antiallergy)
94	Methionine, (liver and oxygenation)
99	Lysine, (anti herpes lesion, nerves)
75	Threonine, (hormonal and energy production)
44	Tyrosine, (thyroid, pit and adrenal function)
89	Valine, (for regulation of blood cells)
76	Glutamine, (supplies energy to brain)
32	Proline, (sulphur disorders)
101	Arginine, (nerves and skin)
103	Glutamic Acid (supplies energy to brain)
96	Aspartic Acid, (nutrasweet, neural destruction)
44	Adenosine, (circulation, energy)
58	Uracil, (for RNA function)
56	Adenine, (for DNA RNA function)
76	Guanine, (for DNA RNA function)
92	Cytosine, (for DNA RNA function)
43	Thymine, (for DNA RNA function)
83	Tryptophan, (mood control, energy, healing, sleep)
104	Taurine, (energy)
104	Citrulline, (thinking)
105	GABA, (thought control)

Double click on an item for individual reactance and resonance
Single click to add to Hold tray, remove extra text and leave product name

Organ: Test Organ For Involvement

Reactance = 148
Resonance = 40
For Oligo use, best Potency: 3 (Double click to treat)
Phenylalanine, Leucine

Treat + Test Hold Items Make Report of edit

Make Report of above

Close


Rectified = 32

RESPIRATORY

Amino Acids are compounds that play many critical roles in the body and are often referred to as the “building blocks of proteins”. They are necessary for vital processes like building proteins and synthesizing hormones and neurotransmitters.

Reactions to Amino Acids may show us reactivities and negative impacts on our nervous, reproductive, immune, and digestive systems.

Treating high reactivity numbers will assist in reducing stress in these areas.

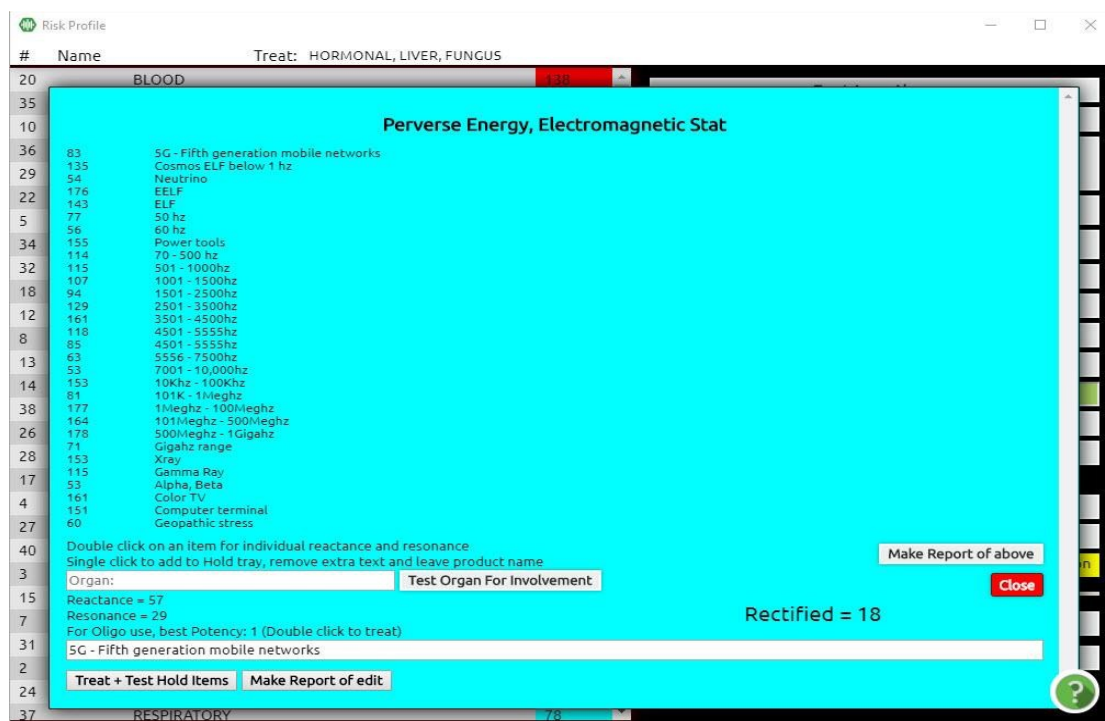

Reactivity values on *Perverse Energy, and ElectroMag Stat* indicates stress in these areas. Discuss with the client about possible exposures.

Due to the high rate of electromagnetic smog exposure, it is always a good idea to go through this panel to reduce your client's stress from radiation exposure.

Risks Profile – Virtual Doctor Test

By choosing any item from this list and then clicking the *Virtual Doctor Test*, you can see a set of additional reactive tri-vector signatures of items in relationship to the “risk” item tested.

This panel will show you the *Primary and Secondary Causes and Cure* related to the stress/risk item you have chosen, as well as the *Possible most stressed organ* (most reactive sarcode tri-vector signature).

You will also have the option to begin training from this panel by clicking on the reactive items that are displayed.

The following Training Protocol can be repeated for any item in the *Risks Profile*.

Click on any high reactivity (red) item to initiate a training protocol. The item that will be tested will appear in the purple tray in the middle of the panel.

Click Virtual Doctor Test

A new screen will appear indicating the Primary and Secondary (energetic) causes and cures (remedies) related to the selected item from the Risk Profile panel, as well as the Possible most stressed organ.



Click on any line for a 3-second training protocol; repeat until **Rectified** is 85 or above.

– 0r –

Click ***Treat All on the Panel*** to begin a training protocol for each item on the screen.

Risks Profile

This program will reveal specific issues that are categorized for your client. Retraining can be performed here. However, before considering the individual test results, it is preferable to build a foundation or framework of the basic issues and priorities. This panel will be able to guide your investigation of the individual reactions.

From the *Main Matrix* screen:

Click **Program** (top menu bar)

Click ***Risks Profile*** (drop-down list)

The screenshot shows the Sensory software interface. On the left, there is a list of items with columns for #, Name, and a numerical value. The list includes items like BLOOD, SENSORY, VIT.DEF OR EXCESS, CHOLESTEROL, EMOTIONAL, NEUROLOGICAL, VIRUS, ALLERGY, FOOD POISONING, BONE, LYMPHATIC, LIVER, TRAUMA, OXIDATION, HORMONAL, RADIATION, FUNGUS, HYPOADRENIA, CONNECTIVE TISSUE, BACTERIA, COGNITION, STRESS, KIDNEY, AMOEBA, TOXICITY, RESPIRATORY, DEGENERATION, and CIRCULATION. The right side of the interface shows a detailed view of a selected item, including a list of symptoms, a list of items, and a list of items with a numerical value. The interface is titled 'SENSORY' and has a menu bar with options like File, Edit, View, and Help.

Double Click the item with the highest reactivity in the left column and it will be added to the top hold tray. Do the same for the next 4 highest reactants (thus adding a maximum of 5 categories to the tray).

The screenshot shows the Sensory software interface. On the left, there is a list of items with columns for #, Name, and a numerical value. The list includes items like BLOOD, SENSORY, VIT.DEF OR EXCESS, CHOLESTEROL, EMOTIONAL, NEUROLOGICAL, VIRUS, ALLERGY, FOOD POISONING, BONE, LYMPHATIC, LIVER, TRAUMA, OXIDATION, HORMONAL, RADIATION, FUNGUS, HYPOADRENIA, CONNECTIVE TISSUE, BACTERIA, COGNITION, STRESS, KIDNEY, AMOEBA, TOXICITY, RESPIRATORY, DEGENERATION, and CIRCULATION. The right side of the interface shows a detailed view of a selected item, including a list of symptoms, a list of items, and a list of items with a numerical value. The interface is titled 'SENSORY' and has a menu bar with options like File, Edit, View, and Help.

Click **Treat Areas Above** to initiate a Training Protocol for these specific Risk factors. Pay particular attention to where your client is responding on the **Selye Bar**, as this will determine how they are reacting to their stressors. Ideally, the client would be in the “**Alarm**” phase, meaning that they are responding to stress with strong reactivity. The “**Adaption**” phase suggests that the client may be getting used to the stress and is not responding with strength today, or possibly symptoms are not appearing.

The “**Exhaustion**” phase can suggest fatigue, or that the client may be depleted and no longer able to respond to stress.

Adaption and Exhaustion phases are the basis of a chronic stress pattern. These phases can be free of medically diagnosable symptoms, as the body has no energy left to respond. A good example is a chronic fatigue.

Short Sarcode Feedback Rx

The Short Sarcode Feedback Rx program is a very useful application giving you many options, depending on what you want to work on.

Note: One of the main reasons why you may want to go into this program first, before doing anything else, is to load into your report basic findings from the testing, which can guide you during your work, while also being piggybacked through the entire session.


From the **Main Matrix** screen:

Click **Program** (top menu bar)

Click **Short Sarcode Feedback Rx** (drop-down list)

The screenshot shows the PLUS GART™ Universal Biofeedback System software. The interface includes a menu bar at the top with options like 'Programs', 'Additional', 'Subspace Programs', 'Reporting', 'Auto Focus Zep', 'Reset', 'Order of Remedies', 'SPR Enhancements', 'Library', 'System Power Settings', 'Audio Player', and 'Help'. The main window is divided into three main sections. The left section is a 'Programs' list with various biofeedback protocols. The central section displays a detailed report for 'John Doe', including a 'VALSOPRE' report, 'Info Report', and 'RETEST' options. The right section shows a 'Data for John Doe' summary with fields for 'Owner', 'User', 'SYS POW', 'Percent Done', 'Time remaining', 'Reactivity / Rectify', and 'Shaping (vsp: 85)'. At the bottom, there are buttons for 'Treat', 'Hold', 'Search', and 'Close'.

The following panel shall appear:

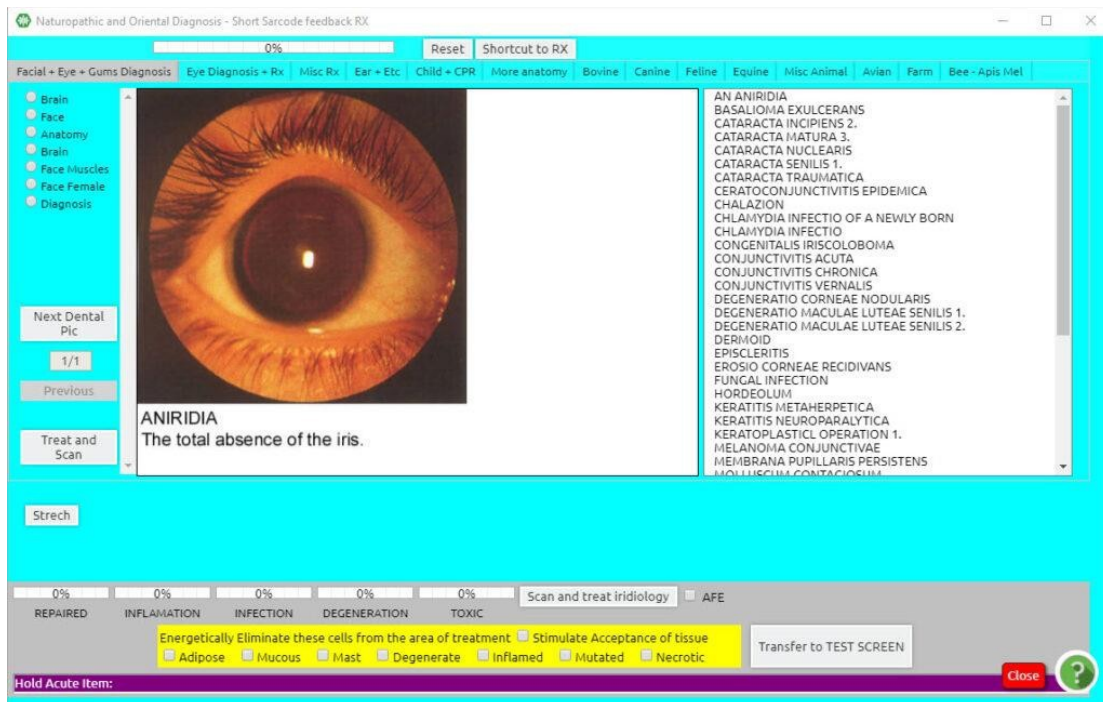


Click *Shortcut to Rx* (upper middle button)

Click *Load* (in the middle of the panel) and take note of all the items.

Click *Close* (in the middle of the panel)

By loading this information into your report, you have *Piggybacked* the information, so that it can be addressed throughout your entire session every time the *Pulsar working* panel is shown.

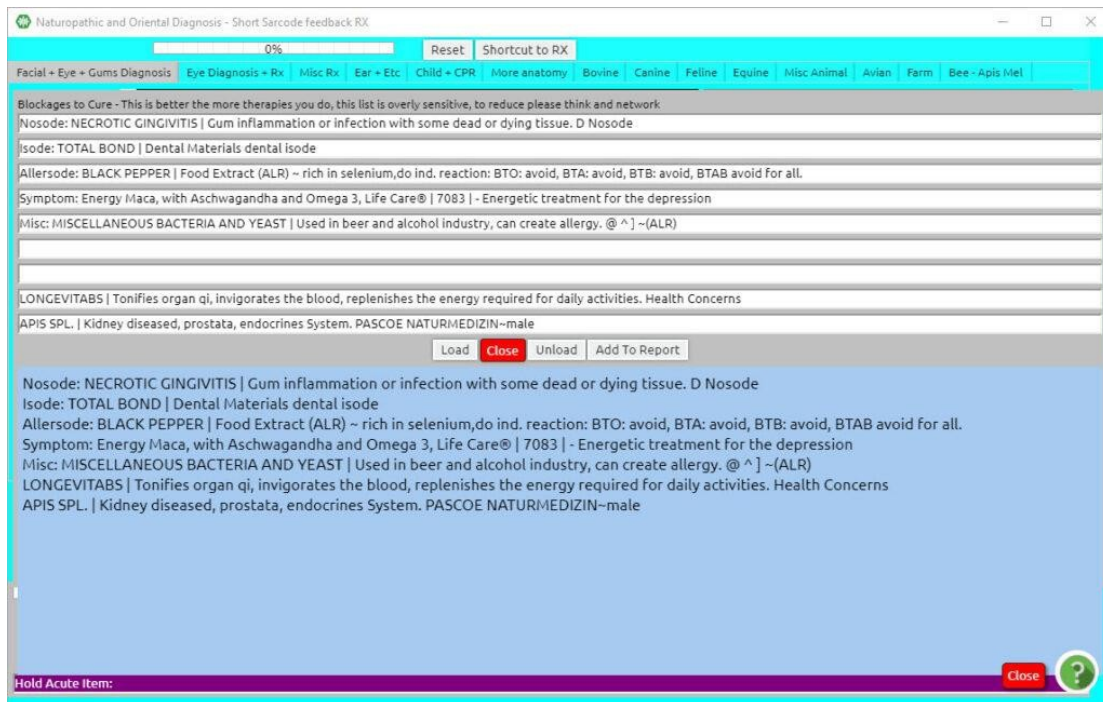

You will also notice a number of tabs on the top part of the Short Sarcode main panel:

- Facial + Eye + Gums Diagnosis
- Eye Diagnosis + Rx
- Misc Rx
- Ear + Etc
- Child + CPR
- More Anatomy
- Bovine
- Canine
- Feline
- Equine
- Misc. Animal
- Avian
- Farm
- Bee – Apis Mel


From the middle part to the right, the tabs refer to animals. This is where you can work animal specific. Let us go into each tab so you see how you can work and apply them to your session.

Note:

Depending on the size of the display of your laptop, you may need to scroll down to see that, at the bottom of each tab, you have an indication for Repaired, Inflammation, Infection, Degeneration, and ***Toxic*** in a percentage bar. When you do a treatment, your objective is to reach at least **85%** on the ***Repaired*** bar, while the other bars (Inflammation, Infection, Degeneration, and Toxic) should be at least 10% lower than the Repaired percentage. For some clients, especially first-time clients, it may take some time to achieve this.

You also have the option to add the following add-ons to the therapy:

- Adipose
- Mucous
- Mast
- Degeneration
- Inflamed
- Infected
- Mutated
- Necrotic
- Click ***Ear + Etc*** (top menu bar).

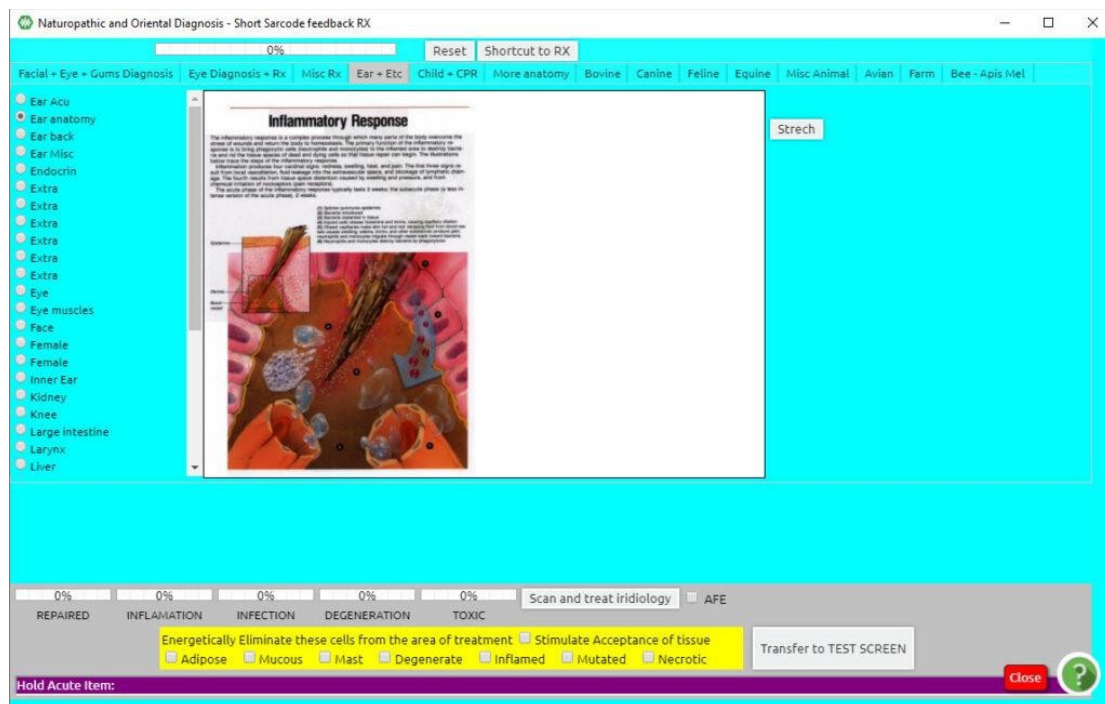


The initial default panel will show ***an Inflammatory Response***. It is always a good idea, when entering this panel, to lower the inflammatory response of your client.

If you select one of the items on the left, you won't be able to return to the Inflammatory Response, which is set as default when entering this panel.

DC/ the grey background of the ***Inflammatory Response*** screen to initiate a training protocol for any stress associated with inflammation. Once initiated, the screen will be scanned, showing that it has been activated. The value indicators at the bottom of the panel provide information based on reactivity.

Once you rectify the Inflammatory Response, you can choose any of the items from the left and proceed with therapies.

For all the training programs in ***Short Sarcode Feedback Rx***, a good rectification value for the Repaired category is 85% or more, or the Repaired category must be at least 10% higher than the other categories (Inflammation, Infection, Degeneration, and Toxic).

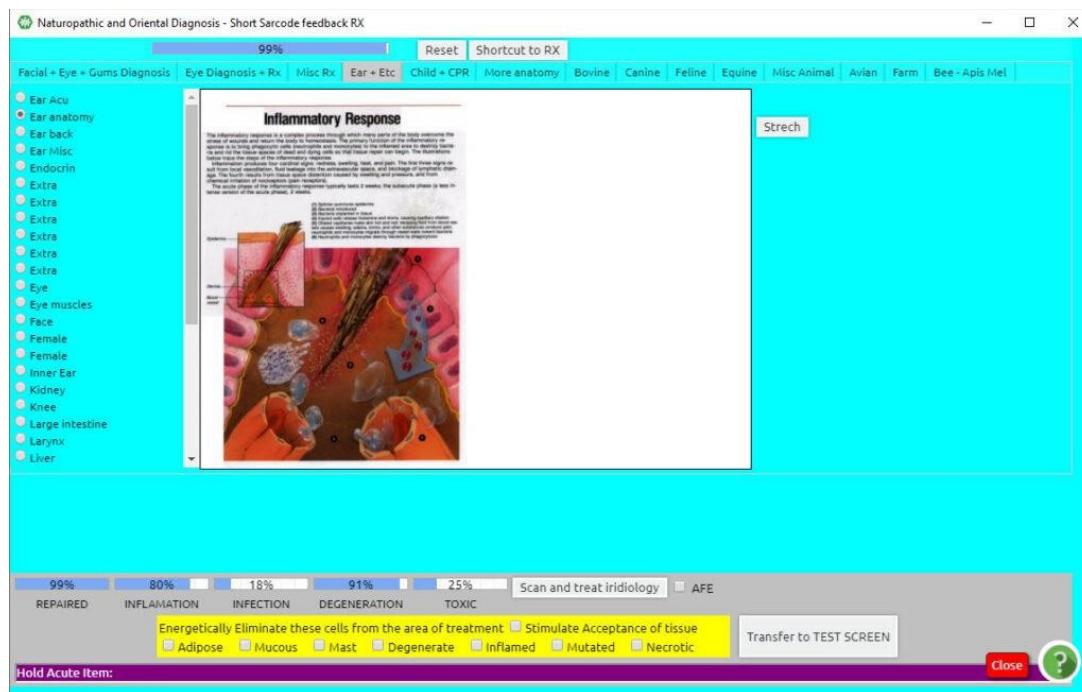

Click ***the Lymph*** button on the left side of the panel.

Double Click on the white background of the screen to initiate a training protocol for all parts of the lymphatic system shown on the panel. The value indicators at the bottom of the panel provide information based on reactivity. Click ***Lungs*** on the left side of the panel with the small radio button.

Double Click on the white background of the screen to initiate a Training Protocol for all parts of the lung shown on the panel. The value indicators at the bottom of the panel provide information based on reactivity.

You can continue working in this way to balance stress associated with all the major organs of elimination.

Consider: Lymph, Lungs Kidney, Small Intestine (including Liver), Male or Female Organs, and **Endocrine** as these will cover the major organs of elimination.

Note:

All screens within the ***Short Sarcode Feedback Rx*** panel work in the following way.


Each Training Protocol is engaged by simply double-clicking anywhere on the white background of the chosen panel. This provides a training balance for all of the areas shown in the panel.

Clicking on the label (or directly on a specific part of the image) initiates additional specialized training for that area.


This is the first tab. It mainly concerns issues of the eyes and dental, with a few other options as well.


The pictures showed depending on what is selected on the list on the right. By default, it is on the first selection of Aniridia.

If you have an eye diagnosis, you may select the issue directly (i.e. cataract).

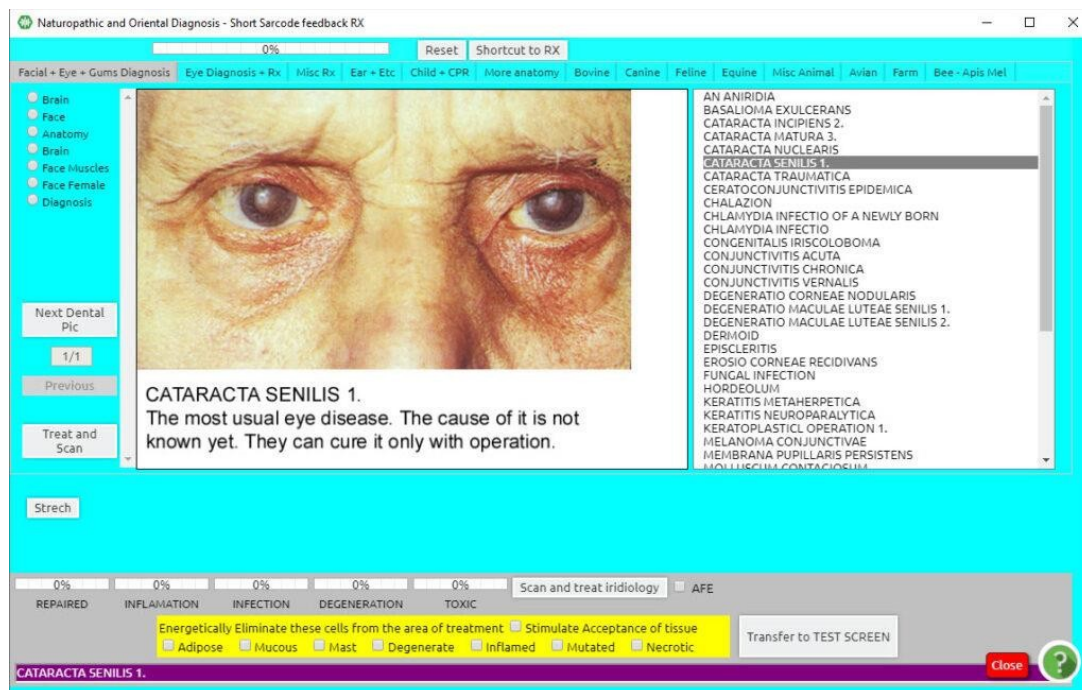

Double Click in the white area to treat, or Click on the button on the mid-left side “Scan and treat”.

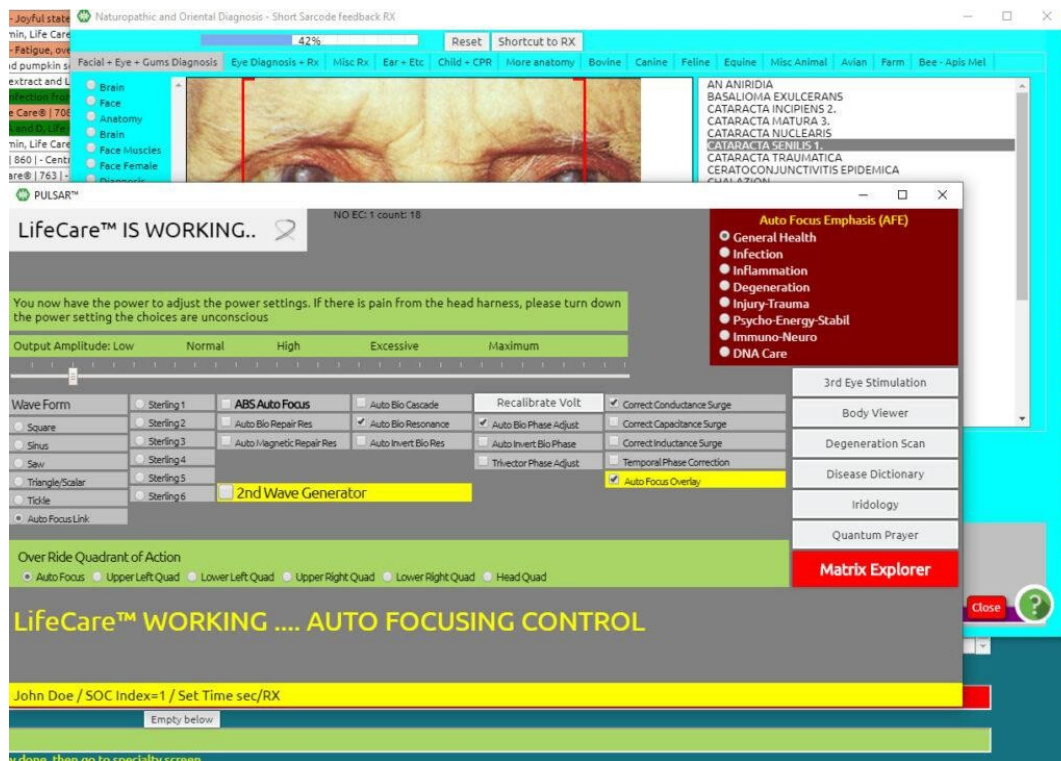


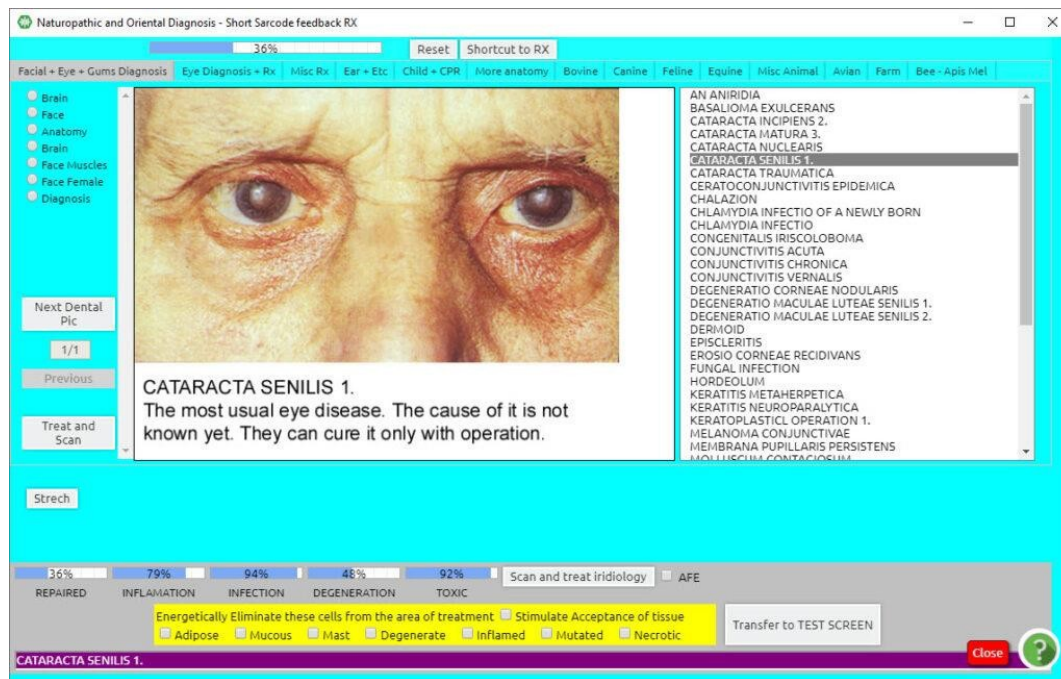

After you *Treat and Scan*, at the bottom of the page you will obtain information concerning the percentage of repaired, inflammation, infection, degeneration, and toxic.

Your objective is always to have at least 85% on the repaired bar, while at the same time all other bars (inflammation, infection, degeneration, and toxic) should be at least 10% lower than the repaired percentage. Repeat as many times as you reach your objective, or depending on the available time.


Furthermore, on this tab, you can also work on gum and Dental issues. By selecting “Next Dental Pic” on the mid-left part of this panel, you will notice pictures of gum/dental issues.

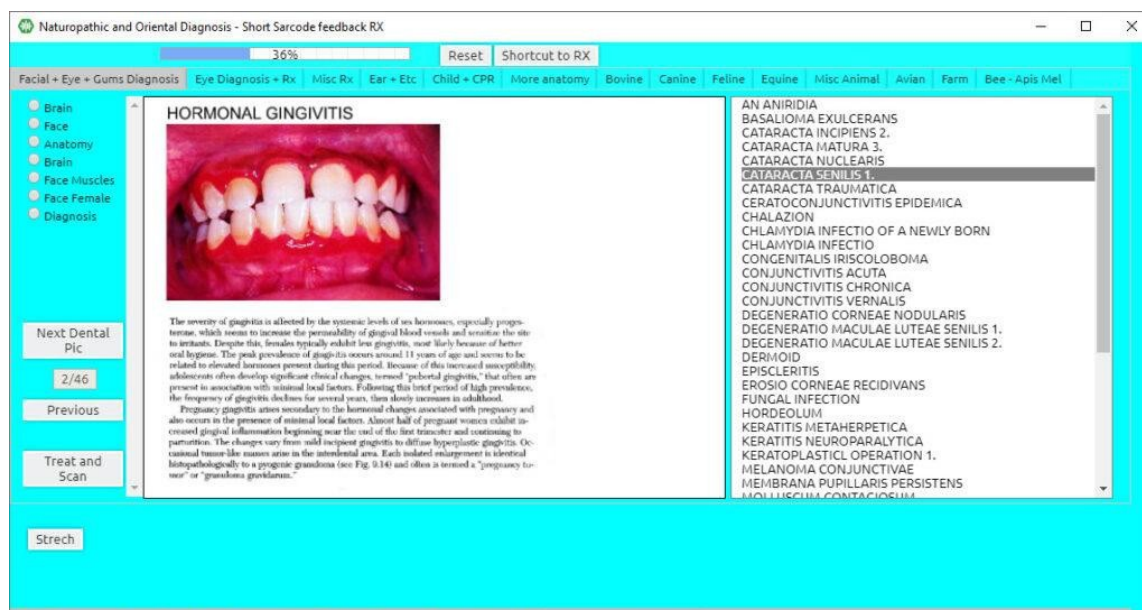

There are 46 pictures, you can go to next or previous till you reach the item you wish to work on.

Other items to work on this panel, on the top left you may select:


- Brain

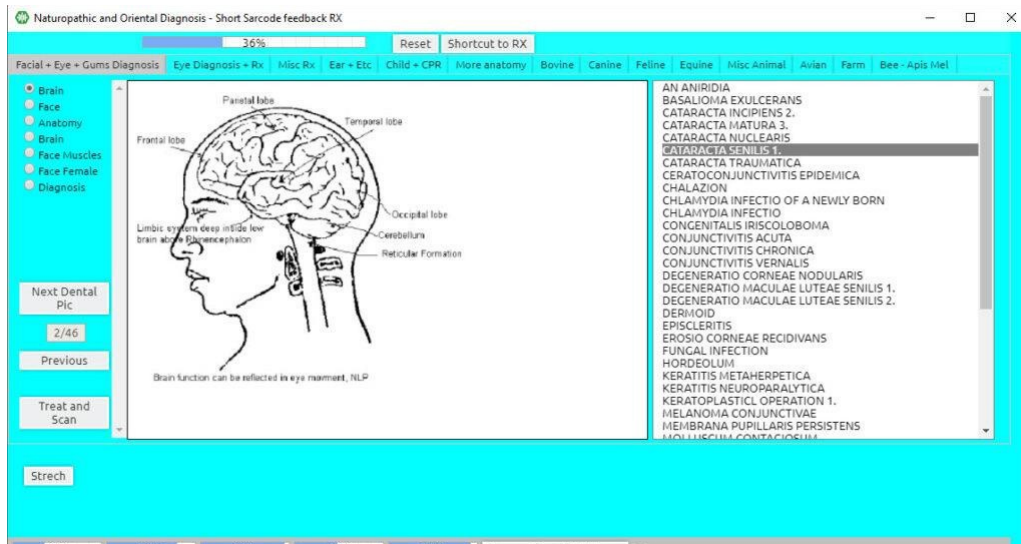

You can either work on a specific part of the brain by double-clicking on the name of the part of the brain, or you may double-click on the white area to work on the entire brain.


Again, please notice the rectification values on the bottom and repeat as needed.

- Face

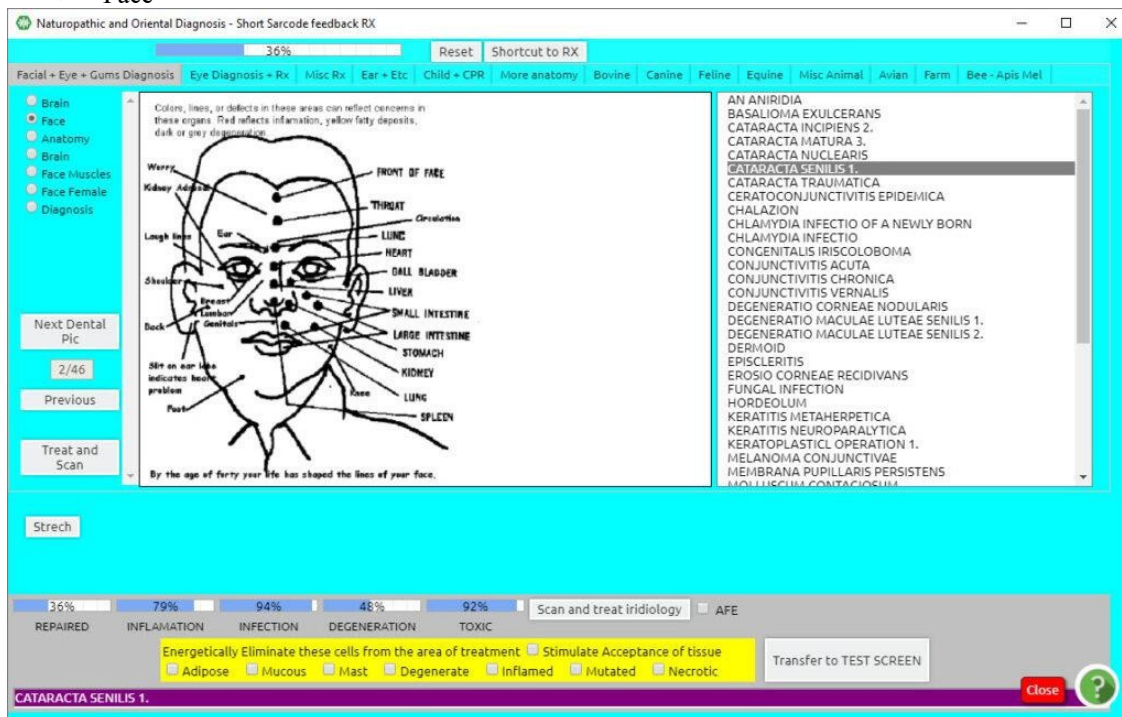


Here you can work on the areas of concern through the face. Our faces show much information through lines, aging, spits, and discolorations, issues of health concern. You can work on specific areas by clicking on the word, or you may double-click in the white area, whereby you can work on the face and health in general. Again, please notice the rectification values on the bottom and repeat as needed.

• Anatomy

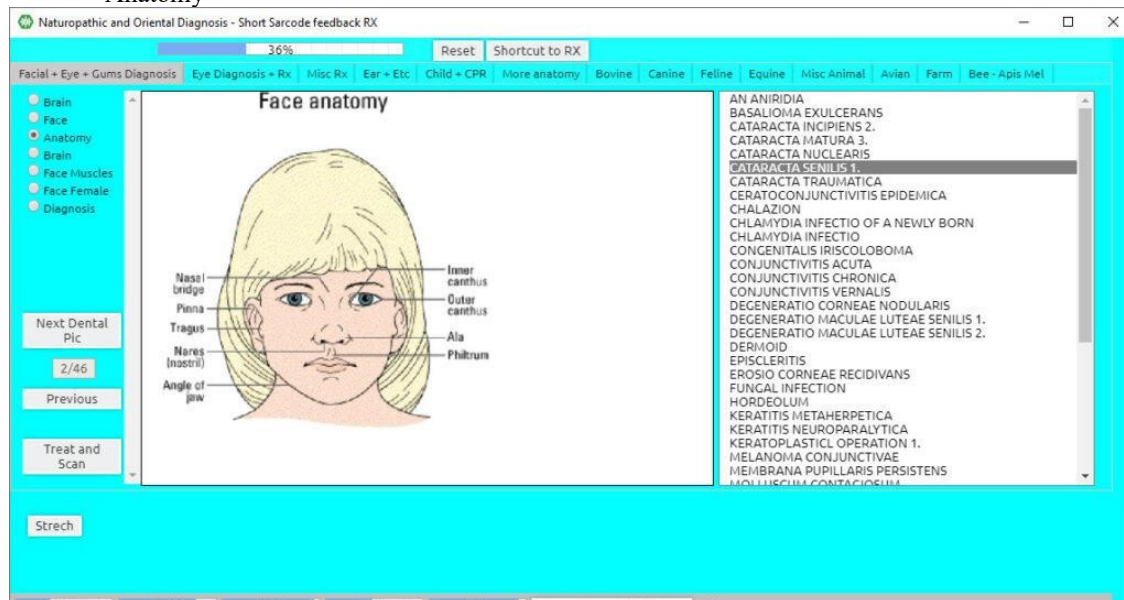

Here you can work on the areas concerning the anatomy of the face. You can work on specific areas by clicking on the word, or you may double-click in the white area, whereby you can work on the facial anatomy in general. Again, please notice the rectification values on the bottom and repeat as needed.


• Brain

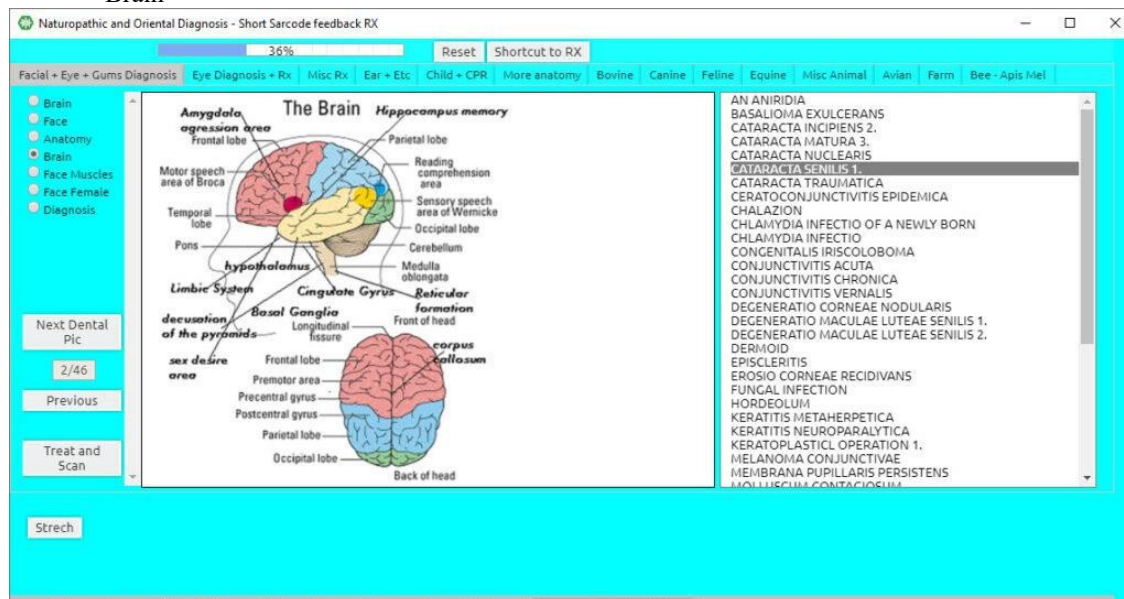

Here is yet another option to work more specifically on the areas concerning the brain. You can work on specific areas by clicking on the word, or you may double-click in the white area, whereby you can work on the brain in general.

Again, please notice the rectification values on the bottom and repeat as needed.


- Face Muscles

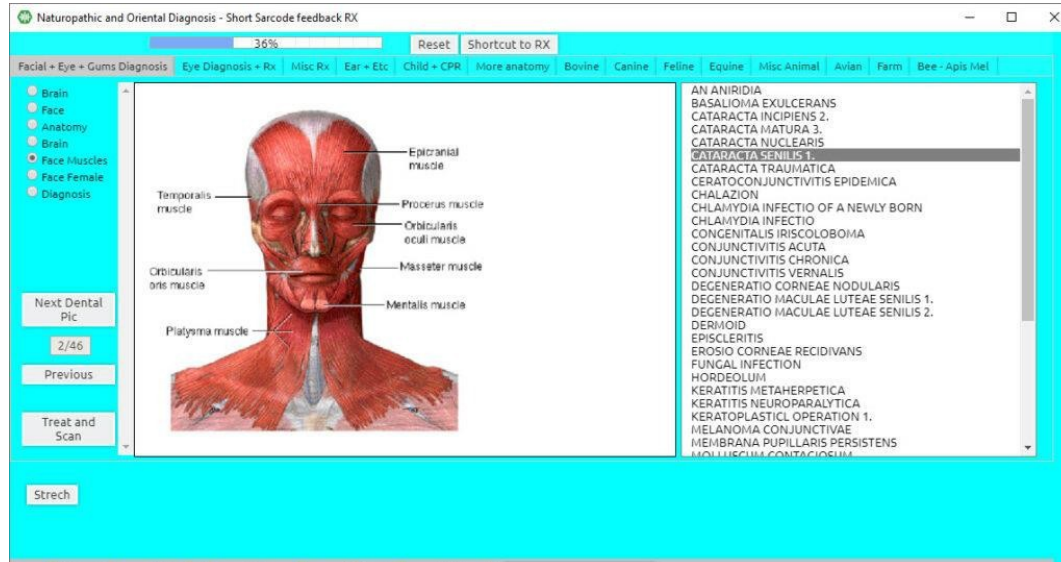

Here, you can work on the facial muscles, smoothing out lines and tension on the face, cranium, or TMJ area. You can work on specific areas by clicking on the word, or you may double-click in the white area, whereby you can work on face muscles in general.

Please check the rectification values on the bottom and repeat as needed.

- Face Female

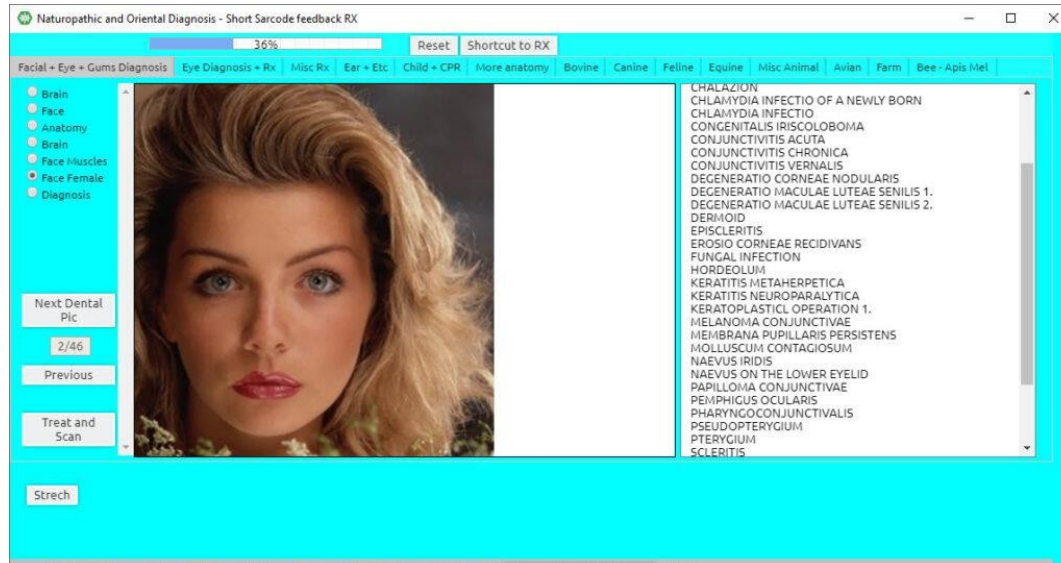

Here you can work on the face in general, areas of concern for anti-aging, spots, skin, or any other issues, smoothing out the issues on the face. You can work on specific areas by clicking on the word, or double click in the white area, whereby you can work on the face in general.


Remember to check rectification values on the bottom and repeat as needed.

Diagnosis

Here you can work on any part of the brain in general or areas of concern. You can work on specific areas by clicking on the word, or double click in the white area, whereby you can work on the brain in general.

Remember to check rectification values on the bottom and repeat as needed.


Click **Eye Diagnosis + Rx** (top menu bar)

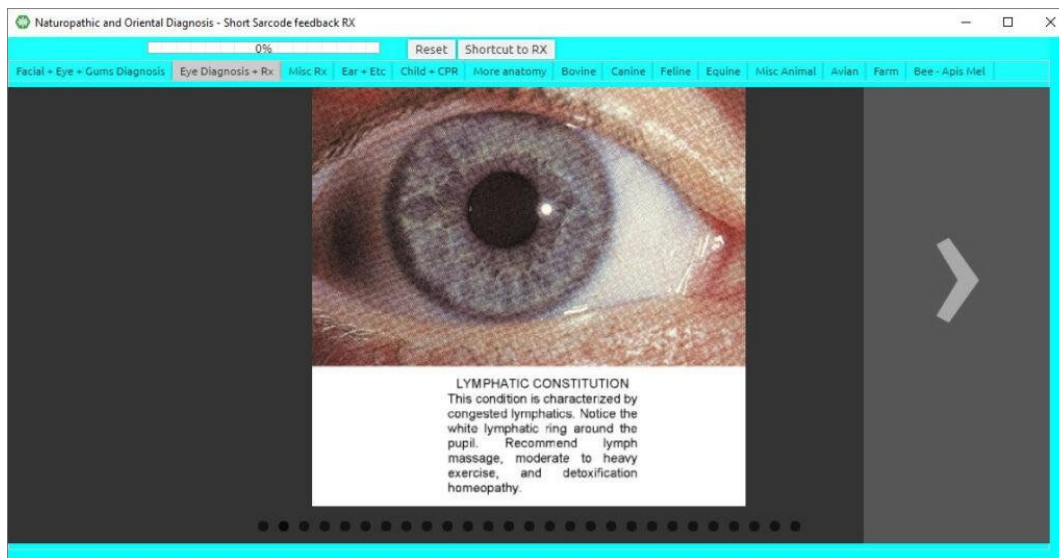

Press the second tab on the top menu of the Short Sarcode panel, and you will enter Eye Diagnosis + Rx.

Within this panel, you may click on the left or right side of the picture where there is an arrow indicating to either move the panels forward or back. Within this panel, you can work on the area of concern depending on the issue.


Remember to scroll down and check rectification values on the bottom and repeat as needed. Click **Misc Rx** (top menu bar)

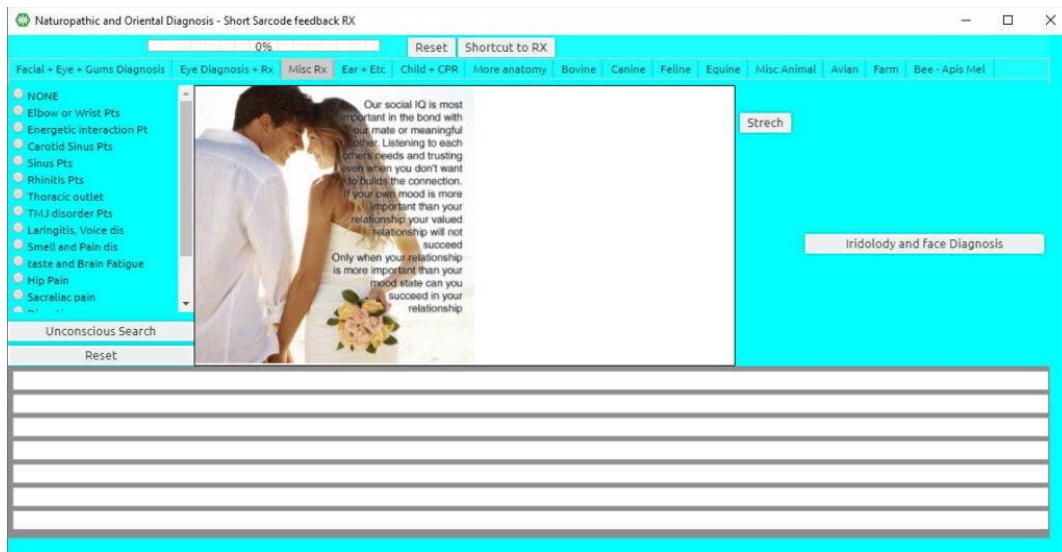

The third tab on the top menu of the Short Sarcode panel is Misc. Rx. You can use this panel in several ways. By clicking **Unconscious Search** several times, you obtain information that has been taken from the testing process in the beginning. This gives you hints on what needs to be worked on during the session, while also loading the items to piggyback for the entire session. Furthermore, by loading these items, we create an area of focus while doing our treatments in Short Sarcodes, and we are also layering in information as we progress in our sessions.




On the left side of the panel, you can see various areas on which you can focus.

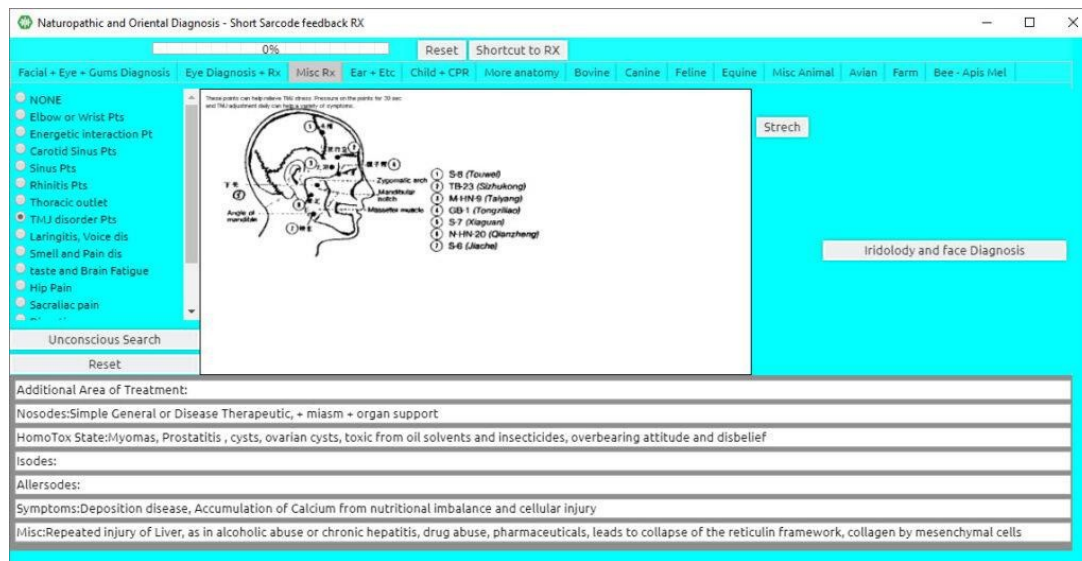

Choose the item you wish to work on, depending on the area of concern (for example TMJ). You can either treat specific areas till rectified or the entire area in general.

Remember to check rectification values on the bottom, and repeat as needed.


On the right side of this panel, you shall notice a button “Iridology and face Diagnosis”. By clicking it once the following appears:

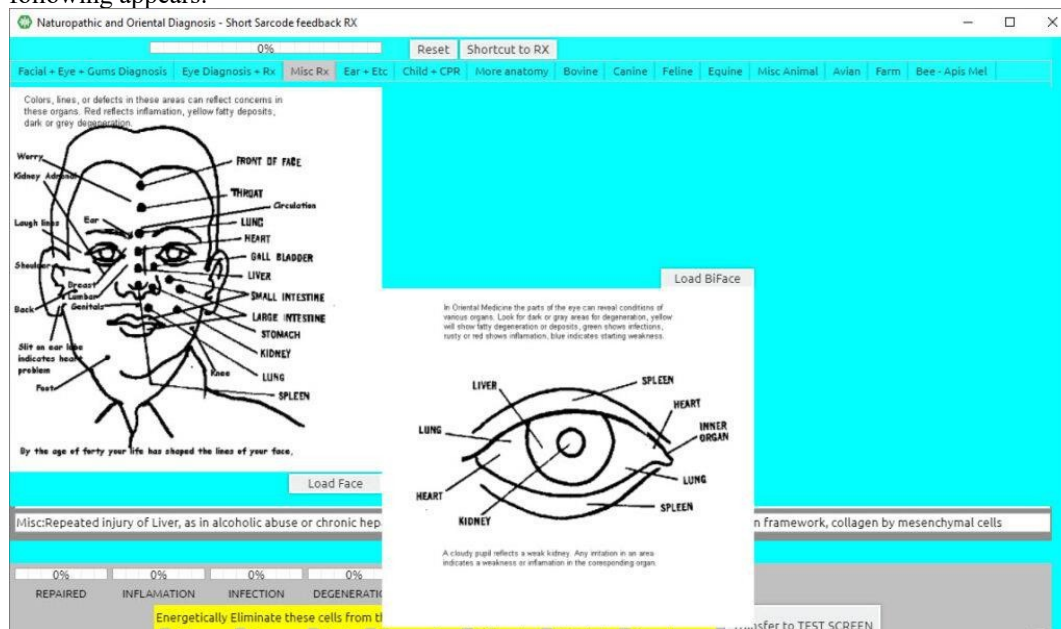

Here, you have the option to work on your client’s health through Facial points, while also being able to load a picture of the face and the Biface. We will get into that in another chapter and video.

Click **Child + CPR** (top menu bar)

Click the fifth tab in the top menu Child + CPR.

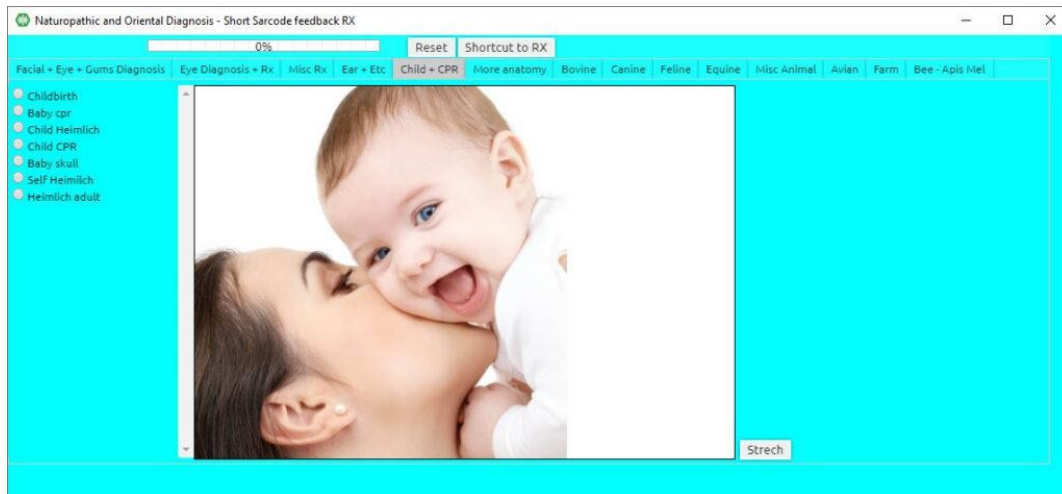


This panel is child-specific, for childbirth, in addition to useful information containing CPR for children while also adults. You can enter the various options on the left side of this panel.


HINT: the various options in this panel can also be used as the treatment of birth trauma, past trauma from choking or CPR, or preparation for birth. However, care must always be taken when treating a pregnant woman. People treating pregnant women must have certification and experience.

Again, please scroll down and notice your rectification values on the bottom and repeat as needed. Click **More Anatomy** (top menu bar)

Click the sixth tab on the top menu bar of the Short Sarcodes, More Anatomy.

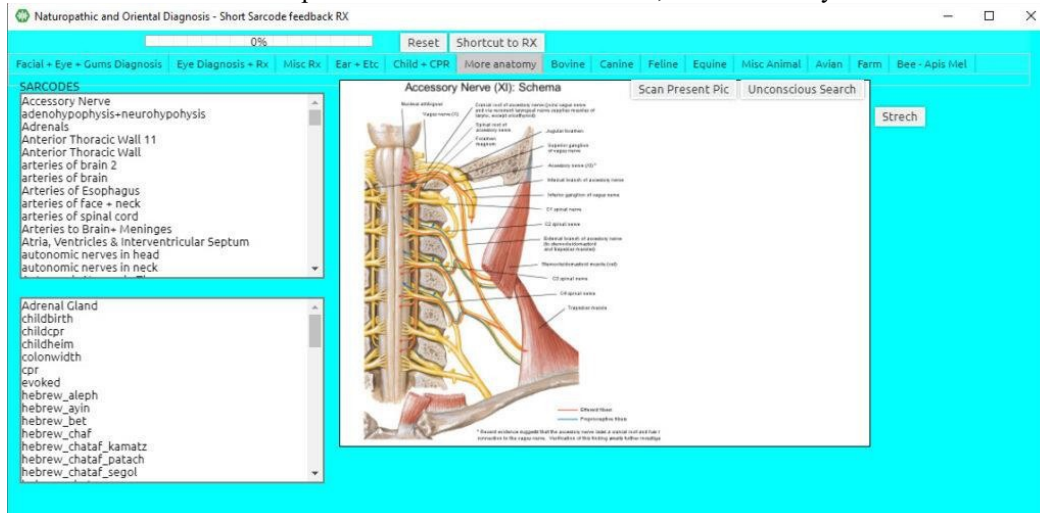

This panel contains full anatomy and is very useful when you want to be specific on your therapy, or you can also do 'Unconscious Search', located in the top right section under the top menu bar.

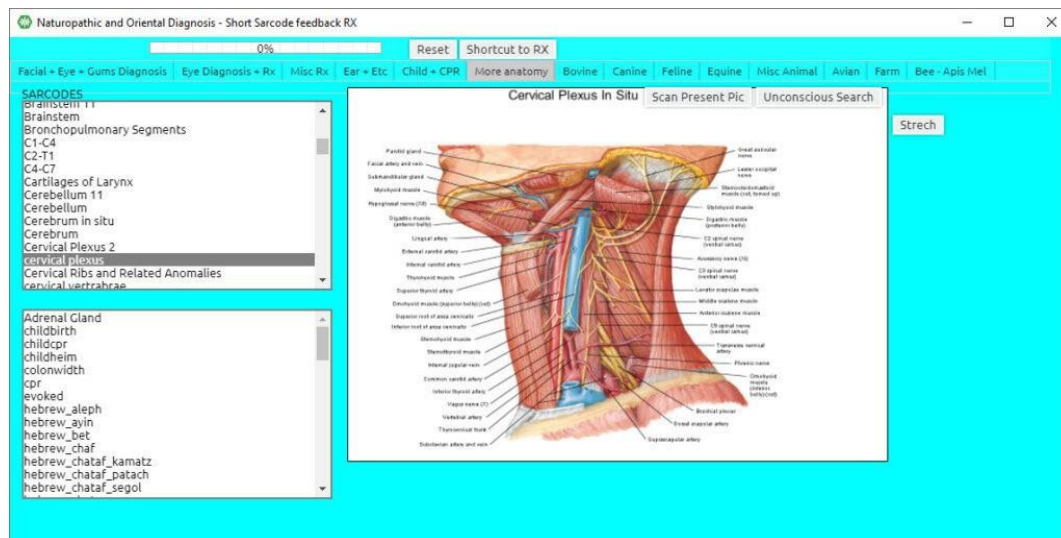


Treatment can be done either specific to the area of concern by clicking the specific word, or in general by clicking the white area.


Remember to check rectification values on the bottom while doing your therapies and repeat as needed.

The remaining tabs are animal-specific and work in the same manner. However, you shall notice that there are two scroll bars on the left side. The top scroll option is organ related whereas the bottom scroll bar is issue related.


Click **Bovine** (top menu bar)

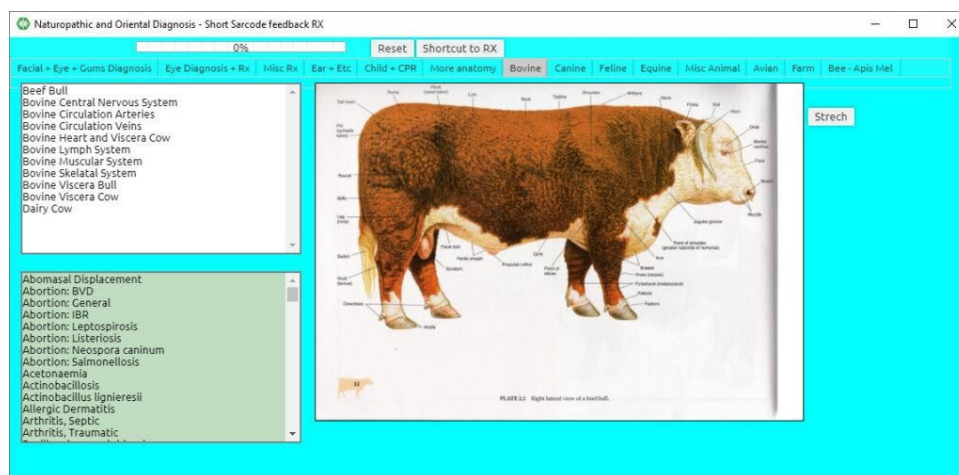


Click **Canine** (top menu bar)

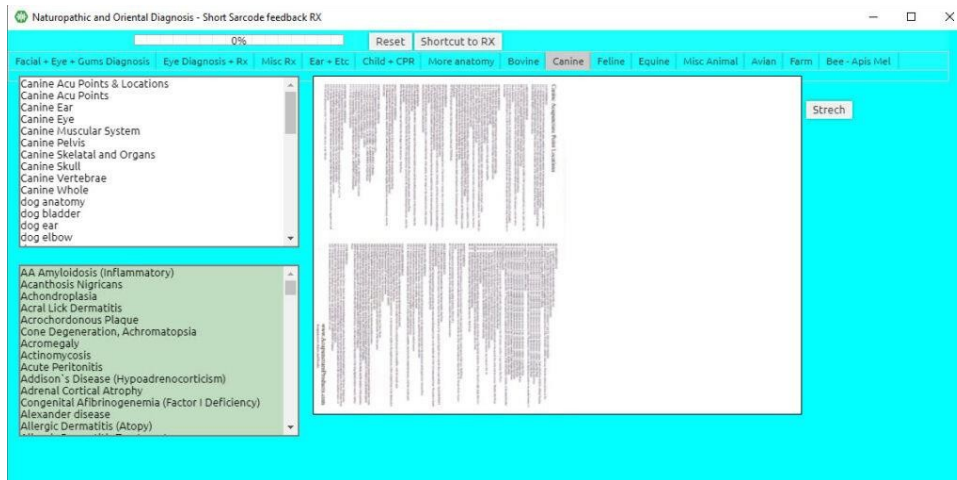


Click **Feline** (top menu bar)

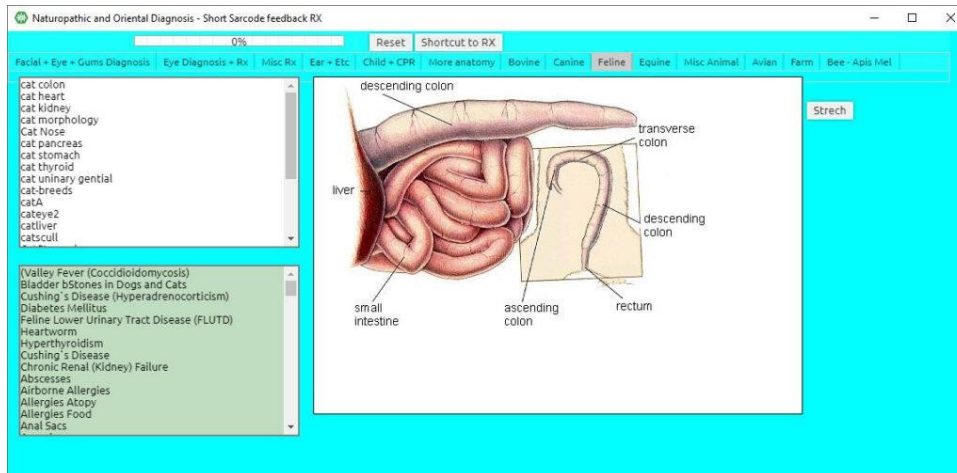


Click **Equine** (top menu bar)

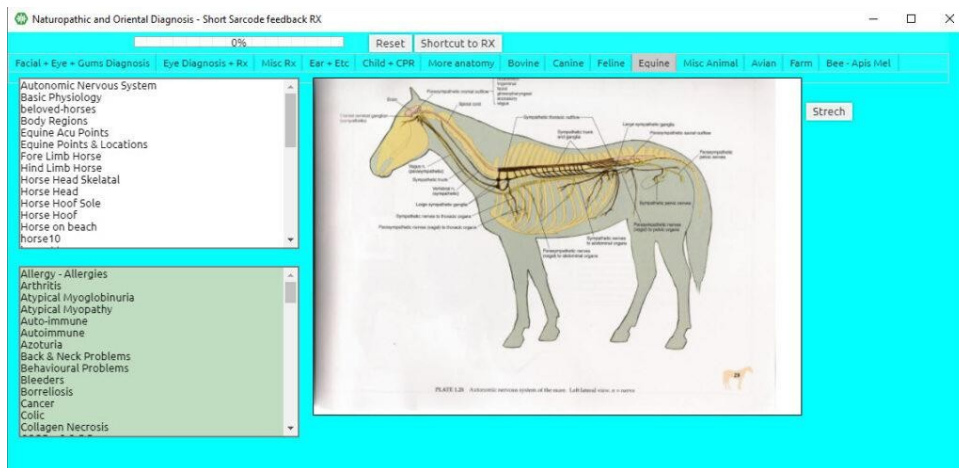


Click **Misc. Animal** (top menu bar)

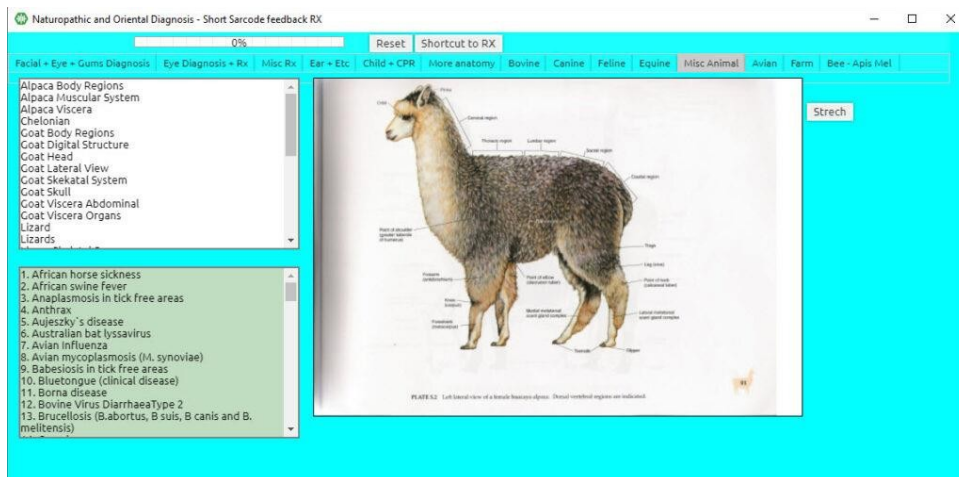


Click **Avian** (top menu bar)

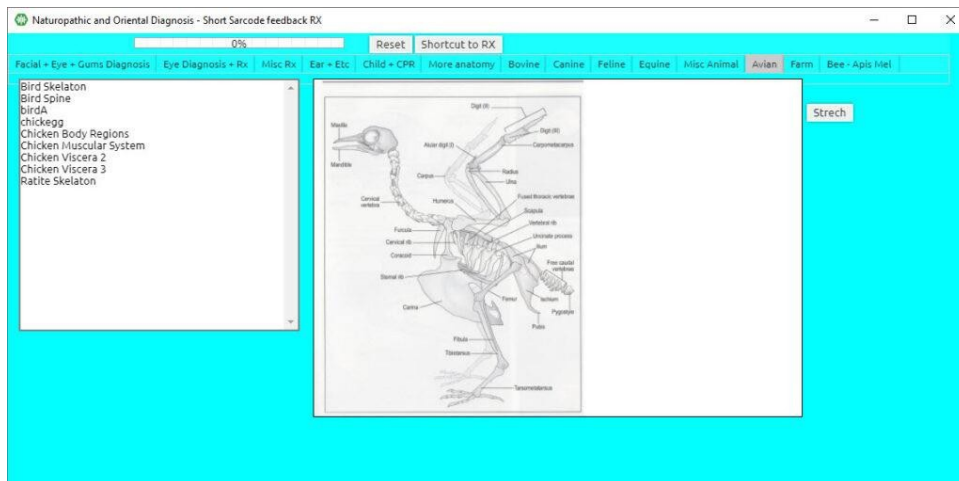

Click **Farm** (top menu bar)

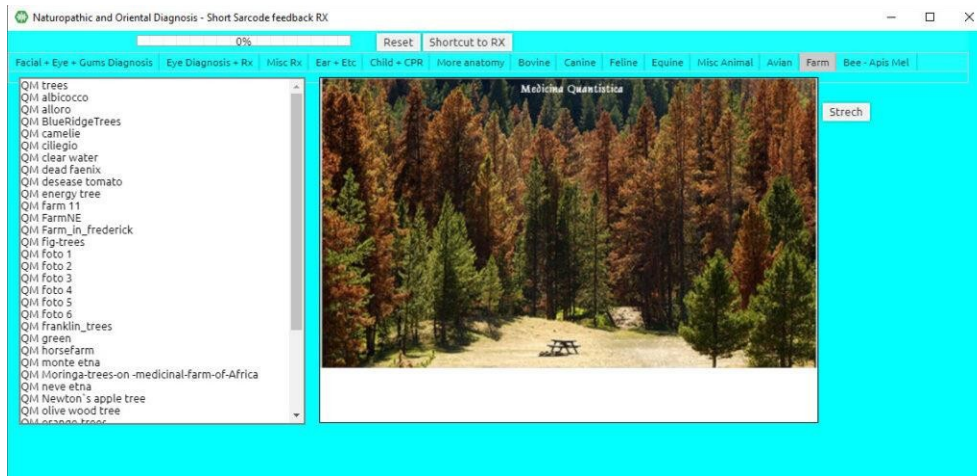

Click **Bee – Apis Mel** (top menu bar)

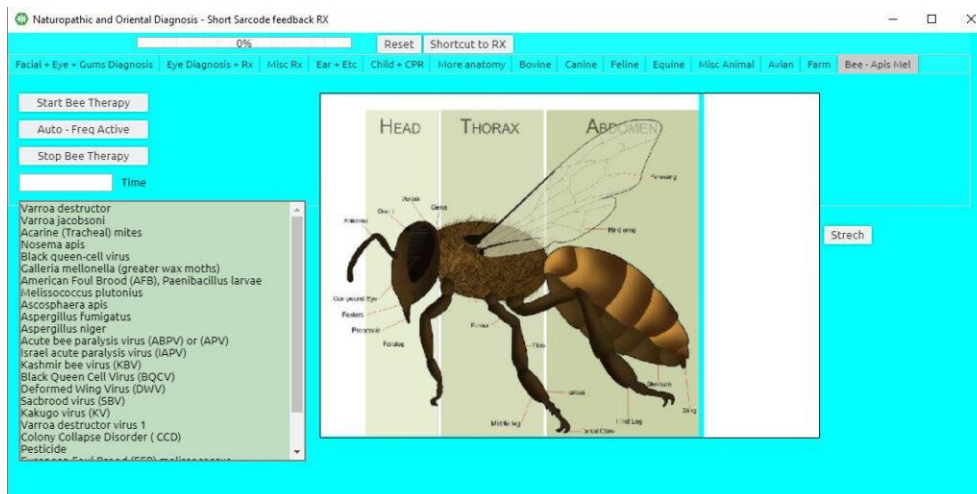

Chapter 4

Bodyviewer Training Program Activation

The **Body Viewer** Program is only available when your device is connected to the computer.

HINT: It is always best to open the Body Viewer and Iridology panels before opening the Disease Dictionary program to prevent the programs from freezing.

From the **Main Matrix** screen:

Click **Subspace Programs** (top menu bar)

Click **Body Viewer** (from the drop-down list)

PULSAR™ Universal Biofeedback System

NEW Programs Additional Subspace Programs Reporting Auto Focus Zap Retest Order of Remedies EPR Enhancements Library System Power Settings Audio Player Help

Body Viewer VARHOPE Info Report [Navigation Icons] Expand

8603	100	Complex of mine	Disease Dictionary (library)	gles, grief, poor vitality, anxiety.
5370	100	Life Impulse® EC	Idiology (library)	
5073	100	Life Impulse® Nu	Personal Evolution	Parkinsonism.
6393	101	Organic Matcha,	Orgone Generator	edema) and high blood pressure. ALO.]
10024	101	Life Impulse® Cc	Quantum Prayer	
4553	101	Life Impulse® Et	Superconscious Choice of Therapies	
5377	101	TurnerCool, wit	3rd Eye Stimulation	ico
4748	100	Life Impulse® Et		
6915	100	ProstaNEW, with pine phytoestrogens and pumpkin seeds, Life Care® 7048 - Generalized edema (dropsy), fluid metabolism problems, diuretic (kidney), gentle digestive strengthener.]		
3644	100	NeuroBooster Complex, with Brahmi extract and L-Tryptophan, Life Care® 7046 - Hysteria with chlorosis, chorea, spinal irritability, dysmenorrhea, bladder tenesmus, hysterical epilepsy, intense.		
3904	100			
4006	100	Reishi Mushroom Extract, 400 mg, Life Care® 7081 - Imbalance in patterns of waking and sleeping-such as insomnia, jet lag, time displacement		
3857	100			
6163	100	TurnerCool, with piperine and curcumin, Life Care® 7012 - Nutrient rich algae. VITAMIN- heavy metal detoxification, energy		
10880	100	Life Impulse® Vitamin C and Calcium 860 - Central Nervous System, Chinese Herbal Formula, anti-carcinogenic formula designed to aid in cellular integrity, immune support and detox		
5625	100	NoAlgero, with black cumin oil, Life Care® 763 -) - listed in d. dictionary program, sugar regulation.		
7370	100	Resveratrol Plus, Life Care® 7490 - Garlic concentrate in a base of vitamin C and chlorophyll to lower cholesterol + Immune. BIOTICS E.N.B.		
7798	100	Life Impulse® Lady's Desire 740 - Chronic fatigue, stomach tonic, NETSMANN-Arzneimittel		
9090	100	OPC Antioxidant, with grape seed extract, Life Care® 7038 - Diabetes mellitus - Sri Ramana Maharshi Natural Spagyric Remedies made for Duplication		
5742	100	Alpha lipoic acid complex, with biotin, Life Care® 7051 -) - listed in d. dictionary program		
5289	100	Stabilized oxygen molecules Life Care® 891 - B15 and other B vitamins from special beer yeast extract-oxygenation, I.Q.		
10307	100	Mineral and vitamin complex, for pregnant women, Life Care® 7071 - Intestinal dysbiosis, digestive disorders & yeast infection		
10626	100	Body cream, with Shea butter and BIO plants, Life Care® 4320 - remedy for cell regeneration and activation / cicatrization, burn trauma, fracture, activation of circulation		
4659	100	GeL, with Snail extract, Life Care® 4046 - Shame in sexual pleasure-helps accept physical self, overcome physical loathing even acne, assists sexual pleasure, acceptance of one's physical body, openminded, hgdxr hgdxr		
9685	100	Life Impulse® ECO Aloe Vera Juice with cranberry flavor 841 - supports detoxification of cigarettes/tobacco/chemicals/ smoking / addiction-		
9016	100	Life Impulse® Super Energy 815 - ASC Plus for male fertility support		
10785	100	Life Impulse® Vitamin C and Calcium 860 - Natural product with anthocyanins la strong antioxidant / Cardiac, vascular, blood circulation, immune system, stomach, intestines and eyes		
1462	105	Stabilized oxygen molecules Life Care® 891 - Methyl donor important for oxygenation, -g vitamin, I.Q.		
1997	100	Organic Pure Moringa, Life Care® 7036 - Bruises, accidents, sprains, bruised feelings, helps prevent strokes.		
3213	100	Organic Pure Moringa, Life Care® 7036 - Bruises, accidents, sprains, bruised feelings, helps prevent strokes.		
3248	100	ProBioLife, with germ complex, Life Care® 8042 - Whooping cough, respiratory affections (lung), laryngeal phthisis, tubercular glands.)		
2748	100	Zinc, Life Care® 7074 - Can id acid alkaline imbalance, electrical sensitivity, mineral-ph		
3126	100	Life Impulse® Vitamins for children 870 - Bones, splinters, imperfect assimilation and defective nutrition; diseases of bone, necrosis; suppurative processes, keloid growth.		
2203	100	DiabForm, with momordica and alpha lipoic acid, Life Care® 7330 - Immune deficiency, manic depressive, diabetes (pancreas), deaf, gout, retardation, sex.)- sugar regulation		
1756	100	Muscle Capsules, Devil's Claw, Life Care® 4101 - Inflammation of the bronchial tree interfering with lung air release.		

– Or –

Wait until the **Body Viewer** icon appears on your taskbar.

Click **Program Start** (top menu bar in the Body Viewer)

Click **Load Unconscious Choices** (start INFO exchange)

Body Viewer

Program Start STOP Animation Organ Systems Chakras Acupuncture Therapies Spinal Psychological Music Misc LifeCare™ Therapy Wipe Therapy Sacred Geometry AFE AutoFocusEmphasis

Load Unconscious Choice (Start INFO Exchange)

Load Individual Movie Program (Start INFO Exchange)

Unload Individual program

Therapist + Unconscious

Chakra Problems

Complete Spinal Therapy

Short Treatment of Disease Path

Treat Disease Path of all Old Trauma

Close Extra Panel

Show Affirmations

Show Mental Conflicts

Show Disease States


Clear Cybernetic movie therapy list

Close


POISONING, SENSORY




Click **Load Individual Movie Program** (from the drop-down list). Allow the graphic to complete.
Click **Therapist + Unconscious** (from the drop-down list)

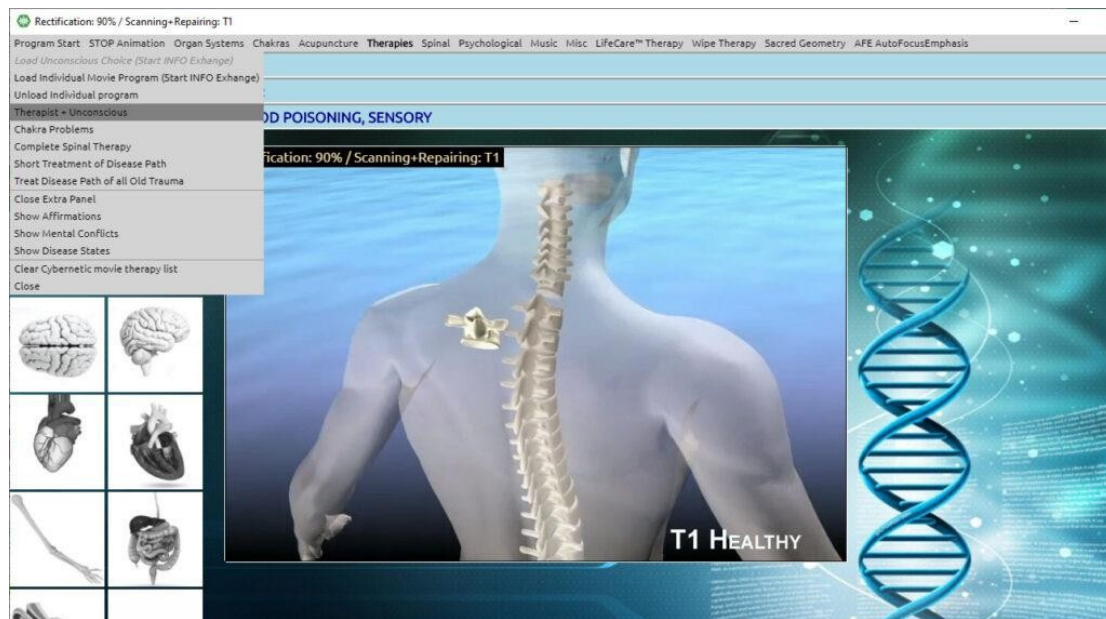

Click **Treat Disease Path of all Old Trauma** (from the drop-down list)

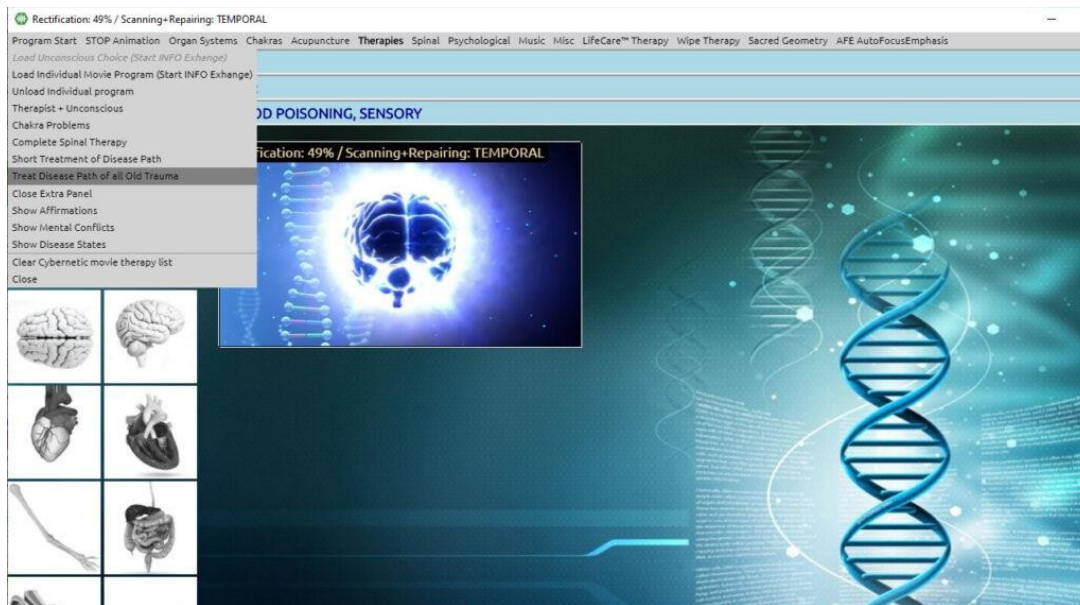


Click the *Consciousness* icon on the bottom task bar to let the program run in the background and go back to the main ED.X and NUCLEUS programs.

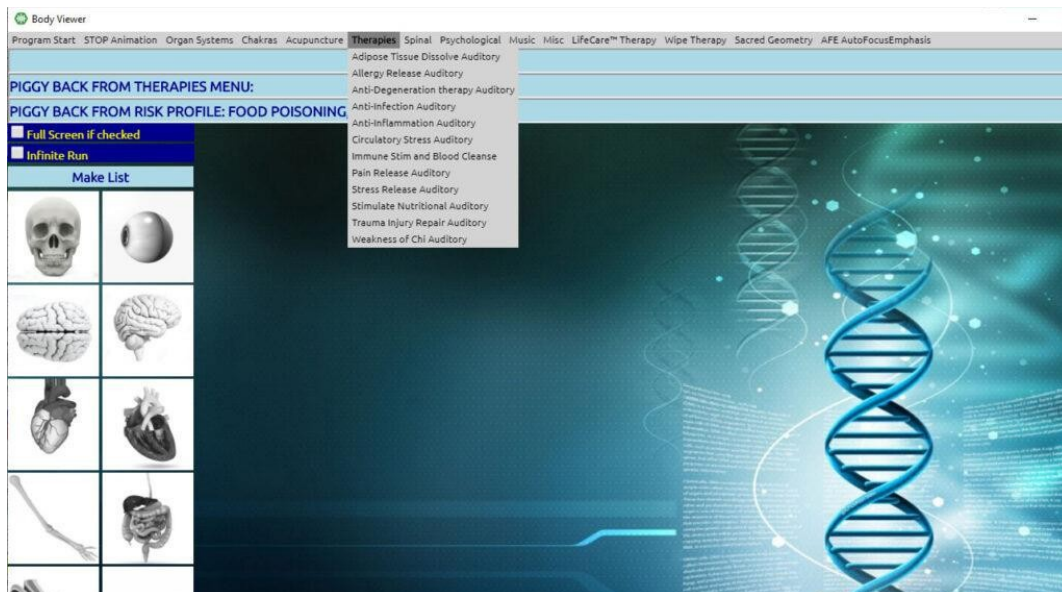

Once you have completed this initial protocol, you are ready to start SELECTING THERAPIES from the Body Viewer.

Organ Systems (top menu bar)

- You have many choices on the *Organ Systems* list, so pick the organs that you consider most relevant for the treatment.
- When you click on any item, it will be loaded into the blue *PIGGY BACK* bars at the top of the screen.
- Consider loading all the major organs of elimination as a beneficial Training Protocol.

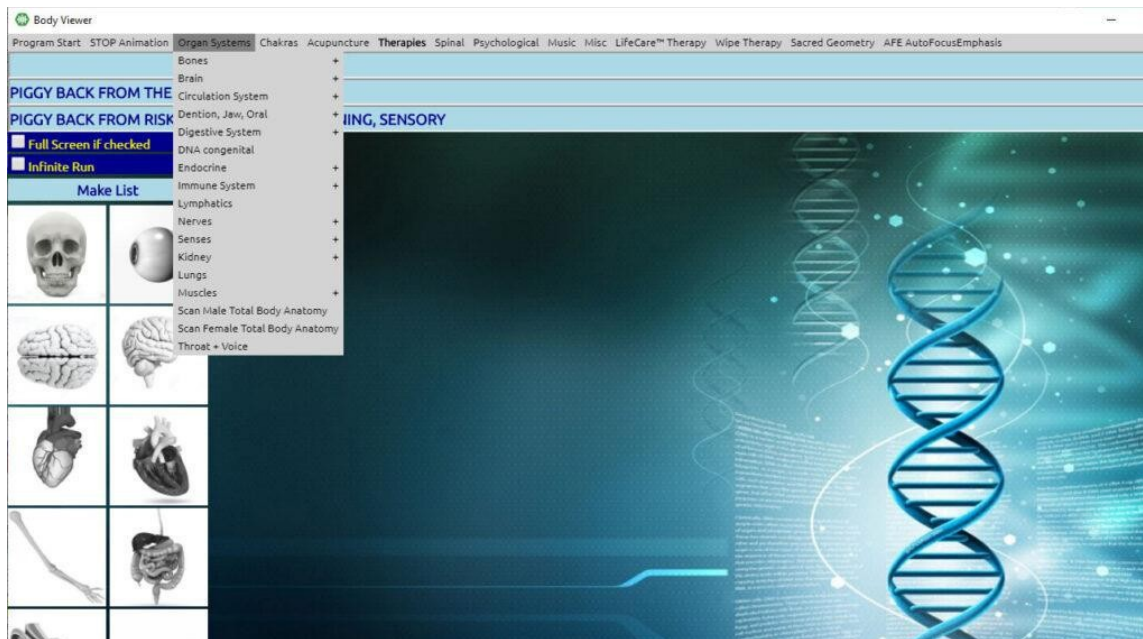

- Click **Chakras** (top menu bar)

Follow the same procedure as above.

Consider selecting **Long Scan Chakra**.

- Click **Acupuncture** (top menu bar)

Follow the same procedure as above.

Consider selecting **Scan Acu Partial System**.

- Click **Therapies** (top menu bar)

Follow the same procedure as above.

- Click **Spinal** (top menu bar)

Follow the same procedure as above.

Consider selecting **Scan Spine**.

- Click **Psychological** (top menu bar)

Follow the same procedure as above.

Consider selecting **Emotions** and **Defence Mechanisms**.

- Click **LifeCare Therapy** (top menu bar)

Follow the same procedure as above.

Consider selecting **Aura Cleanse** or **Foci Detection Automatic**.

- Click **Sacred Geometry** (top menu bar)

Follow the same procedure as above.

Consider selecting **Unconscious Choice**.

- Click **AFE Autofocus Emphasis** (top menu bar)

Follow the same procedure as above.

Consider selecting therapy in conjunction with AFE from the main PULSAR Software to achieve maximum effect.

Iridology

The **Iridology Program** is only available when your device is connected to your computer.

To open the program from the **Main Matrix** screen:

Click **Special Adjunct ED.X and NUCLEUS Programs** (top menu bar)

Click **Iridology** (drop-down list)


– Or –

During any treatment, on the “**PULSAR is Working**” panel:

Click **Iridology** (on the bottom right)


The screenshot displays the PULSAR Universal Biofeedback System interface. The top menu bar includes options like 'NEW', 'Programs', 'Additional', 'Subsequent Programs', 'Reporting', 'Auto Focus Zap', 'Reset', 'Order of Remedies', 'EPR Enhancements', 'Library', 'System Power Settings', 'Audio Player', and 'Help'. Below the menu bar, there's a 'Body Viewer' section with a list of programs. The central panel shows a list of items, each with a number and a description. The right-hand panel features a 'Data for John Doe' section with a table showing 'Owner', 'User', 'Regulatory dysfunction', 'Explore risk of degeneration', 'Reactivity/normal', 'Return to Main Matrix', 'Allersodes', 'Amino Acids', 'Ayurvedic', 'Bach Flowers', 'Bacteria', 'Bio War', and 'Blood'. Below this is a 'HOLD' section with a table showing 'Treat', 'Indiv. Reaction', 'Test item in Hold Tray', and 'Both Trays'. At the bottom right, there's a 'Search for Items' section with a search bar and a 'Clear Search' button.

Iridology program will appear with the conscious button on the bottom taskbar, along with the **Body Viewer**, and will open automatically.

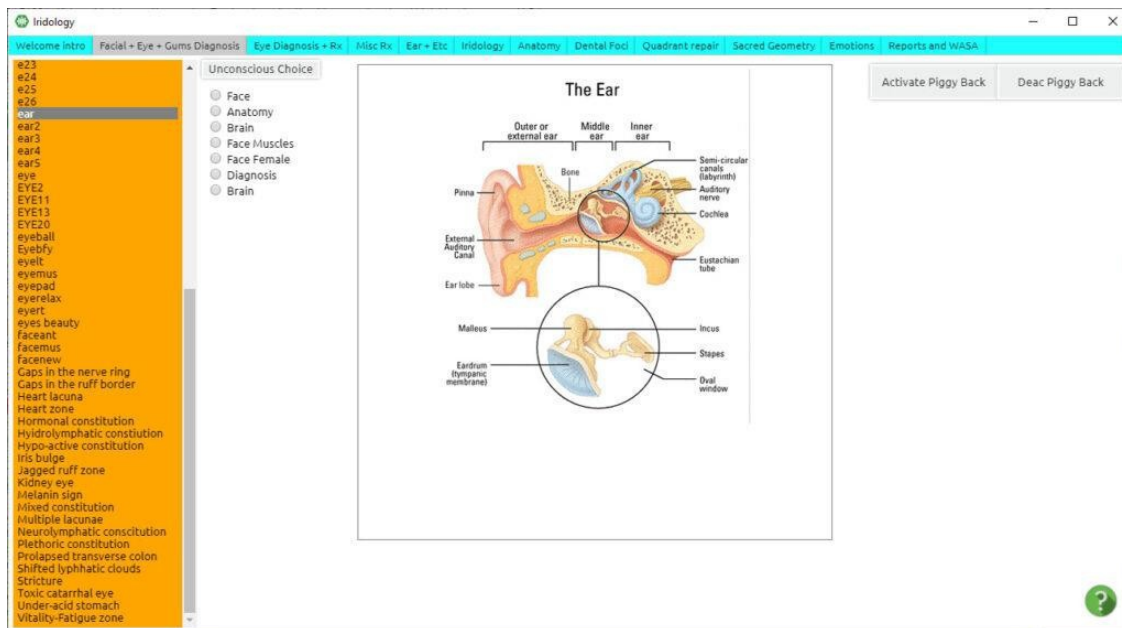


Click **Activate Piggyback** (on the upper right of the screen)

Click **Anatomy** (top menu bar)

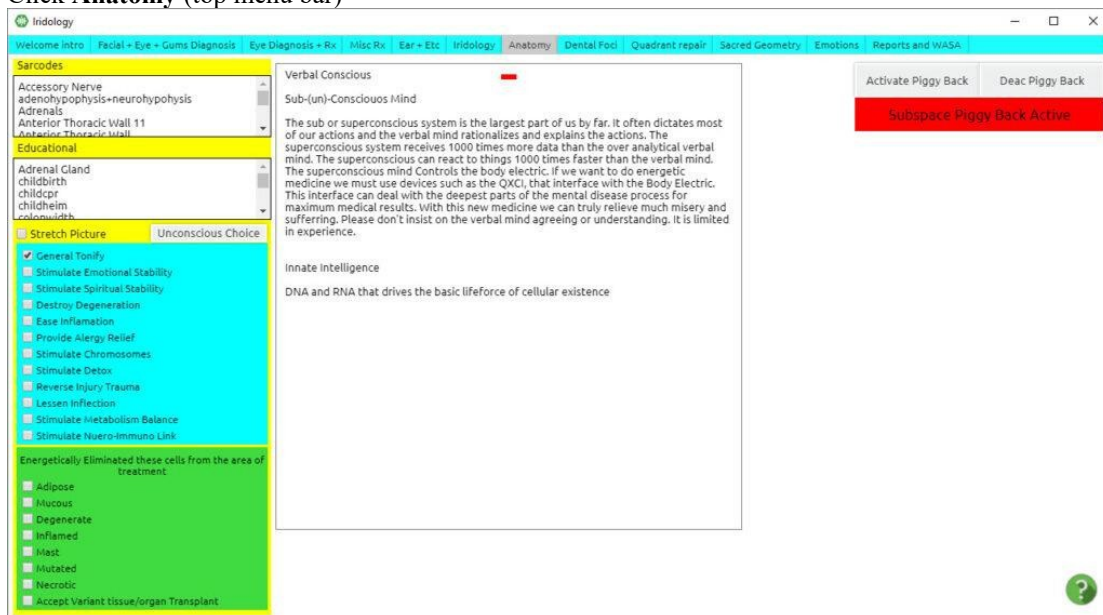

Click **Unconscious Choice** (on the left side, in the middle of the screen)

Click **Activate Piggyback** in the upper right of the screen

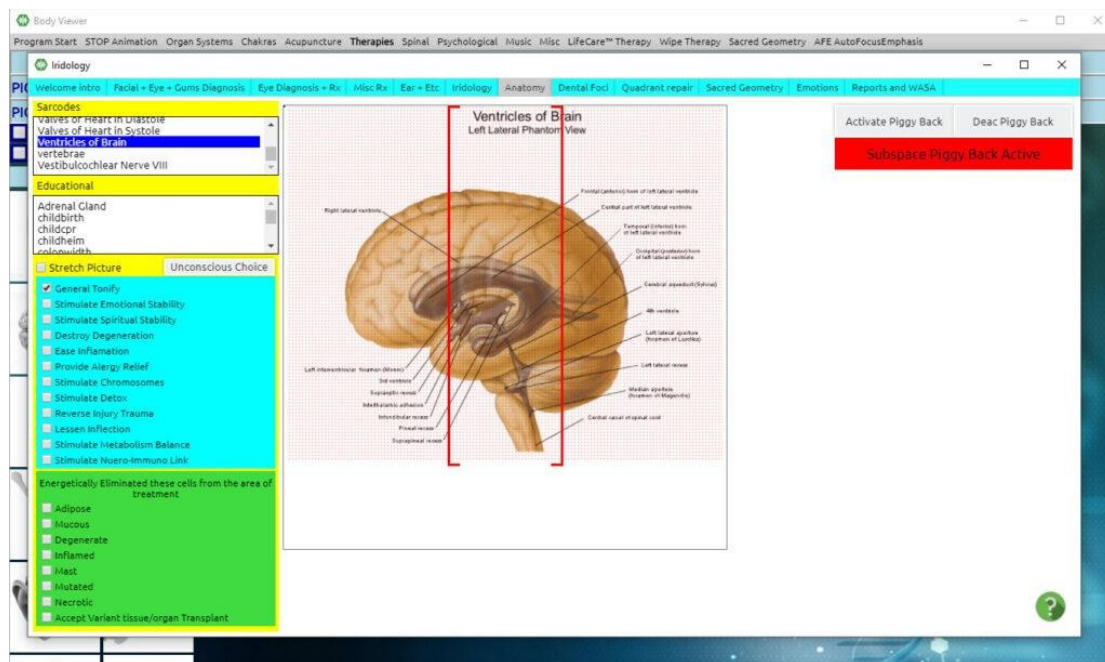

Click **Facial + Eye + Gums Diagnosis** (top menu bar)

Click **Auto Search** (left, middle of the screen)

Click **Activate Piggyback** (upper right of the screen)

Click **Eye Diagnosis + Rx** (top menu bar)

Click **Eye Anatomy** (upper left of the screen) or any item you choose
necessary Click **Activate Piggyback** (upper right of the screen)

Click **Misc Rx** (top menu bar)

Click **Chose any item** (from the left column)

Choosing **SKIN** is a major organ of elimination.

Click **Unconscious Search** (mid-left of the screen) in order to load unconscious information at bottom of the page

Click **Activate Piggyback** (upper right of the screen)

Click **Ear + etc** (top menu bar)

Click **Lymph** (on the left of the screen) or any choice you feel

appropriate Click **Activate Piggyback** (upper right of the screen)

Click **Iridology** (top menu bar)

Click **Iridology Choice Visible** (upper left of the screen)

Click **Lymph1** (from the blue column on the left)

Click **Activate Piggyback** (in the upper right of the screen)

Click **Dental Foci** (top menu bar)

Click **Auto Search** (on the left side, in the middle of the screen)

Click **Activate Piggyback** (in the upper right of the screen)

Click **Quadrant Repair** (top menu bar)

Select piggyback choices from the blue bars on the left of the screen

Double Click the area of the diagram you would like the piggyback signals to be directed as an area of additional focus.

Click **Sacred Geometry** (top menu bar)

Click **Activate Piggyback** in the upper right of the screen

Click the **Consciousness icon** on the bottom task bar to return to the main PULSAR program and allow the Iridology program to run in the background.

Chapter 5

Disease Probability and Dictionary

Disease Probability and Dictionary Training Program

Hint: *It is always best to open up Bodyviewer and Iridology panels first and thereafter open the Disease Dictionary program to prevent the programs from freezing.*

To open the Disease Dictionary from the Main Matrix screen:

Click *Subspace Programs* (top menu bar)

Click *Disease Dictionary* (library)

– Or –

During any treatment, on the “*LifeCare is Working*” panel:

Click *Disease Dictionary* (on the upper right)

PULSAR™ Universal Biofeedback System


NEW Programs Additional Subspace Programs Reporting Auto Focus Zap Retest Order of Remedies EPR Enhancements Library System Power Settings Audio Player Help

Body Viewer

VARHOPE Info Report

8603 100 Complex of mind
5370 100 Life Impulse® E
5073 100 Life Impulse® N
6393 100 Organic Matcha
10024 100 Life Impulse® C
4553 100 Life Impulse® E
5377 100 TurmeriCool, wit
4748 100 Life Impulse® E
6915 100 ProstaHEV, with pine phytosterols and pumpkin seeds, Life Care® | 7048 | - Generalized edema (dropsy), fluid metabolism problems, diuretic (kidney), gentle digestive strengthener.
3644 100 NeuroBooster Complex, with Brahmi extract and L-Tryptophan, Life Care® | 7046 | - Hysteria with chlorosis, chorea, spinal irritability, dysmenorrhea, bladder tenesmus, hysterical epilepsy, intense.
3904 100 A. K. K. Pure Tregado Oil | 4053 | - Infection from use of steroids or exposure to antibiotics.
4806 100 Reishi Mushroom Extract, 400 mg, Life Care® | 7051 | - Imbalance in patterns of waking and sleeping—such as insomnia, jet lag, time displacement
3857 100 Pure Ester C, with Omega-3, vitamin E and D, Life Care® | 7611 | - Fatty acid that helps form, on finger growth or deficiency of fatty acids.
6163 100 TurmeriCool, with piperine and curcumin, Life Care® | 7012 | - Nutrient rich algae. VITAMIN® heavy metal detoxification, energy
10880 100 Life Impulse® Vitamin C and Calcium | 860 | - Central Nervous System, Chinese Herbal Formula, anti-carcinogenic formula designed to aid in cellular integrity, immune support and detox
5625 100 Nio-Alergo, with black cumin oil, Life Care® | 763 | -) - listed in d. dictionary program, sugar regulation.
7370 100 Resveratrol Plus, Life Care® | 7490 | - Garlic concentrate in a base of vitamin C and chlorophyll to lower cholesterol + immune. BIOTICS E.N.B.
7798 100 Life Impulse® Lady's Desire | 740 | - Chronic fatigue, stomach tonic, NETS/HAN-Arzneimittel
9090 100 OPC Antioxidant, with grape seed extract, Life Care® | 7038 | - Diabetes mellitus - Sri Ramana Maharshi Natural Spagyric Remedies made for Duplication
5742 100 Alpha lipoic acid complex, with biotin, Life Care® | 7051 | -) - listed in d. dictionary program
5289 100 Stabilized oxygen molecules Life Care® | 891 | - B15 and other B vitamins from special beer yeast extract—oxygenation, I.Q.
10307 100 Mineral and vitamin complex, for pregnant women, Life Care® | 7071 | - Intestinal dysbiosis, digestive disorders & yeast infection
10626 100 Body cream, with Shea butter and BIO plants, Life Care® | 4320 | - remedy for cell regeneration and activation / cicatrization, burn trauma, fracture, activation of circulation
4659 100 Gel, with Snail extract, Life Care® | 4046 | - Shame in sexual pleasure—helps accept physical self, overcome physical loathing even acne, assists sexual pleasure, acceptance of one's physical body, openminded, hgdrr hgd2
9685 100 Life Impulse® ECO Aloe Vera Juice with cranberry flavor | 841 | - supports detoxification of cigarettes/tobacco/chemicals/ smoking / addiction—
9016 100 Life Impulse® Super Energy | 815 | - ASC Plus for male fertility support
10785 100 Life Impulse® Vitamin C and Calcium | 860 | - Natural product with anthocyanins is strong antioxidant / Cardiac, vascular, blood circulation, immune system, stomach, intestines and eyes
▶ 1462 105 Stabilized oxygen molecules Life Care® | 891 | - Methyl donor important for oxygenation—g vitamin, I.Q.
1997 100 Organic Cholesterol, Life Care® | 7071 | - Lipid-lowering, disease form of cholesterol, must bind to affect lungs, &
3213 100 Organic Pure Moringa, Life Care® | 7036 | - Bruises, accidents, sprains, bruised feelings, helps prevent strokes.
3248 100 ProBioLife, with germ complex, Life Care® | 8042 | - Whooping cough, respiratory affections (lung), laryngeal phthisis, tubercular glands.)
2748 100 Zinc, Life Care® | 7074 | - Can id acid alkaline imbalance, electrical sensitivity, mineral—pH
3126 100 Life Impulse® Vitamins for children | 870 | - Bones, splinters, imperfect assimilation and defective nutrition, diseases of bone, necrosis; suppurative processes, keloid growth.
2203 100 DiabForm, with momordica and alpha lipoic acid, Life Care® | 7530 | - Immune deficiency, manic depressive, diabetes (pancreas), deaf, gout, retardation, sex,]- sugar regulation
1756 100 Muscle Capsules, Devil's Claw, Life Care® | 4101 | - Inflammation of the bronchial tree interfering with lung air release.

Click *Disease Dictionary* (drop-down list)

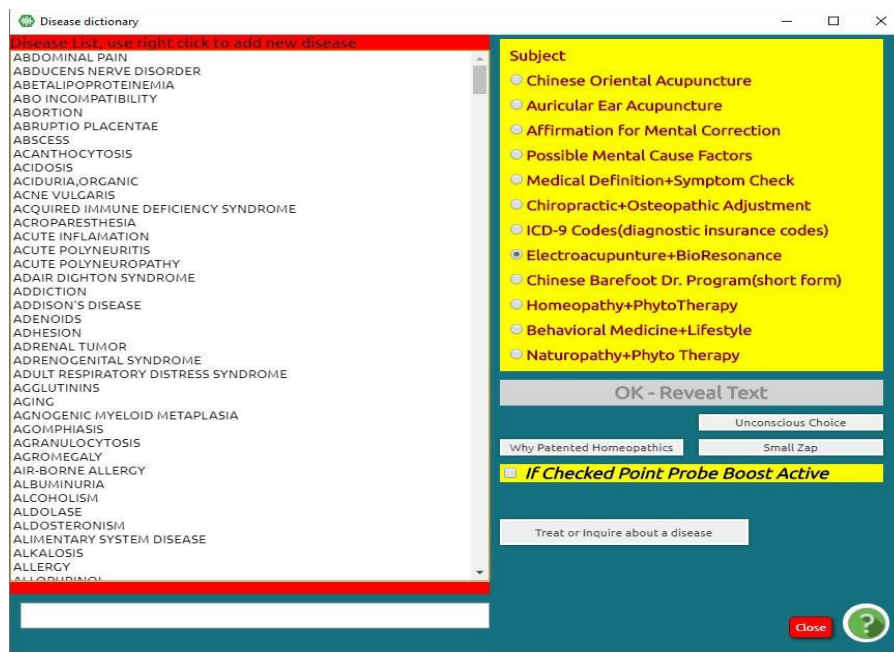


There are two ways to choose the items from the ***Disease Dictionary***:

- Search for the item of your choice in the alphabetic list on the left (for example “*Stress*” and click on it in order to select it.
- Click on the ***Unconscious Choice*** button and the program will choose the item according to the patient’s needs. The unconscious choice item will appear at the bottom of the list and will also be highlighted in the list.

Now you can select the type of therapy you wish to do in the yellow box on the right-hand side of the panel.

It is automatically set on ***Chinese Oriental Acupuncture***, but it is recommended to choose ***ElectroAcupuncture+BioResonance***.

Click ***OK – Reveal Text***

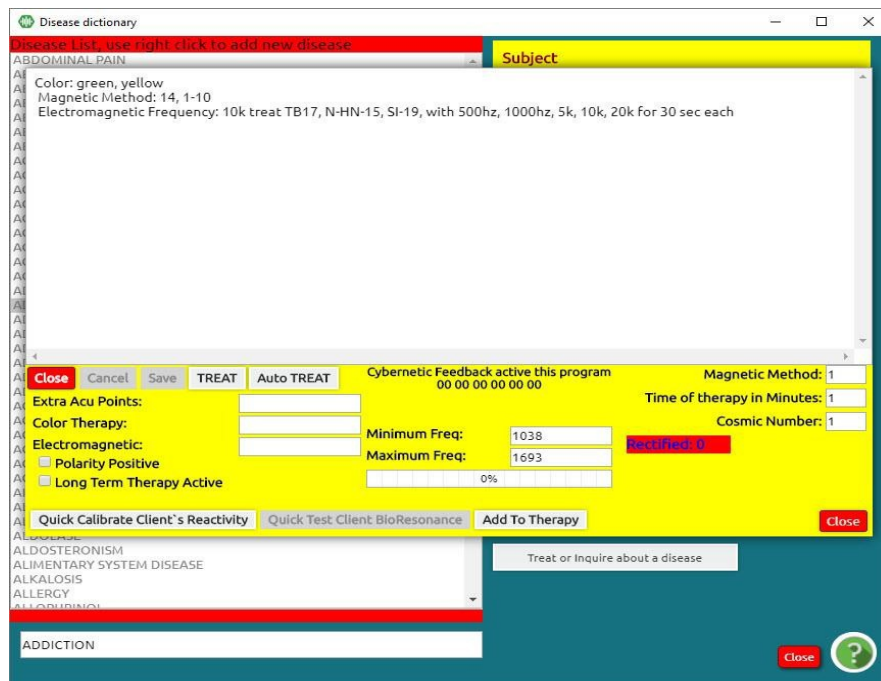


Click Quick Calibrate Client's Reactivity

This will do a quick calibration on the item of choice.

Click Quick Test Client BioResonance

This will do a quick test on the item of choice to show the necessary time and frequencies.

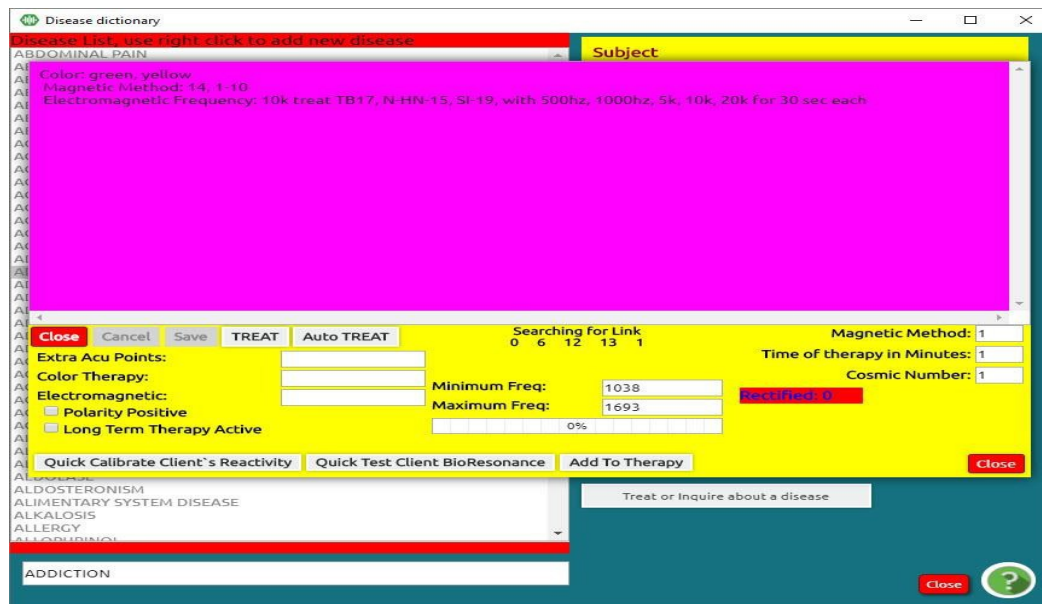

Click *Add to Therapy*

Here you can load additional information to add to your therapy. Load and close. Click *Auto Treat*

Click the *Consciousness* icon on the bottom task bar to let the *Disease Dictionary* training protocol run in the background.

Chapter 6

Close Out of Programs

Step 1.

- Close any opened program and return to the *Main Matrix* screen.
- Close the *Main Matrix* screen.
- This will bring you to the *Main Menu* screen

Step 2.

On your *Bottom Task Bar*, *Right-Click* and *Close*:

1. *Body Viewer*
2. *Iridology*
3. *Disease Dictionary* and any *Subspace Program* you may have opened

Note: Do not close *Consciousness* on your Task Bar, as this is the main ED.X and NUCLEUS program.

Shut Down Your Device

Click **OK**


To turn off your **ED.X and NUCLEUS** press the **On/Off Button** located at the front of your device.
The lights will go out.

- **Never remove the USB cable while the power is still on.**

Check Current Rectifications



Click **Check Current Rectifications** (on the lower right side of the Main Menu Screen) Click **OK** on the pop-up box
The blue Rectification panel appears.

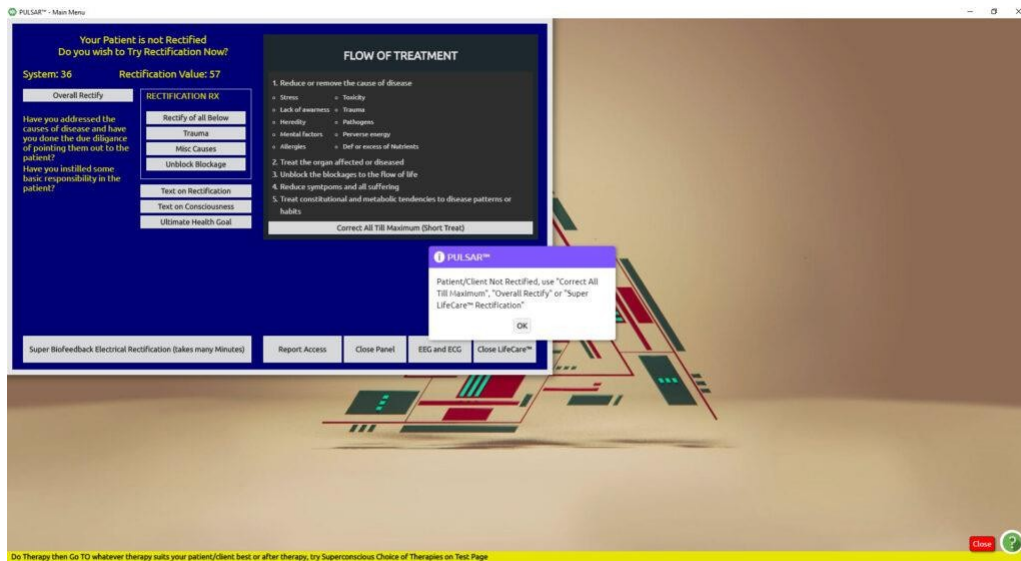

Click Correct All till Maximum (Short Treat)
The panel will change color. When complete, it will turn blue again.

Close Out of ED.X and NUCLEUS SOFTWARE

Right Click and close any items open on the bottom task bar (such as Iridology, Body Viewer, or Disease Probability and Dictionary)

Click **Close ED.X and NUCLEUS** on the lower right corner of the Check Rectification panel

Click **Yes** to *Are you sure you want to leave?*

Closing a session with your client

This is optional but a good idea when closing and finalizing your session with your client. Before you shut down your software and your device, from your Main Panel, after rectifying your session, click on TEST to return to your main matrix panel.



From your main matrix panel, go to the drop-down menu “Retest” and choose the options Test VARHOPE Improvement.

The screenshot shows the PULSAR™ Universal Biofeedback System interface. The main matrix panel displays a list of remedies with columns for 'no', 'value', and 'name'. The 'Retest' dropdown menu is open, showing options like 'VARHOPE', 'Info Report', and 'Expand Name'. The 'VARHOPE' option is selected, and the 'Info Report' option is also visible. The 'Expand Name' option is also visible. The 'Retest' dropdown menu is open, showing options like 'VARHOPE', 'Info Report', and 'Expand Name'. The 'VARHOPE' option is selected, and the 'Info Report' option is also visible. The 'Expand Name' option is also visible.

This process takes 1.5 minutes and basically recalibrates your client to get a rectification value on your initial VARHOPE numbers.

The screenshot shows the PULSAR™ Universal Biofeedback System interface during the RETEST process. The 'RETTEST' dropdown menu is open, showing options like 'VARHOPE', 'Info Report', and 'Expand Name'. The 'VARHOPE' option is selected, and the 'Info Report' option is also visible. The 'Expand Name' option is also visible. The 'RETTEST' dropdown menu is open, showing options like 'VARHOPE', 'Info Report', and 'Expand Name'. The 'VARHOPE' option is selected, and the 'Info Report' option is also visible. The 'Expand Name' option is also visible. The 'RETTEST' dropdown menu is open, showing options like 'VARHOPE', 'Info Report', and 'Expand Name'. The 'VARHOPE' option is selected, and the 'Info Report' option is also visible. The 'Expand Name' option is also visible.

As soon as the retesting process is done, a number appears in the midsection of the VARHOPE panel ‘Today’s Improvement Percent’.

PHAROS™ Universal Feedback System

new Programs Additional Subspace Programs Reporting Auto Focus Zap Reset Order of Remedies EPR Enhancements Library System Power Settings Audio Player Help

Personal message to patient started in red

How long How Often to treat

Original Items, and Image Visible

Make LifeCare™ Medicine Report

Add Report to Info File

Compare to Previous Results

Unconscious selection of next visit

Possible cause: allopathic drugs, or lack of Awareness

Possible most stressed organ: Parathyroid

Possible Blockage of Flow: Perverse Energy Exposure

Possible Diagnosis: Nicotine toxicity

Suggested Rx Plan: Heavy Metals Detox, or imponderable for Perverse Energy

Constitutional Type: Cerebral Venous, metabolic organ: Adrenal

Immediate Healing Goal:

Long Term Goal: take responsibility for my disease and let healing flow into my body

Treatment Plan: Addiction Therapy

Referral to Medical and other Professionals Plan: tell patients should consult with their medical doctor or ask for help

Patient Awareness to today's therapy: minimal verbal awareness

Nutritional dysfunction: Deficiency of: VITAMIN A, vit/mined with: Kidney

Behavioral Change: Exercise is essential for healthy life, some form of exercise can be always arranged, John Doe

Affirmation: I Love And Respect The Natural Forces Of The Universe.

Close

1268 95 Phosphorus For Digestion, Life Care® 1738 - Digestive system's enzyme, for digestion

412 95 Calcium and Magnesium Coral Complex, Life Care® 1774 - Detox remedy for industrial

1289 95 Stabilized oxygen molecules Life Care® 891 - B15 and other B-vitamins from special beer yeast extract-oxygenation, IQ.

1988 100 Life Impulse® Coral calcium 885 - Hormonal Balance Support

1280 100 Ginkgo Biloba, with lecithin, Life Care® 1781 - Menopause Hormone Balance / Female -feminine

10870 100 Life Impulse® Vitamin C and Calcium 860 - Lymphatic System, Chinese Herbal Formula, anti-carcinogenic formula designed to aid in cellular integrity, immune support and detox.

1766 100 Schizandra, Energizing, Life Care® 1780 - Hypertension, hypertony, NET/tyrosine/tyramine/serase "coffee" by blood types, could be reactive because of homeopathic reasons: anxiety, nervousness, insomnia, etc.

8152 100 Schizandra, Energizing, Life Care® 1780 - Nerves, psyche, sleep, pain, migraine, dizziness, PASCOS NATUROS-EDZIN-insomnia

8755 100 Life Impulse® Ace 852 - Protein maldigestion, high potency proteases, detox, allergies

ENDORPHINS - Excessive pain, inability to deal with emotional or physical pain--high

Click this button to find the Chiropractic connection for this compound on the Spinal Screen

Click this button to find the tissue connection for this compound on the Homeostatology Screens

Load Report Item

Empty below

Harmonic Coupling of 2 Items:

When Done with data Go to Therapy, or if therapy done, then go to specify screen

Info Report

Expand Name

RETEST

Data for John Doe

Soc Index: 1

Owner: DEMO

SYS POW: Normal

Percent Done: 0%

Regulatory dysfunction in EEG

Explore risk of energy management dysfunction

Explore risk of degeneration

Reactivity / Rectify

Reactivity: normal

Return to Main Matrix

Gut Foods

Guna Remedies

Herbaplus

Hormones

Horny Enzymes

Imponderables

Life Care™

Shaping (iron 85)

54%

Treat

Indu. Reaction

Test Item in Hold Tray above

Both Trays

Hold

Get/Enter

Load Report Item

Test Item in Hold Tray above

Get/Enter

Search for Item:

Search Above Tray

Clear Search

Close

Situation Hold:

Click on the third selection on the top right of this panel, Make LifeCare Medicine Report.

PHAROS™ Universal Feedback System

new Programs Additional Subspace Programs Reporting Auto Focus Zap Reset Order of Remedies EPR Enhancements Library System Power Settings Audio Player Help

Personal message to patient started in red

How long How Often to treat

Original Items, and Image Visible

Make LifeCare™ Medicine Report

Add Report to Info File

Compare to Previous Results

Unconscious selection of next visit

Possible cause: allopathic drugs, or lack of Awareness

Possible most stressed organ: Parathyroid

Possible Blockage of Flow: Perverse Energy Exposure

Possible Diagnosis: Nicotine toxicity

Suggested Rx Plan: Heavy Metals Detox, or imponderable for Perverse Energy

Constitutional Type: Cerebral Venous, metabolic organ: Adrenal

Immediate Healing Goal:

Long Term Goal: take responsibility for my disease and let healing flow into my body

Treatment Plan: Addiction Therapy

Referral to Medical and other Professionals Plan: tell patients should consult with their medical doctor or ask for help

Patient Awareness to today's therapy: minimal verbal awareness

Nutritional dysfunction: Deficiency of: VITAMIN A, vit/mined with: Kidney

Behavioral Change: Exercise is essential for healthy life, some form of exercise can be always arranged, John Doe

Affirmation: I Love And Respect The Natural Forces Of The Universe.

Close

1268 95 Phosphorus For Digestion, Life Care® 1738 - Digestive system's enzyme, for digestion

412 95 Calcium and Magnesium Coral Complex, Life Care® 1774 - Detox remedy for industrial

1289 95 Stabilized oxygen molecules Life Care® 891 - B15 and other B-vitamins from special beer yeast extract-oxygenation, IQ.

1988 100 Life Impulse® Coral calcium 885 - Hormonal Balance Support

1280 100 Ginkgo Biloba, with lecithin, Life Care® 1781 - Menopause Hormone Balance / Female -feminine

10870 100 Life Impulse® Vitamin C and Calcium 860 - Lymphatic System, Chinese Herbal Formula, anti-carcinogenic formula designed to aid in cellular integrity, immune support and detox.

1766 100 Schizandra, Energizing, Life Care® 1780 - Hypertension, hypertony, NET/tyrosine/tyramine/serase "coffee" by blood types, could be reactive because of homeopathic reasons: anxiety, nervousness, insomnia, etc.

8152 100 Schizandra, Energizing, Life Care® 1780 - Nerves, psyche, sleep, pain, migraine, dizziness, PASCOS NATUROS-EDZIN-insomnia

8755 100 Life Impulse® Ace 852 - Protein maldigestion, high potency proteases, detox, allergies

ENDORPHINS - Excessive pain, inability to deal with emotional or physical pain--high

Click this button to find the Chiropractic connection for this compound on the Spinal Screen

Click this button to find the tissue connection for this compound on the Homeostatology Screens

Load Report Item

Empty below

Harmonic Coupling of 2 Items:

When Done with data Go to Therapy, or if therapy done, then go to specify screen

Info Report

Expand Name

RETEST

Data for John Doe

Soc Index: 1

Owner: DEMO

SYS POW: Normal

Percent Done: 0%

Regulatory dysfunction in EEG

Explore risk of energy management dysfunction

Explore risk of degeneration

Reactivity / Rectify

Reactivity: normal

Return to Main Matrix

Gut Foods

Guna Remedies

Herbaplus

Hormones

Horny Enzymes

Imponderables

Life Care™

Shaping (iron 85)

54%

Treat

Indu. Reaction

Test Item in Hold Tray above

Both Trays

Hold

Get/Enter

Load Report Item

Test Item in Hold Tray above

Get/Enter

Search for Item:

Search Above Tray


Clear Search

Close

Situation Hold:


By doing this it will load information in the white lines of this panel.

NOTE: doing this process at the end after “Testing Improvement”, gives you information on what is still needed for the client to do, so you can discuss with your client what may be needed in addition to programming your next visit. If you do this process at the beginning of your session, it will give you an idea of what the client needs in the present session.

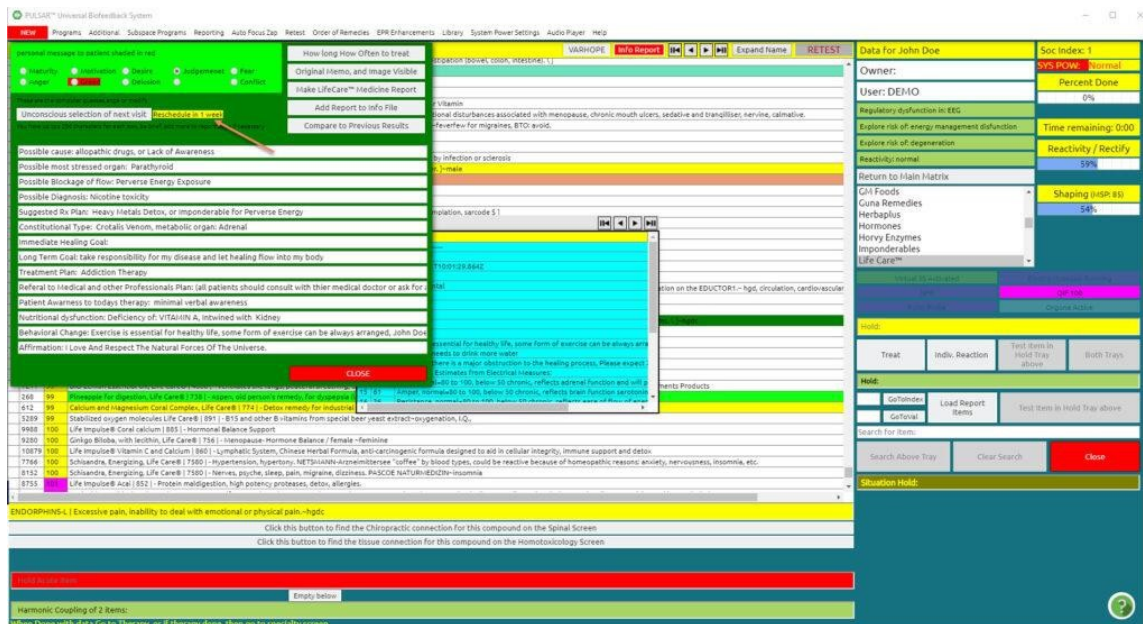

By pressing the “Unconscious selection of next visit”, the device estimates when the client should repeat their therapy for optimum results. Our example above, which gives us 1 week for the next visit, gives us information that the therapy done on them shall take approximately one week to absorb, and for optimum results should be repeated within 1 week again.

Other statistics on this panel is

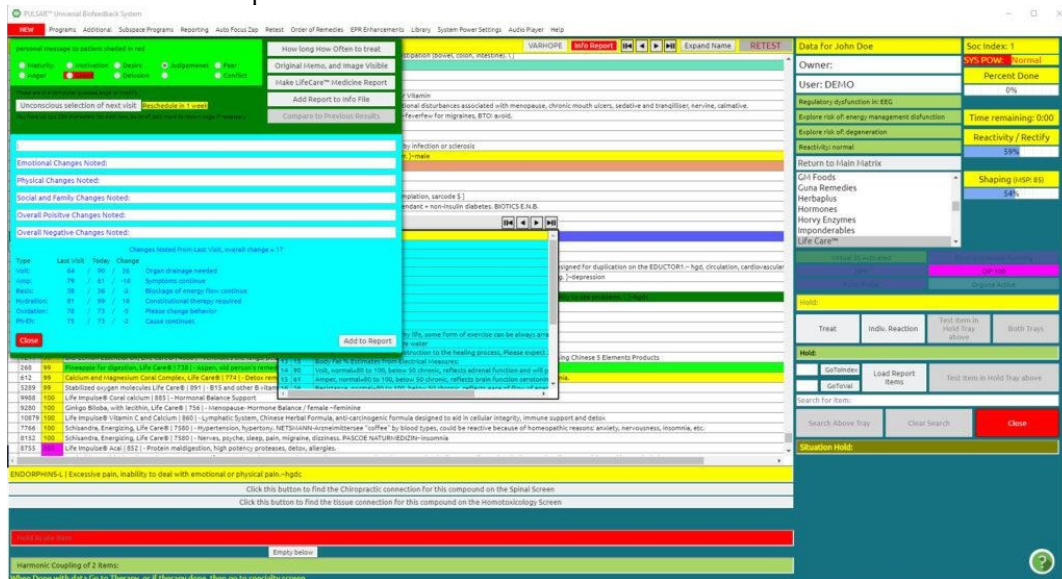

Compare to Previous Results, which is the last option on the list. This option compares VARHOPE values from the last visit.

Continue to close the software as previously explained in this chapter.

Chapter 7

Modifying the SOC Index of a Previous Client

Upon selecting your previous client as explained above, Calculate SOC will appear in bold.

Demographics

Patients Data Birth Time and Place Backup and Restore Demographics Reporting Focus on Symptoms not Diagnosis Help

Patients

Today is: 2022/10/14 Visit # 17 Number of Patients: 2

Name: John Doe Use Gender specific history: General complaint First visit:

Date of Birth: 10/14/1975 Male Female

Place of Birth: Amsterdam Inherited Disorder: Amount of negative: 0-100 0

Country: Netherlands

Address

City/Country/State: Zip Code: Insurance Company: New Patient Modify Patient Data Show Current Patient Delete Patient Cancel Hide Patients

Show Patients

Calculate SOC

Load Patient Picture

Delete Patient Picture

Add more Data to Report

Report Access

Bio rhythms

New Symptom Profile

Check Body Fat Percentage

Print Interview Help

Interview Inventory

Animal testing

Close

SOC INDEX: 1	John Doe to be Tested
GENICE 1000	
What part of your body do you do you dislike most?	
What possible things does your disease for you?	
How would your life change if your disease was gone?	

Click Calculate SOC

The SOC Index panel appears. You may alter any part of the questionnaire. When you finish modifying, by clicking CLOSE on the right, your modifications are automatically saved.

In order to load a previous client, after the password

From the **Main Menu** screen:

Click **Demographics**

As the client is already saved in the software database, you do **not** need to type the client's name in the **Patient Name** area, nor add any information

As we have already mentioned, for GDPR standards, the list of clients is hidden. If you wish to view the clients you can click on “**show patients**”

The screenshot shows the 'Demographics' screen of a software application. The top navigation bar includes 'Patient Data', 'Birth Time and Place', 'Backup and Restore', 'Demographics', 'Reporting', 'Focus on Symptoms not Diagnosis', and 'Help'. The main area is divided into a left sidebar with 'Patients' and 'Number of Patients: 1', and a central form area. The form contains fields for 'Name', 'Date of Birth', 'Place of Birth', 'Country', 'Address', 'City/Country/State', 'Street', 'Zip Code', 'Insurance Company', and 'Test Date'. A 'Show Patients' button is located at the bottom of the form. On the right side, there is a vertical menu with buttons for 'Calculate SOC', 'Load Patient Picture', 'Delete Patient Picture', 'Add more Data to Report', 'Report access', 'Biometrics', 'Add Symptom Profile', 'Check Body Fat Percentage', 'Print Interview Help', 'Interview Inventory', 'Animal testing', and a 'Close' button at the bottom.

Click Search for Past Client

This screenshot shows the same 'Demographics' screen, but with the search results displayed. A blue box highlights the 'Search for past Client' button. Below the form, a table lists the search results. The table has columns for 'Name', 'Birth Date', 'Birth Place', 'Country', 'Sex', 'Exam Count', 'SOC', 'Birth Time', 'Test Count', and 'Creation Date'. The first row shows a client named 'John Doe' with a birth date of '1980-10-20' and a creation date of '2023-10-11T10:01:29.864Z'. Below the table, there is a section for the 'SOC INDEX' and a 'John Doe to be Tested' section. The right sidebar remains the same as in the previous screenshot.

Name	Birth Date	Birth Place	Country	Sex	Exam Count	SOC	Birth Time	Test Count	Creation Date
John Doe	1980-10-20	Amsterdam	Netherlands	Male	1	1	1	0	2023-10-11T10:01:29.864Z

A prompt will come up telling you to “Please type the name of the client, then press Enter”

please type the name of the client and press enter


The details of your selected client will appear.

It is **recommended** that you review the **SOC Index** questionnaire at each session to record any changes that may have occurred since the previous session.

To record the new changes, please see **Modifying the SOC Index of a Previous Client**.

Modifying Birth Date

After you have selected the previous client, and in the case you need to edit either date of birth or place of birth Click **Modify Patient (Client)**

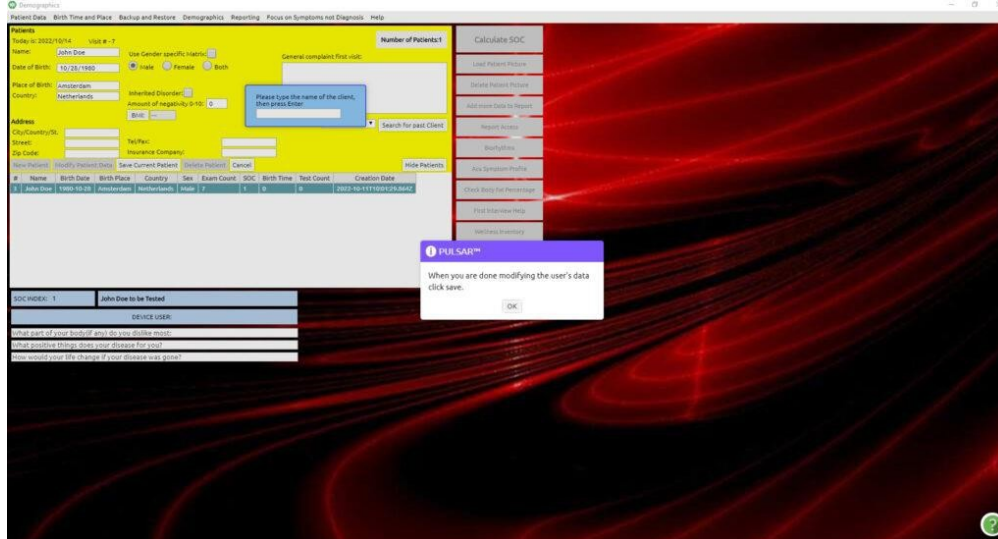

It shall give you a prompt “when you are done modifying the user’s data, click save” and press **OK**. You can now edit the birth information in the **Patient Data (Client)** panel. Click **Save Current Patient (Client)**.



Click **OK** on the questions.

Click **Close**.

Chapter 8

Auto Frequency Therapy / Auto VARHOPE

Step 1

From the Main Matrix screen:

Click Program

Click Universal Biofeedback Therapy

Click Auto Freq Therapy

Click Start

This may end up automatically switching to the Auto VARHOPE screen.

If so, continue with Step 2. If not, continue with Step 3.

Step 2

Click all of the items in the lower part of the panel that is relevant to the client.

This will put check marks in the yellow boxes.

(Some will have been selected automatically).

Click Emotion Stress > Electrical Imbalance (until Repaired for Today appears)

Click Auto Aggressive Immune Dysfunction (until Repaired for Today appears)

Click Cellular Vitality Disorder (until Repaired for Today appears)


Click Phase Aberration Present (until Repaired for Today appears)

Click Start Injury Repair (until the message reports Injury Repair Complete for Today)

Click Close

This will bring you back to the Auto Frequency Panel to complete the initial Training Protocol.

Note: Additional information may be typed in the tray next to Describe Injury Fully.

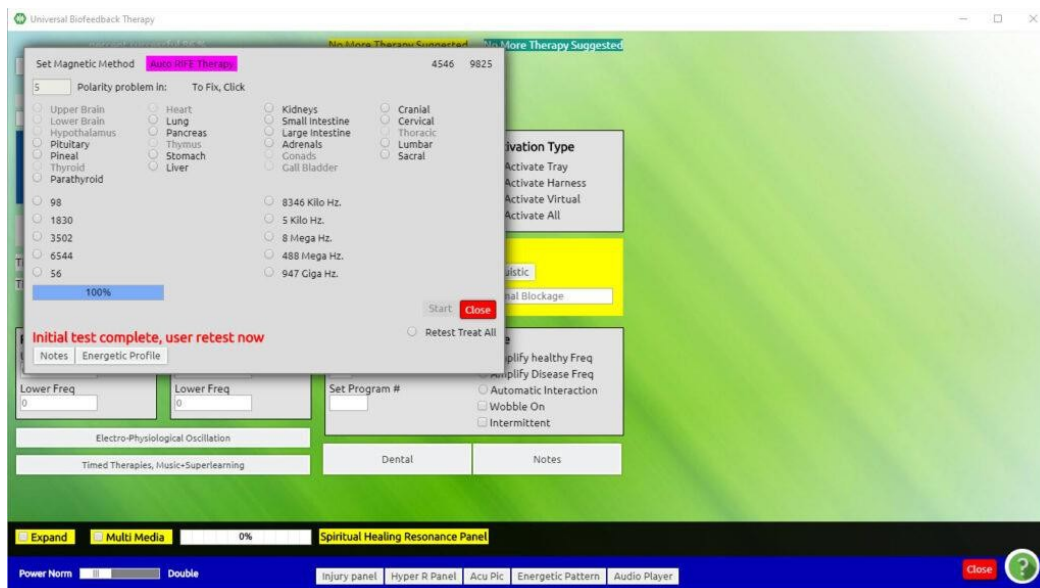

Note:

There is a possibility of up to four buttons that may appear in the mid-right section of the panel. Not all four buttons show up for every client, as the program chooses only what is appropriate. At times, only 2 or 3 of these will be available.

Step 3

Click **Retest Treat (Train) All**


Wait until the retest is complete.

Items will show up as either corrected or not corrected.

Click on any individual still **not corrected**. This gives another opportunity to attempt correction.

Click **Close**

Click **Auto VARHOPE** (and follow **Step 2**)

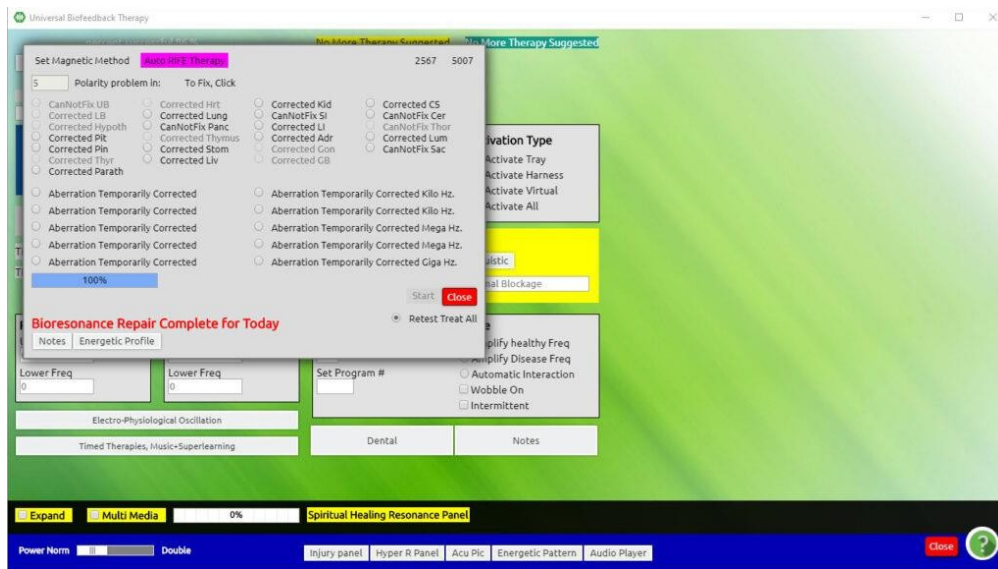

Auto Meridian / Hyper-Reactants

From the *Main Matrix* Screen:

Click **Program** (top menu bar)

Click **Universal Biofeedback Therapy** (drop-down list)

Click **Auto Meridian/Hyper-Reactants**

Click **Start Meridian Therapy** (bold black print indicates energetically blocked meridians)

This begins an automatic training protocol for any energetic disturbances in the acupuncture meridians.

Once complete:

Click **Retest All**. This will repeat an energetic flow through the meridians and initiate further training. A purple screen (the Acu Diag) will appear with information and indications of areas requiring further training.

Click **Close on the Acu Diag Panel**

Click **Close**

– Or –

If you wish to continue additional training protocols

Click **Acupuncture Pictures and Consciousness Form Test**

Click **Superconscious Search for Point** (in the lower left corner of this screen)

Double Click on the background of the picture or diagram of the panel to begin the training process. Continue clicking until the pop-up box shows **This Acu-Point or Points Corrected Energetically**.

You may choose to select and run the training process for several of the Superconscious Searches.

Click Close



Start Meridian Therapy (bold black print indicates energetically blocked meridians)
This begins an automatic training protocol for any energetic disturbances in the acupuncture meridians.



Once complete:
Retest All. This will repeat an energetic flow through the meridians and initiate further training. A purple screen (the Acu Diag) will appear with information and indications of areas requiring further training.



Close on the Acu Diag Panel

Close

Or –

If you wish to continue additional training protocols...

Acupuncture Pictures and Consciousness Form Test

Superconscious Search for Point (in the lower left corner of this screen)

DBL on the background of the picture or diagram of the panel to begin the training process. Continue clicking until the pop-up box shows This Acu-Point or Points Corrected Energetically.

You may choose to select and run the training process for several of the Superconscious Searches.

Close

Auto Scalar

This program is providing training protocols for the subtle energy centers known as the Chakra system. Like acupoints, chakras are specialized energy centers throughout our body where a unique form of subtle life-force energy is absorbed and distributed to our cells, organs, and tissues. The flow of this energy is greatly affected by our personality structure, emotions, and state of spiritual development.

From the **Main Matrix** screen:

Click **Program**

Click **Universal Biofeedback Therapy**

Click **Auto Scalar**

Click **Treat (Train) All**

This may take several minutes to train all the **Chakras**. When it is complete the panel should show the message **System in Correct Phase**.

Click **Close**



Auto Trivector

From the **Main Matrix** screen:

Click **Program** (top menu bar)

Click **Universal Biofeedback Therapy** (drop-down list)

Click **Auto Trivector**

Click **Use Electro Chelation Boost** (this is optional)

Click **Start**

Click **OK** on the message in the pop-up box. If complete, a box will appear **Trivector Complete for Today**

If the process is not complete a message stating so will appear. In order to initiate further Training, do the following:

Click **OK** on the message

Click **Prepare for Special Treatment (Training)**. Wait for the **Start** key to become bold again.

Click **Start**

When the process is complete:

Click **Close**



Note: It is possible a message may come up:

“Warning, There is a POSSIBLE Personality Oppression Spirit Superimposition Energetically Detected”.

There will be three options:


Continue

Correct Personality Aberration

Quit and exit Therapy

Click the **Continue** option in the message box (more in-depth information will be given on this in class)

Spiritual Resonance

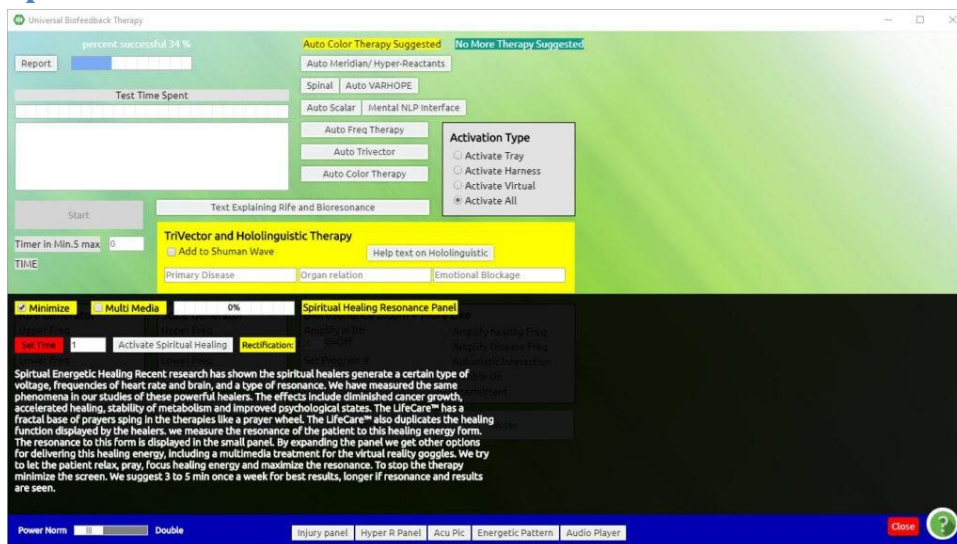

This is a beautiful training protocol that may encourage one’s awareness of the innate healing vibration potential of the body. It will also create a foundation whereby the client can now tolerate a higher rate of frequency. Notice that

the power (bottom left corner of the *Universal Biofeedback Therapy* panel) has now shifted to double after completing this training and will remain double for all the Auto training programs.

From the *Main Matrix* Screen:

Click **Program** (top menu bar)

Click *Universal Biofeedback Therapy* (drop-down list)

On the bottom left side of the *Universal Biofeedback Therapy* panel, in bright yellow:

Click **Expand** and it then changes to:

Click **Energetic Link**

The *Spiritual Resonance Training* panel appears.

Scroll the time up to one minute (or simply erase 0 and type the numeral 1) in the **Set Time** box and the **Activate Spiritual Healing** button appears.

Click **Activate Spiritual Healing (Training)**

The bar at the top of this panel shows the progression of the training as a percentage of completion.

Once the training has finished, the percentage now represents the rectification value.

If the program completes with anything less than 90% please repeat:

Click **Activate Spiritual Healing (Training)**

Trivector and Hololinguistic Therapy



Trivector and Hololinguistic Therapy

☐ Add to Shuman Wave Help text on Hololinguistic

Primary Disease Organ relation Emotional Blockage

From the *Main Matrix* screen:

Click **Program** (top menu bar)

Click *Universal Biofeedback Therapy* (drop-down list)

At the center of the screen, you will see the green

Therapy panel:

Type in the name of the known stressor(s) in the field of **Primary Disease**.

For example **Stress** (if you type in more than one is sure to double space between each one so the system will recognize them as separate items).

Type in which organs are related to that stress in the section **Organ Relation**. For example *All Organs and Tissues*

Type in the **Emotional Blockage**.

For example *Mental Retention of Toxins*

Click **Add to Schuman Wave**

This program runs throughout the session.

Two Suggested Training Programs

At the top of the screen, you will find two suggested training programs.

These are the suggested Superconscious choices. They are determined by the client's energy field (through Biofeedback).

The **yellow** one is the *primary Superconscious Choice*


Green is the *secondary Superconscious Choice*.

These colored buttons cannot be clicked on. To select one of the suggested training programs, you must click the appropriate button from the list below.

We will explore some of the main Training Programs accessible through this panel.

Note: In many of the training programs on this page, a particular area will appear highlighted with a specific choice, selection, or suggestion. These are often referred to as ***Superconscious Choices***. They are determined by the client's energy field (through biofeedback reactivity) and will guide you in training program options.


Auto Colour Therapy

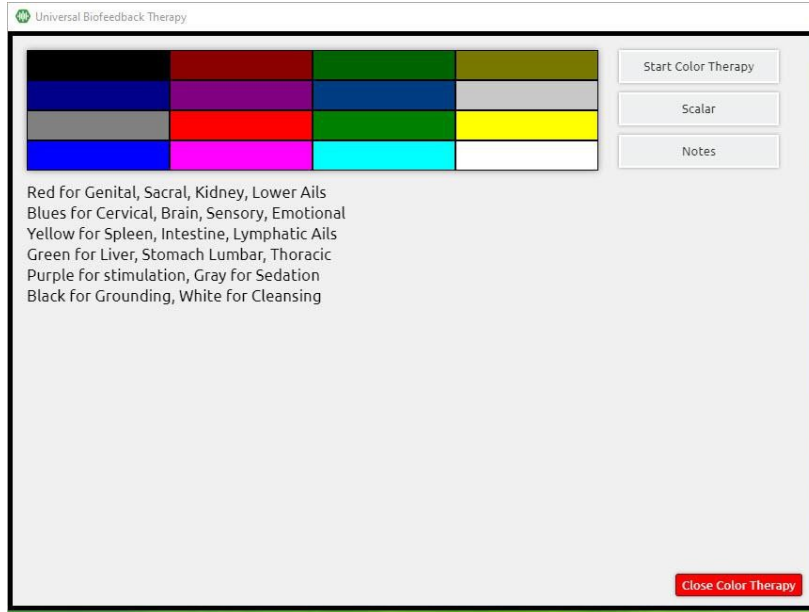

This program will flash through 256 different energetically harmonizing colors, eventually stopping on the most beneficial color.


From the *Main Matrix* screen:

- Click Program (top menu bar)
- Click Universal Biofeedback Therapy
- Click Auto Color Therapy
- Click Start Color Therapy
- Click Close Color Therapy

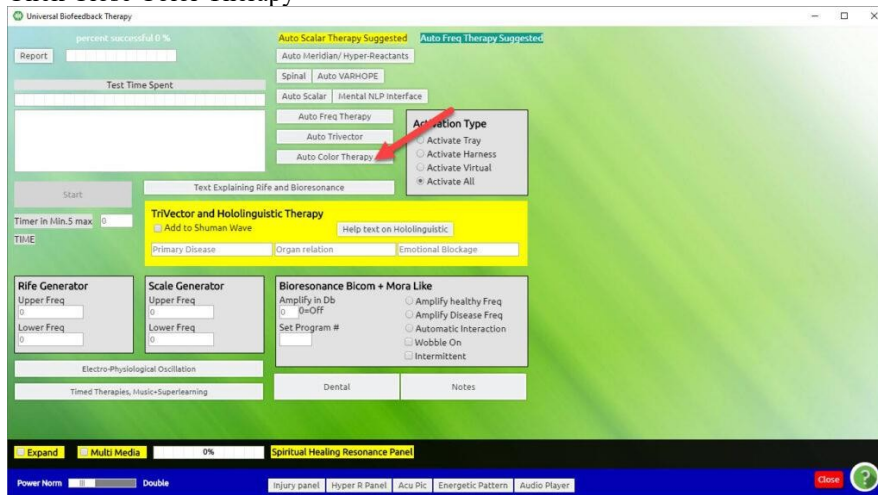

Chapter 9

Superlearning

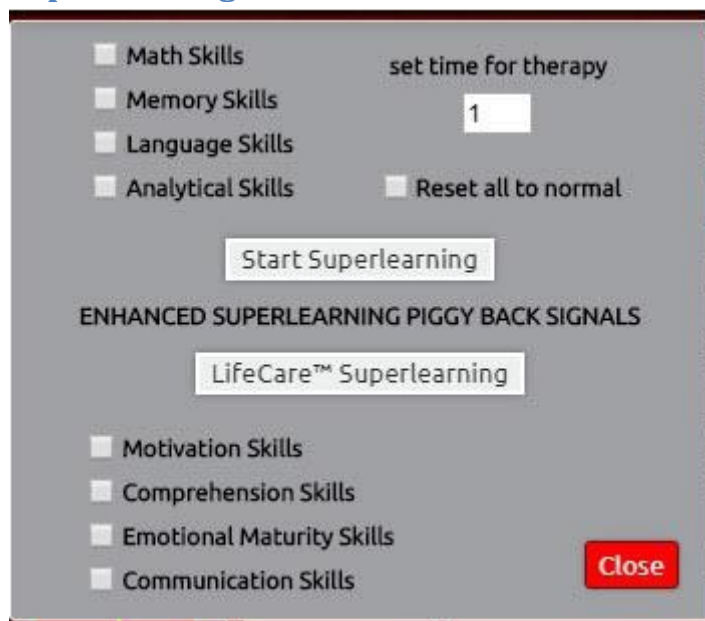

From the *Main Matrix* screen:

Click **Program** (top menu bar)

Click **Timed Cybernetic Feedback Therapy, Music, Superlearning** (drop-down list)

Click **Quick Check Phase**

Click **Enhanced Learning**

A grey panel appears.

Click appropriate selections

Set Time to 1 minute

(Time can be set for a maximum of 30 minutes).

Click **Start Superlearning**

Click **Harmonic Therapy (optional)**

Click **Homeopathic Superimposition** to add a specific learning frequency or suggestion/affirmation to be piggybacked along with the process.

Consider typing in the index number of some of the following Main Matrix items (making sure to DBL space between the numbers):

670: Memory (NV)

3325: LM 26 Memory

7011: 14 Essential Brain Nutrients

946: Remedy for Learning Disability

947: Remedy for Brain Circulation

Click **Quick Check of Phase** (with Harness, optional)

– **Or** –

Click **Virtual Quick Check of Phase** (without Harness using subspace, optional)

Click **Harmonic Therapy** (with Harness, optional)

– **Or** –

Click **Virtual Harmonic Therapy** (without Harness using subspace, optional)

Repeat this training protocol until all the bars on the graph are aligned to the left. Click ***ED.X and NUCLEUS Superlearning (optional)***

Note:

The ***Special ED.X and NUCLEUS or AFE Training Programs*** are found in different parts of the ***ED.X and NUCLEUS*** Software. They can be recognized by the word ED.X and NUCLEUS or AFE in the name.

These are only available when your ED.X and NUCLEUS device is connected to your computer.

They are different from other Training Protocols as they need to be activated from ***within*** the ***ED.X and NUCLEUS IS WORKING*** Box.

Click ***Enhanced Learning***

Click ***ED.X and NUCLEUS Superlearning***

Click and drag the ***ED.X and NUCLEUS IS WORKING*** box up, so you can see the full panel.

(To do this, you will need to click on the top border of the ***ED.X and NUCLEUS IS WORKING*** box, and, without releasing the click, drag it up so that the whole panel is fully visible).

Click ***Start Superlearning*** (in the lower right panel of the ***ED.X and NUCLEUS IS WORKING*** box)

To the bottom right you will notice 000. When it begins to rectify, these numbers will change. At first, you will see a rectification of 60 and it will then proceed to 100. This program is auto-focusing (self-selects the number of minutes needed) so you can expect it to take anywhere from a few to several minutes before rectifying.

Click ***End Superlearning*** (in the lower right panel of the ***ED.X and NUCLEUS IS WORKING*** box)

A pop-up box will appear for you to read.

Click ***OK***.

Click ***Close***.

Click ***OK*** on the pop-up box again.

Click ***Close***.


Timed Therapies

From the ***Main Matrix*** screen:


Program (top menu bar)

Timed Cybernetic Feedback Therapy, Music, Superlearning (drop-down)

Timed Treatments (Training)

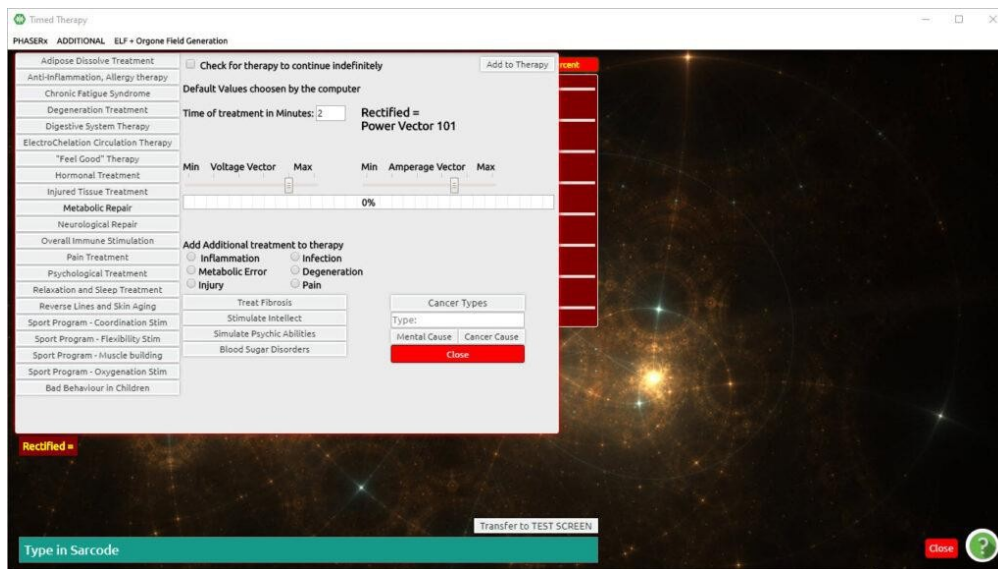

To adjust the time in minutes: click, hold and drag the slider.
 -the highlighted item in bold text (*Superconscious Choice*)
 By clicking on this any of the grey buttons, the training protocol will begin.
 Repeat any choice three times, rectifying to 85 or better.

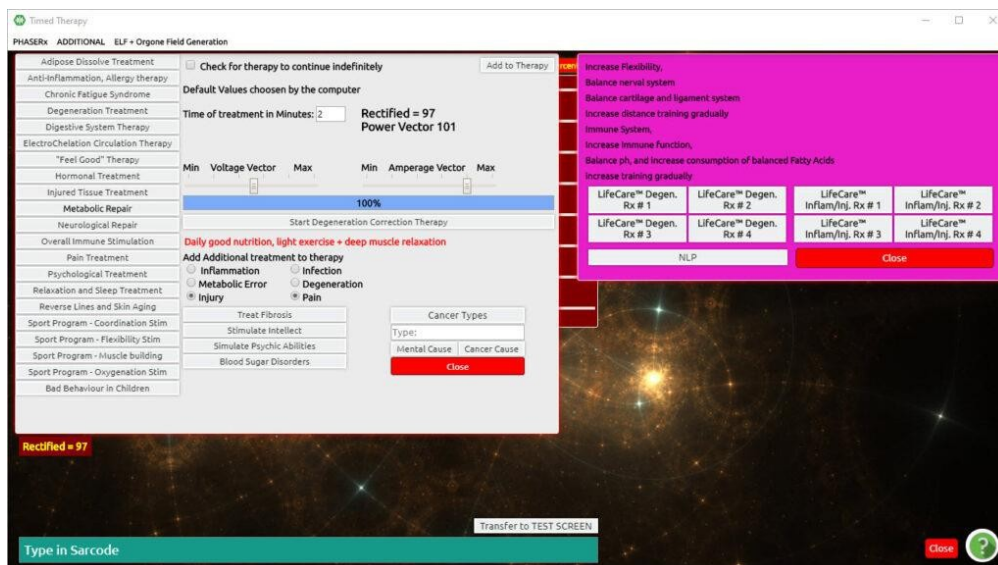


Chapter 10

ECG Fast Fourier Transform

From the *Main Matrix* screen:
 Click **Program** (top menu bar)
 Click **Universal Biofeedback Therapy** (drop-down list)

Click ***Electro-Physiological Oscillation***
The ***ECG Fast Fourier Transform Panel*** now appears.

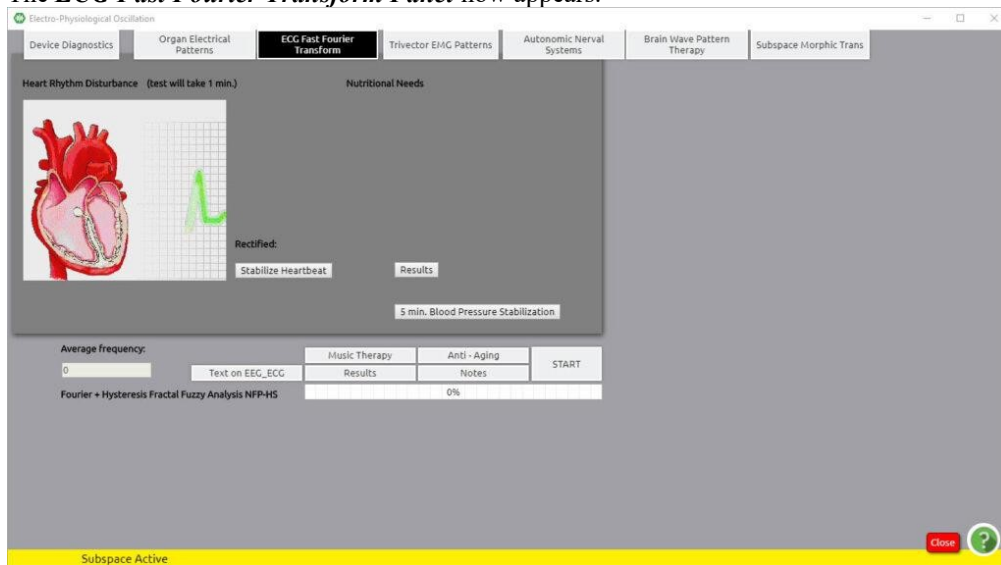

Click Start
Click OK
Click OK (finished)
Nutritional and organ-related energetic reactions are listed in a grey box.
Double Click the background of this grey box to remove it.
Click Stabilize Heart Beat – Repeat up to three times, attempting to rectify to 85 or better
Click Close


Superconscious Reduction Panel

From the ***Main Matrix*** Screen:

Click ***EPR Enhancements*** (top menu bar)
Click ***Client Superconscious Reduction Panel*** (drop-down list)



In the **List Disease** box, type in any concern of the client
 For example, **Resistance to change** or **Stress**
 Click **Start Superconscious Link**
 Click **Treat (Train) All**

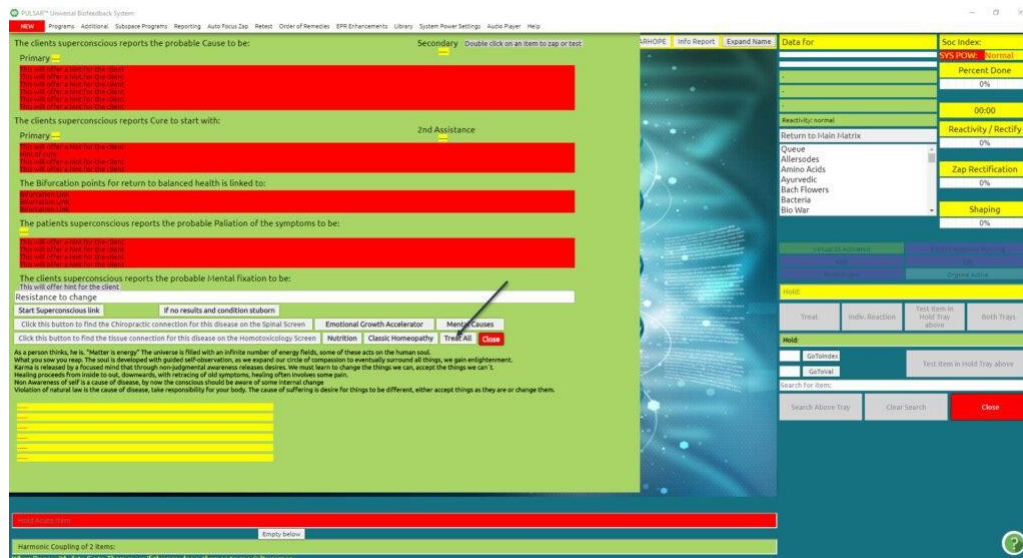

Autonomic Nerval System


From the **Main Matrix** screen:

Program (top menu bar)

Universal Biofeedback Therapy (drop-down list)

Electro-Physiological Oscillation

Autonomic Nerval System tab (top menu bar)
The **Autonomic Nerval System** screen appears.

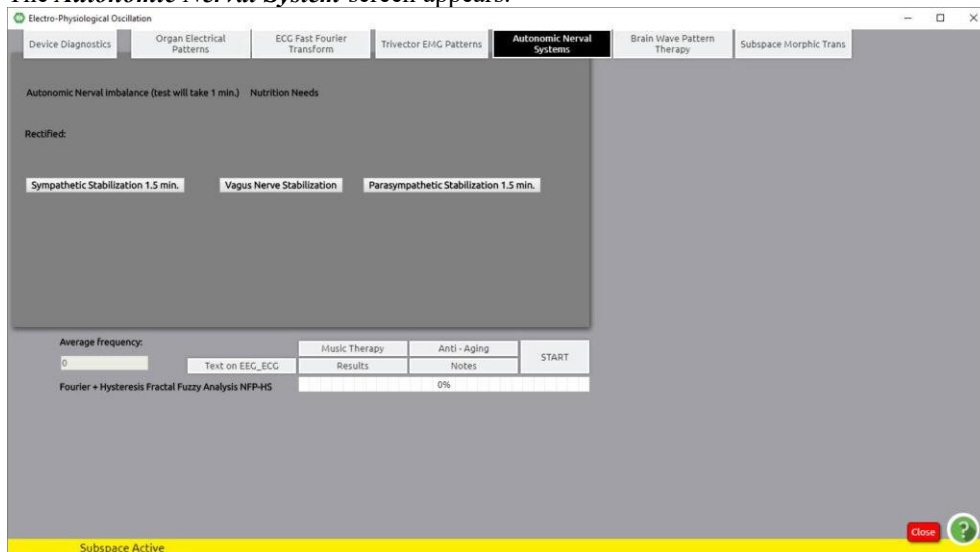

Start (will test for one minute)

Train Vagus Nerve Parasympathetic, or Sympathetic Stabilization

Depending on which is indicated in the grey report panel. Consider balancing each nervous system (for a complete relaxation experience). If the client is in good health and has no history of Vasovagal episodes

Vagus Nerve Stabilization (rectify to 85 or above)

Parasympathetic Stabilization (rectify to 85 or above)

Sympathetic Stabilization (rectify to 85 or above)

Close

CAUTION Possible Undesirable Side Effects:

Before training on the Autonomic Nerval System it is advised to ask your client: Do you have a history of Vasovagal Episodes, fainting, syncope, and/or generalized heart weakness? If the answer is yes, you need to proceed with caution, and be aware of the symptom of Vasovagal Episodes:

- A slight feeling of light-headedness after a session. Please note that this may be due to the client's own health such as low blood pressure, or suddenly standing up after being in a relaxed, reclined position during the session. If this feeling seems abnormal in any way, or of concern, then the client should see their primary healthcare practitioner.
- Through biofeedback electrical nerve stimulation; interaction with the CNS the ED.X and NUCLEUS may induce an autonomic nervous system cascade that can lead to a **Vasovagal** crisis. The client might sweat excessively or report nausea and dizziness. If this occurs, place a cool wet rag over their eyes gently with light pressure, tell the client to relax and breathe deeply, and wait about 5 to 10 minutes for the Vasovagal Storm to pass. In extreme cases, there might be **syncope (fainting)**.

VASOVAGAL Reaction

Biofeedback systems using skin resistance all use the stimulation of electricity.

Some clients may have a Vasovagal reaction when their system tries to shift from sympathetic nervous dominance to a more relaxing para-sympathetic dominance. If this happens too fast there can be a Vasovagal episode.

A Vasovagal episode. A vasovagal response or Vasovagal attack (also called Neurocardiogenic Syncope) is mediated by the Vagus nerve. When it leads to syncope (fainting) it is called Vasovagal Syncope, which is the most common form of fainting.

Symptoms

In clients with Vasovagal episodes, the episodes are typically recurrent, usually happening when the client is exposed to a specific trigger. The initial episode often occurs when the person is a teenager, and then recurs in clusters throughout their life. Prior to losing consciousness, the individual frequently experiences warning symptoms

such as lightheadedness, nausea, sweating, ringing in the ears (tinnitus), uncomfortable feeling in the heart, weakness, and visual disturbances such as lights seeming too bright, fuzzy, or tunnel vision. These last for at least a few seconds before consciousness is lost (if it is lost), which typically happens when the person is sitting up or standing. When sufferers pass out, they fall down (unless this is impeded); allowing the person to wake up. The autonomic nervous system's physiologic state leading to loss of consciousness may persist for several minutes, so:

- If sufferers try to sit or stand when they wake up, they may pass out again;
- The client may be nauseated, pale, and sweaty for several minutes.

Vasovagal syncope is rarely life-threatening in and of itself but is mostly associated with injuries from falling while having an episode.

The Autonomic Nervous System has two parts 1) the Sympathetic for stress (fight-or-flight) also known as the Thoracic Lumbar, and 2) the Parasympathetic for immunity and digestion, also known as the Cranial Sacral Nerval system. To stimulate the Parasympathetic system apply cold compresses for five minutes to the upper neck and sacrum, and then alternate with warm compresses for five minutes.

What to Do

The method of dealing with *Vasovagal syncope* focuses on avoidance of triggers, restoring blood flow to the brain, and measures that interrupt the process; such as cold and hot compresses. Also to avoid physical injury due to fainting.

Because *Vasovagal syncope* causes a decrease in blood pressure, relaxing the entire body as a mode of avoidance isn't favorable. The client can cross their legs and tighten leg muscles to keep blood pressure from dropping.

Before known triggering events, the client may increase the consumption of salt and fluids to increase blood volume. Sports and energy drinks may be particularly helpful.

Clients should be educated on how to respond to further episodes of syncope, especially if they experience pre-warning signs: they should lie down and raise their legs, or at least lower their head to increase blood flow to the brain.

If the client has lost consciousness, they should be laid down with their head turned to the side. Tight clothing should be loosened. If the inciting factor is known, it should be removed if possible (for instance, the cause of the pain).

Wearing graded compression stockings may be helpful.

There are certain orthostatic training exercises that have been proven to improve symptoms in people with recurrent Vasovagal syncope.

How do I use my Biofeedback device with a client who has a history of Vasovagal episodes?

Do NOT use maximum settings for too long and stay at the safe calibrated levels. This will reduce the possibility of a Vasovagal episode. Always ask if there is a history of Vasovagal episodes, most often fainting or cold sweats with heart palpitations. If so always use low-calibrated settings. Do not try to push therapy after therapy on the client with a Vasovagal episode history. Please always keep a watchful eye on the client who is getting therapy on electrical stimulation. Be prepared to respond. Respond by cold water on the face gently and/or gentle pressure over the eyes with a cold rag. Push with your knuckle on the acupuncture emergency spot above their upper lip in the cleft under their nose (note, the point is on the bone so, therefore, push hard to apply pressure to the bone). Lie the client down, and be cautious that they do not stand up too fast. Reduce stress and wait for ten minutes. It will most likely pass. Do not let them leave until they are better.

Call 911 if they pass out for more than a minute, have pain or vomiting, or anytime your better judgment tells you to call 911. Report the fainting spell to the 911 crew.

Chapter 11

Bones, Nerves, Circulation, Lymph, Digestion

These training programs can all be operated in the same way as *the Muscles* above.

From the *Main Matrix* screen:

Click **Program** (top menu bar)

Click **Spinal Muscular Re-Education, Sarcode Feedback**

Click **Bones, Nerves, Circulation, Lymph, Digestion, Eye, or Sinus**

index	value	name
1218	96	ACID PHOSPHATASE Ids problems with blood levels of acid phosphatase, ids bone liver and or kidney (bladder, urethra, ureter) disease.]
8692	97	Vitamin E Forte, soft capsules, Life Care® 7063 - Poor calcium absorption metabolism, bone changes.
10320	97	Recovery gel, „Knitbone“, with comfrey and BIO plants, Life Care® 4450 - Bone fractures, slow healing bone fractures & bone pain
5078	97	Bone and Teeth Homeopathy Japan]
3560	98	CRANIAL BONES FACIAL PALANTINE sarcode \$]
3732	108	Life Impulse® Joint Elixir 880 - The gatherer stone, gathers and absorbs vital nutrient energy, helps teeth and bones, crystal
3558	110	CRANIAL BONES FACIAL LACRIMAL sarcode \$]
5804	111	MARFAN'S SYNDROME A hereditary condition of connective tissue, bones, muscles, ligaments & skeletal structures. patients are taller than average for age and family,]- listed in d. dictionary program

Broken Bone and Fracture Repair

● Treat Circulation overall ● Lower Right ● Upper Right ● Upper Torso ● Heart
● Lower Left ● Upper Left ● Low Torso ● Head + Neck

Add an Additional therapy superimposed on the treatment

Set time in minutes: 1

Start treatment

Rectified : 000

Quadrant Tension: 000

Sarcode Resonance: 000

Tissue Disturbance: 000

Pain Types

- ☒ Heat Improves
- ☐ Cold Improves
- ☒ Touch Improves
- ☒ No Touch Improves
- ☒ Motion Improves
- ☐ No Motion Improves
- ☐ Acid Pain
- ☐ Alkaline Pain
- ☐ Electromagnetic Pain
- ☐ Toxic Pain
- ☒ Referred pain
- ☐ Geographical Pain
- ☒ Emotional or Psychic

Close

Click on an item to highlight the line

Click in the little box next to each additional item that would be appropriate to add as a piggyback signal—[These are found in the section *Add an additional therapy superimposed on Treatment (Training)*].

Click **Start Treatment (Training)**

– Or –

Choose a special focus by clicking the desired area's radio button (Lower Right, Upper Right, Upper Torso, Heart, Low Torso, Lower Left, Or Head and Neck). Only one choice can be selected at a time.

The timer can be left on the default of 0 or set for 1 minute.

Click **Start Treatment (Training)**

Click **OK** in the *Treatment (Training) Over* the box

Repeat up to three times, attempting to rectify to 85 or better

Click **Close**

Variant Treatment

index	value	name
2545	95	Thyroidaltonin Bone and calcium regulator, lds with thyroid or calcium diseases.
1218	96	ACID PHOSPHATASE lds problems with blood levels of acid phosphotase, lds bone liver and or kidney (bladder, urethra, ureter) disease.]
8692	97	Vitamin E Forte, soft capsules, Life Care® 7063 - Poor calcium absorption metabolism, bone changes.
10320	97	Recovery gel, „Knitbone“, with comfrey and BIO plants, Life Care® 4450 - Bone fractures, slow healing bone fractures & bone pain
5078	97	Bone and Teeth Homeopathy Japan]
3560	98	CRANIAL BONES FACIAL PALANTINE sarcode \$]
3732	108	Life Impulse® Joint Elixir 880 - The gatherer stone, gathers and absorbs vital nutrient energy, helps teeth and bones, crystal
3558	110	CRANIAL BONES FACIAL LACRIMAL sarcode \$]

Broken Bone and Fracture Repair
Show picture

Treat Circulation overall
Lower Left

Lower Right
Upper Right

Upper Left
Upper Torso

Head + Neck
Heart

Add an Additional therapy superimposed on the treatment:
Set time in minutes: 1
Start treatment

Rectified : 50
Quadrant Tension: 25
Sarcode Resonance: 3
Tissue Disturbance: 3

Pain Types
☒ Heat Improves
☐ Cold Improves
☐ Touch Improves
☐ No Touch Improves
☐ Motion Improves
☐ No Motion Improves
☐ Acid Pain
☐ Alkaline Pain
☐ Electromagnetic Pain
☒ Toxic Pain
☒ Referred pain
☐ Geographical Pain
☐ Emotional or Psychic

☒ Treat Facia in area
☐ Treat Inflammation
☐ Treat Metabolic dysfunction
☐ Treat Infection
☐ Treat Clonus or spasm in area
☒ Treat Injury in Area
☒ Treat Degeneration
☐ Treat Pain in area
☒ Treat Fibrous tissue in the area
☒ Use Electro Chelation Boost Therapy for circulation
☒ Vascular Endothelial Growth Factor, energy that stimulates growth of new blood vessels

Close

Spinal Muscular Re-Education, Sarcode Feedback

Spinal Muscular Re-Education, Sarcode Feedback

Please Adjust the Spine before Testing, the system will show remaining Foci
Reset Colors
Dbl Click on Dermatome Area to

Cervical
C1
C2
C3
C4
C5
C6
C7

Thoracic
T1
T2
T3
T4
T5
T6
T7
T8
T9
T10
T11
T12

Lumbar
L1
L2
L3
L4
L5

Cranial Sacral
Cranial
Sacral
Occiput

Sensitivity
Maximum
Normal
Minimal

Test + Treat Trivector Energy Flow
Neuro Emotional Complex NEC
Original Sarcode Results
Nutrition Profile
Cranial Sacral
Timed therapies
Sport Programs
Biorhythm
Pain
Dental
Muscles
Bones
Nerves
Circulation
Lymph
Digestion
Eye
Sinus
NLP
Allersodes
Make Report of Spinal Program
VARHOPE GRAPH
Transfer to Test Screen

Cervical
Thoracic
Lumbar

Cervical
Thoracic
Lumbar

Close Spinal


From the **Main Matrix** screen:

Click **Program** (top menu bar)

Click **Spinal Muscular Re-Education, Sarcode**

Feedback Click **Test and Treat (Train) Trivector Energy**

Flow Repeat until all vertebrae are marked as corrected.

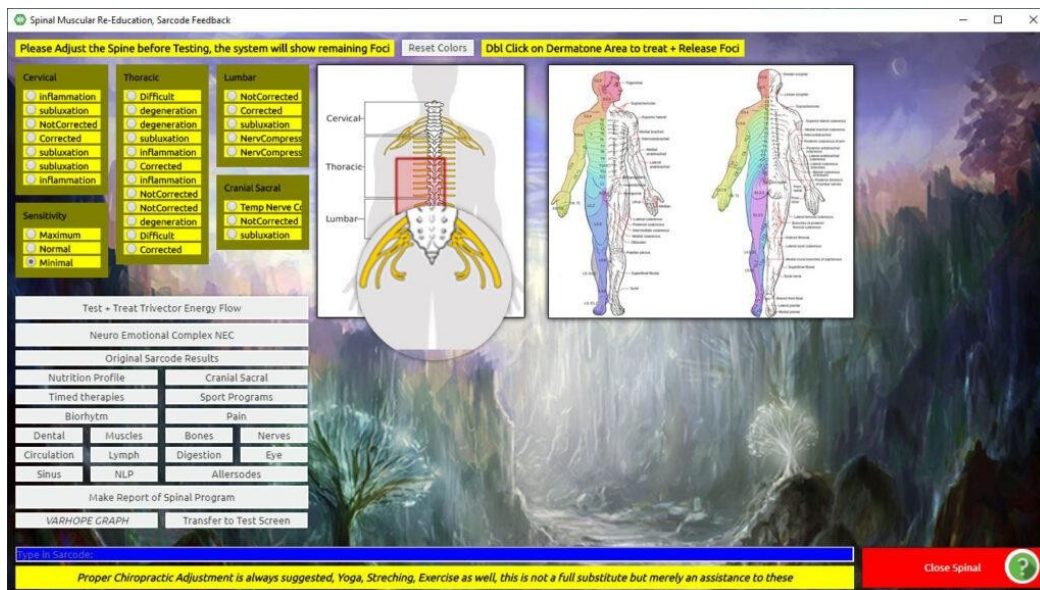


Click **OK, YES** or **No** (to any pop-up boxes that may come up during the process)

From the **Main Matrix** screen:

Click **Program** (top menu bar)

Click **Spinal Muscular Re-Education, Sarcode Feedback**

Click **Muscles**

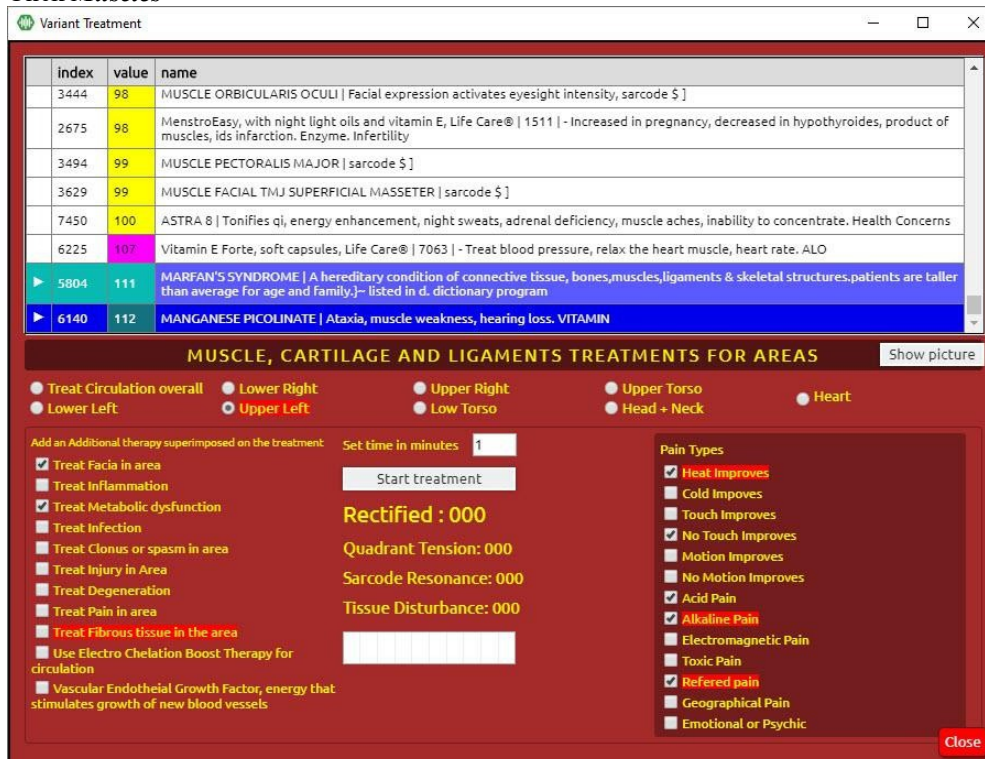

The program will isolate a particular muscle sarcode, ligament, or tendon based on reactivity; this will be indicated by a small arrow in the upper panel listing the matrix items.

Click on this item to highlight the line

Click in the little box next to each additional item that would be appropriate to add as a piggyback signal—[These are found in the section *Add an additional therapy superimposed on Treatment (Training)*].

Click *Start Treatment (Training)*

– **Or** –

Choose a special focus by clicking the desired area's radio button (Lower Right, Upper Right, Upper Torso, Heart, Low Torso, Lower Left, Or Head and Neck). Only one choice can be selected at a time.

The timer can be left on the default of 0 or set for 1 minute.

Click *Start Treatment (Training)*

Click *OK* in the *Treatment (Training) Over* the box

Repeat up to three times, attempting to rectify to 85 or better

Click *Close*

– **Or** –

Click anywhere on the upper panel and begin to scroll up or down, using the up/down arrow keys on your computer keyboard until you find the specific muscle sarcode you are seeking.

Click each additional item that would be an appropriate piggyback signal

Click to highlight the muscle you want to train

Click *Start Treatment (Training)*

Click *OK* in the *Treatment (Training) Over* the box

Repeat up to three times, attempting to rectify to over 85 or better

Click *Close*.

Timed Biofeedback Therapies

From the *Main Matrix* screen:

Click *Program* (top menu bar)

Click *Spinal Muscular Re-Education, Sarcode Feedback*

Click *Timed Therapies*

Often there is a **Superconscious** choice highlighted in red on this panel. If so:

Click the red highlighted item (the timer default is set to one minute)

Click *Start* (in the lower right of the screen-rectify to 85 or above)

If there is not a red highlighted item, unclick any items that may have been checked in the yellow boxes as piggybacks and

Click *Repeat Unconscious Search*

Once the red highlight shows:

Click the red highlighted item

Click *Start* (rectify to 85 or above)

Spinal Muscular Re-Education, Sarcode Feedback

☐ Adrenal dysfunction Stim	☐ Feel good therapy	☐ Memory Stim	☐ Pituitary
☒ Brain Stim	☐ Frontal Lobe of Brain	☐ Metabolic Stim	☐ Restore Prostate Health
☐ Cerebellum	☐ Heart Stabilizer	☐ Mucous Membranes	☐ Reticular Formation
☐ Cerebral spinal fluid	☐ Hypothalamus	☐ Muscle Tone	☒ Sex System Stim
☐ Circulation Stimulation	☐ Immune Stim	☐ Nerval pulse Stim	☐ Sinus
☐ Congested Lymph Stim	☐ Intestinal Stim	☐ Occipital	☐ Skin healing pulse
☐ Connective tissue Stim	☐ Kidney	☐ Spleen	☐ Sluggish Liver
☐ Digestive Stim	☐ Larynx	☐ Parathyroid	☐ Stimulate fertility Male
☐ Ears	☐ Limbic System	☐ Parietal	☐ Teething Pain
☐ Eyes	☐ Lung therapy	☐ Pineal	☐ Temporal
			☐ Thymus Gland
			☐ Thyroid dysfunction Stim
			☐ Vagus Nerve

☐ Build Eye-Hand Coordinate	☐ Reduce Paranoia-Distrust	☐ Treat Pain in the Area
☐ Build Eye-Foot Coordinate	☒ Reduce Hormone Imbalance	☐ Relieve Inflammation in Area
☐ Treat Dryness in the Area	☐ Increase Mental Focus	☐ Treat Metabolic dysfunction
☐ Treat Damp in the Area	☐ Release Emotional Stress	☐ Treat Injury in the Area
☐ Treat Cold in the Area	☐ Treat Phlegm in the Area	☒ Treat Infection in the Area
☐ Treat Heat in the Area	☐ Treat Energy Stagnation	☐ Treat Degeneration in the Area
☐ Treat Wind in the Area	☐ Treat Weak Chi	☒ Treat Fibrous Tissue
Wholistic Organ Health Restorer	☐ Adjust Trivector Disorders	☐ Send Detox Pulse into Area

Enter Additional Data for Therapy

Link to: A Reticular Formation Primary Connection: Nerval ,

Isode:

Allersode: NETTLE | Pollen Extract (ALR), POLLEN AMARANTHUS (ALR) | Allergy sensitivity or toxic reaction, avoid and desensitize, PAINT SOLVENTS | phenc

Sarcode: Link to: A Parathyroid

Denial of Fear

Energetically Eliminate these cells from the area of treatment

☐ Atopose ☐ Mucous ☐ Mast ☐ Degenerate ☐ Inflamed ☐ Mutated ☐ Necrotic ☐ Stimulate Acceptance of Tissue

1 minutes

20 min

40 min

Rectified: 52



Chapter 12

Hormonal Training Program

From the **Nutrition** panel:

Click **Male** or **Female Problems** (depending on the client, the button will show as male or female)

Click **Do you wish to use the Questions? Yes or No**

Selecting **Yes** will provide a questionnaire that will fine-tune the training protocol.

Click **Hormone Balance** – rectify to 85 or above

Click **Ovaries Balance** or **Prostate Balance** – rectify to 85 or above

Click **Close**.

QED TRIVECTOR SPINAL FOCI ANALYSIS - Nutrition Profile

Notes Report

index	value	name
2706	91	Ornithine Decarboxylase Ids with pituitary disorder. Enzyme.]
2692	91	Histaminase Enzyme for histamine regulation, Ids histamine concern
1423	92	HISTIDASE Enzyme involved in anti-inflammation, Ids inflammation c
599	92	Organic Spirulina, Life Care® 7013 - Ids with some anemia or vitam
9313	92	Co-enzyme CoQ 10 Carnitine Vitamin E / for energy regulation - HZ
2701	92	Lytic Enzymes Enzyme used in blood and spleen to destroy red bloo
796	93	PHOSPHATYDL CHOLINE Lecithin for brain function and cholinergic
2666	93	Choline Acetyltransferase Ids autonomic (parasympathetic neural sy
940	93	PULMO LIQUITROPHIC (DR) Combo remedy to assist in lung repair.
8806	93	LIPID DEPOSITION DISEASE Gaucher's disease (glucoceramide), depc
8807	94	LIPID DEPOSITION DISEASE Niemann-Pick disease (sphingomyelin), c
245	94	PROPHYRA Liver, sea plant, liver remedy, Ids amino acid deficiency,
377	94	APIOL (4X) Phenol compound used for allergy desensitization.
2732	94	Xanthine Oxidase Ids milk allergies or toxic bowel (colon, intestine), r

Hormonal Imbalances = Sugar, Guilt, Vaccination, T3, Anger, Oxytocin,
Disease factors = Unmanaged emotional stress
Other Conflicts with : Thyroid, Mental Factors,
Intwined:
Intwined with Circulation
Intwined with Allergy
Intwined with Bacteria
Intwined with Environment
Therapy: Prostate Therapy Stress Reduction, Prostate Therapy

Sort Table by: Value Name Index

More Info On Natural Cooking on the NET Constitutional

Increase Enzyme Therapy Stabilize Enzyme Therapy Increase Hormone Therapy Stabilize Hormone Therapy Nutri/Enzy/Horm Stabilize

Close ?

Internal Enzymes

QED TRIVECTOR SPINAL FOCI ANALYSIS - Nutrition Profile

Notes Report

119	Vitamin A family	108	Protein Digestion
74	Vitamin B complex family	106	Carbohydrate Digestion
87	Vitamin G complex family	76	Fat Digestion
60	Vitamin C family	68	General Digestion
74	Vitamin D family	99	Internal Enzymes
70	Vitamin E family		Some Data on Risk menu
111	Vitamin F (fatty acids)	82	Minerals
72	Vitamin K family	96	Amino Acids
81	Vitamin U family (cq Q)	79	Oriental Herbal Combos
76	Biochemic Tissue salts		

Calculate, strike this key first: Go to Progressive Homeopathy list

Basic Homeopathy From Test Scan Make master nutritional report

Access Reports + Texts Good foods Virtual Doctor

Insert Disease:

Male problems Stress and Cortisol

Transfer To Test Screen

Intwined:
N/A
N/A
N/A
N/A
N/A

More Info On Natural Cooking on the NET Constitutional

Increase Enzyme Therapy Stabilize Enzyme Therapy Increase Hormone Therapy Stabilize Hormone Therapy Nutri/Enzy/Horm Stabilize

Close ?

From the **Nutrition** panel:

Internal Enzymes

-on the highlighted item – rectify to 85 or above

Close

Nutrienergy Balance (from the panel behind – rectify to 85 or above)
Close.

Nutrition and Homeopathy Feedback

QED TRIVECTOR SPINAL FOCI ANALYSIS - Nutrition Profile

Notes Report

119	Vitamin A Family	108	Protein Digestion
74	Vitamin B complex family	106	Carbohydrate Digestion
87	Vitamin G complex Family	76	Fat Digestion
60	Vitamin C Family	68	General Digestion
74	Vitamin D Family	99	Internal Enzymes
70	Vitamin E Family		Some Data on Risk menu
111	Vitamin F (fatty acids)	82	Minerals
72	Vitamin K Family	96	Amino Acids
81	Vitamin U Family (cq Q)	79	Oriental Herbal Combos
76	Biochemic Tissue salts		

Calculate, strike this key first

Go to Progressive Homeopathy list

Basic Homeopathy From Test Scan

Make master nutritional report

Access Reports + Texts

Good foods

Virtual Doctor

Insert Disease:

Male problems

Stress and Cortisol

Transfer To Test Screen

Intwined:

N/A

N/A

N/A

N/A

N/A

More Info On Natural Cooking on the NET

Constitutionals

Increase Enzyme Therapy

Stabilize Enzyme Therapy

Increase Hormone Therapy

Stabilize Hormone Therapy

Nutri/Enzy/Horm Stabilize

Close

From the **Main Matrix** screen:

Program (top menu bar)

Nutrition and Homeopathy Feedback (drop-down list)

Calculate Strike This Key First (values for various vitamins etc. appear)

OK in the pop-up box

DBL C is any specific Vitamin name to reveal individual items of that group in the degree of reactivity. Anything over 100 is significant.

Nutrienergy Balance

This will provide a training protocol for the whole group of items.

Rectify any choices to 85 or above

Stress and Cortisol

QED TRIVECTOR SPINAL FOCI ANALYSIS - Nutrition Profile

Notes Report

119 64 Adrenalin index value name

74 55 Prostaglandin 2689 99 Glutathione Peroxidase | Enzyme used in energy and immune function

87 22 Testosterone 8824 99 DEPOSITION DISEASE | Accumulation of lipofuscins in ageing cells we

60 113 Cortisol 587 99 D-XYLOSE | Wood, sugar, non fermentable pentose, undigestible by h

74 115 Histamine 9274 99 Proactazyme | Vegetarian Enzymes / Sunshine Tool Works

70 120 Pineal Pituitary Hypothalamus 2677 99 Cytochrome B | lds oxygen and energy disease. Enzyme.]

111 19 Emotional Allergy tendency 8966 100 Fish Oil and Co-enzyme Q10 (Co-Q10) /Omega3 fatty acids- BTO: avoi

72 112 Corticotropin Releasing Factor 7655 100 HAD COW DISEASE REMEDY - BOVINE SPONGIFORM ENCEPHALOPATHY

81 122 ACTH 2716 100 Phospholipase A | Enzyme for phospholipid regulation, lds phospholip

76 Dbl Click name to treat 2654 100 NeuroBooster Complex, with Brahmi extract and L-Tryptophan, Life C

Adrenal Balance 6106 100 REFSUM'S DISEASE | Phytanic acid storage disease see diet listed in d.

Make report Set Time min: 1 6078 100 CERVIX FEMALE DISEASE |]-"feminine"- listed in d. dictionary progra

2508 100 GASTRIN | FRAGMENTS | lds gastrointestinal disease, (bowel, colon, in

785 100 Organic Matcha, Life Care® | 7100 | - For anti inflammation enzyme ar

1171 100 PYRUVATE OXIDASE | Enzyme used in vitamin B utilization, lds regulat

Hormonal Imbalances = Mental Factors, Nicotine, Adrenalin, Pituitary, Mental Factors, Mental Factors,
Disease factors = Unmanaged emotional stress
Other Conflicts with - Guilt, Guilt,
Intwined:
Intwined with Circulation
Intwined with Allergy
Intwined with Bacteria
Intwined with Environment
Therapy: Herbal, Stress Reduction, Herbal,

Sort Table by: Value Name index

More Info On Natural Cooking on the NET Constitutionals

Increase Enzyme Therapy Stabilize Enzyme Therapy Increase Hormone Therapy Stabilize Hormone Therapy Nutri/Enzy/Horm Stabilize

Close ?

From the *Nutrition* panel:

1. *Stress and Cortisol*

2. *Do you wish to use the Questions? Yes or No*

Selecting *Yes* will provide a questionnaire that will fine-tune the training process.

3. *Adrenal Balance* – rectify to 85 or above

4. *Close*.

Chapter 13

Stimulate Emotional Growth Stress Reduction

From the *Main Matrix* screen:

Click *Programs* (top menu bar)

Click *NLP Emotional Growth Stress Reduction* (drop-down list)

(Rectify all chosen Training Protocols to 85 or above)


Click *Stimulate Emotional Growth* (stabilizes defense mechanisms)

Click *Therapies* (top menu bar)

Click *Unconscious choice of therapy* (bottom of drop-down list) – after half a minute or so, an item will appear in a blue bar at the bottom left of the screen (right below where Rectified: is).

Double Click the item posted in the blue bar.

Click *Close*.

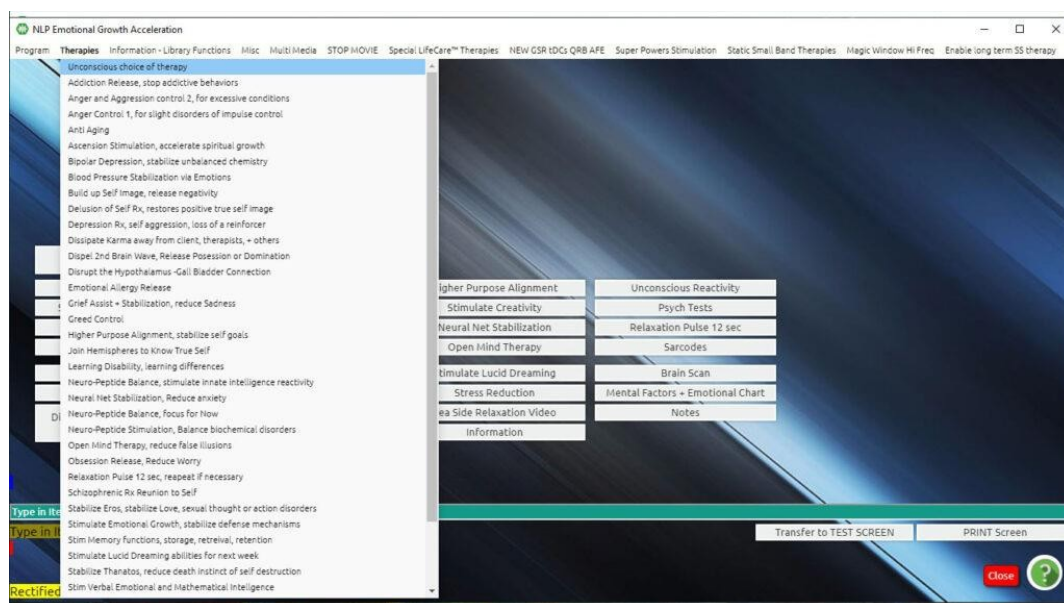

– Or –

Consider adding the following options:

Click Dissipate Karma

Caution: This is a very powerful selection and it is possible to expect that life may rapidly produce situations presenting opportunities to change old patterns of thought and behavior that are no longer conducive to well-being. A possibility would be to follow up with Ascension Stimulation and Higher Purpose Alignment.

Click Ascension Stimulation – for accelerating spiritual growth

Click Higher Purpose Alignment – for self goals

Click Stimulate Insight – for self-awareness

Click Stimulate Creativity – for inner potential

Click Neural – Net Stabilization – for anxiety reduction

Click Neuro-Peptide Stimulation – for biochemical disorders

Click Neuro-Peptide Balance – for focus and “get-up-and-go”!

Click Stim Memory – for memory functions of storage, retrieval, and retention

Click Stim Intelligence – for verbal, emotional, and mathematical intelligence

Click Open Mind Therapy – for reducing false illusions

Click Disrupt the Hypothalamus Gall Bladder Connection – for obsessive-compulsive behaviors

Click Close (when complete with choices).


Chapter 14


Other Biofeedback Choices

Any choice can be made from the drop-down list of items – watch the rectified value in the top menu bar.

Click **End Therapy** -This will end the program.

Close.

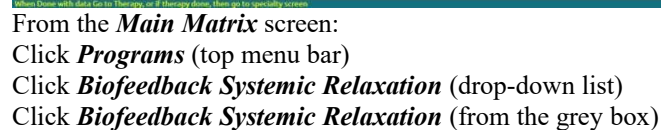

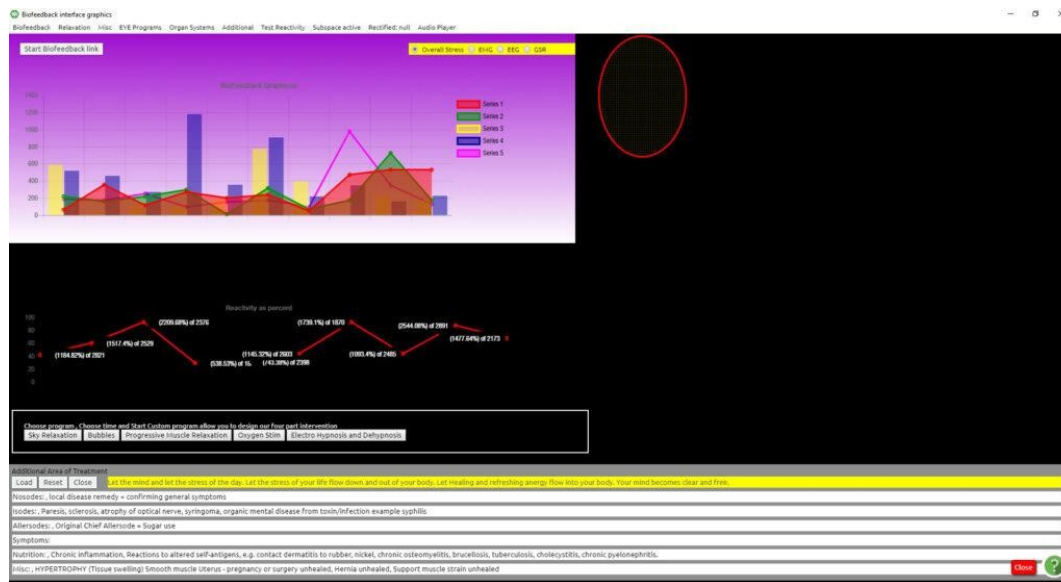

Click **Biofeedback** (top menu bar of the Biofeedback

Panel) Click **Reduce Stress in Nerves** (drop-down list)

Watch the “Rectified” value in the top menu bar on the right side. When the Training Protocol is complete, the rectification score will stop changing. For the Biofeedback program, it is important to realize that you must stop it yourself or it will continue indefinitely. The panel will continue to flash until you:

Click **End Reduce Stress in Nerves Therapy** – This will end the program.

Any choice can be made from the drop-down list of items – watch the rectified value in the top menu bar.

Click **Close**

Chapter 15


Body Scan and Face Therapy

From the **Main Matrix** screen:

Programs (top menu bar)

Body Scan and Face Therapy (drop-down list)

Scan All (top center of the screen) – let this complete

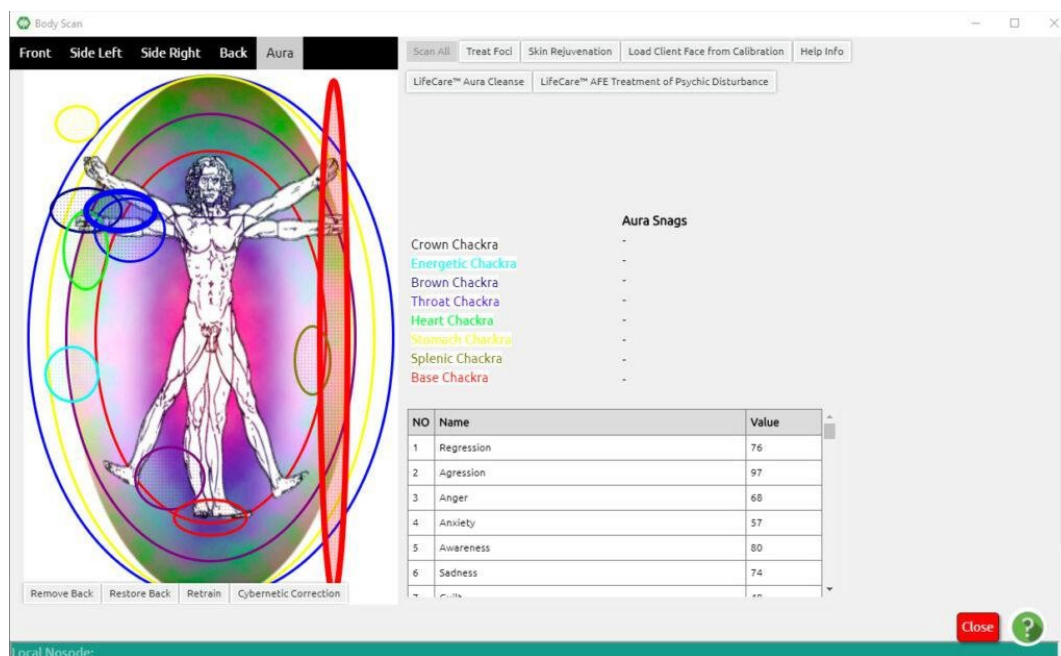

Front (upper left of the screen)

Treat Foci until there are no longer any red or purple (Foci detection) nets on the diagram (Do not mistake the Crown or Base Chakra for Foci as they are shown in the same color)

Side Left (upper left of the screen)

Treat Foci (as above)

Side Right

Treat Foci (as above)

Back

Appendix A

EDUCATIONAL LANGUAGE

Example #1:

A client reports: "I have a sore throat. What can I do about it?"

A *diagnostic* response would be: "Your sore throat is probably strep because that is what has been going around".

The practitioner is diagnosing what is wrong, by telling the client why she/he has a sore throat.

A *prescriptive* response would be: "Take vitamin C three times a day for a week, and drink slippery elm tea to soothe your throat and you will feel better".

An *educational* response would be: "When I feel like that, I take extra vitamin C, and drink tea made from slippery elm bark. If you look up and read material on these two substances you can learn about their relationship to building the immune system and how the bark of the slippery elm has been used to soothe mucous membranes throughout the body and ease the soreness of my throat. I find that it also tastes good".

Example #2:

A client reports: "My knee is really bothering me. It hurts every once in a while, ever since a high school football injury. Lately, it just hurts all the time. I'm tired of taking these painkillers constantly, I'm sure it can't be good for me. What else can I do?"

A diagnostic response would be: “You probably have arthritis in your knee”. Or, “Arthritis usually starts several years after an injury”.

A prescriptive response would be: “Try taking Glucosamine three times a day, to decrease the swelling in your knee, and the pain will subside”.

An educational response would be: “Some of my other clients have told me that some dietary supplements are reported to reduce the swelling and decrease the pain in their joints. You might want to read up on the effects of adding Glucosamine to a daily routine”.

Using educational language verbally throughout your session, and using only educational and “client” language throughout all of your written forms, notes, and materials in your record-keeping is important both legally and ethically. It is not illegal to inform and educate about wellness and health, yet the interpretation of how the information is conveyed can make the difference between being viewed by some as “practicing medicine without a license” or not. Educating and empowering the client is the ultimate goal of the practitioner. Practicing with high professional integrity and competency is imperative.

INCORRECT TERM	CORRECT TERM	COMMENTS
Cure	Re-educate and Relax	Cure: to correct, heal, mend, regenerate, rehabilitate, relieve, remediate, repair, restore, treat or otherwise help a person recover from any medical or psychological condition. It is illegal for any unlicensed professional to cure or treat any condition, disease, or disorder. It is illegal for any unlicensed professional to cure a medical or psychological condition, disease or disorder. Re-Educate: to educate again, as for new purposes; to re-develop or re-train. Relax: to make less tense, rigid, or firm, make lax; to relax the muscles; to release or bring relief from the effects of nervousness and tension, and anxiety; to become less intense.
Diagnose	Relaxation training programs manage the stress related to the reactivity of a tri-vector signature.	Diagnose: to analyze, appraise, assay, assess, calculate, evaluate, estimate, examine, experiment, gauge, guess, investigate, judge, probe, rank, rate or reckon any medical or psychological condition, disease or disorder. It is illegal for any unlicensed professional to diagnose a medical or psychological condition, disease or disorder. Relax: to make less tense, rigid, or firm, make lax; to relax the muscles; to release or bring relief from the effects of nervousness and tension, and anxiety; to become less intense. Biofeedback is a practice that enables an individual to learn how to relax and change some physiological activities for the purpose of improving health.

INCORRECT TERM	CORRECT TERM	COMMENTS
Diseased	Stress	Disease: a disordered or incorrectly functioning organ, part, structure, or system of the body resulting from the effect of genetic or developmental errors, infection, poisons, nutritional deficiency or imbalance, toxicity, or unfavorable environmental factors; illness; sickness; ailment. Only licensed healthcare professionals can treat a disease. Stress: a biological reaction to physical, emotional, mental, psychological, or spiritual stimuli, whereby normal functioning is compromised.
Patient	Client	Patient: the term “Patient” is applied to clients under the care of a licensed healthcare professional. Certified Biofeedback Technicians and Specialists use the term “Client”. Client: a customer; a person that engages the professional services of another.
Prescribe	Educate and	Prescribe: to advise, apply therapies, coerce, counsel, manipulate, promote,

	Coach	recommend, require, suggest, or otherwise treat a client. Educate: to train by formal instruction and supervised practice; to provide information. Coach: one who instructs trains. As Certified Biofeedback Technicians and Specialists, we educate and coach the clients to increase their awareness.
--	-------	---

INCORRECT TERM	CORRECT TERM	COMMENTS
Test	Electro Dermal Response Stress Reaction	Test: a particular process or method for trying and assessing; to examine or evaluate. Electro Dermal Response (EDR): a technique by which sensors monitor the skin's electrical resistance. Biofeedback devices measure the voltage, current, and resistance of the skin. Certified Biofeedback Technicians and Specialists are trained to use biofeedback equipment to measure the Electro Dermal Response stress reactions and to teach their clients how to relax and re-educate their muscles, so they can manage their stress and pain.
Therapies	Training Programs	Therapies: the treatment of disease or disorders, as by some remedial, rehabilitating, or curative process; a curative power or quality. Training Programs: for the purpose of correct verbiage, we will refer to the biofeedback programs as training programs Trained practitioners can teach clients how to relax, and re-educate their muscles in order to manage their stress and pain.
Therapist	Certified Biofeedback Technician or Specialist	Therapist: one who specializes in the provision of a particular therapy; one specializing in therapy, especially a person trained in methods of treatment. Certified Biofeedback Technician or Specialist: a person trained to use biofeedback equipment to measure the biological activities of their clients, in order to teach their clients how to relax and re-educate their muscles, so they can manage their stress and pain.

INCORRECT TERM	CORRECT TERM	COMMENTS
Therapy	Relaxation Training	Therapy: the treatment of disease or disorders, as by some remedial, rehabilitating, or curative process; a curative power of quality. Relaxation Training: to help make less tense, rigid, or firm, make lax; to help relax the muscles; to help release or bring relief from the effects of nervousness and tension, and anxiety; to help become less intense. Replace the word therapy, treat, or treatment with Train or Training Program. Biofeedback is the process of helping clients measure and changes their biological activities by relaxing and re-educating their muscles in order to help them manage their stress and pain.
Treat	Train	Treat: to deal with a disease or a patient in order to relieve or cure; to care for or deal with medically or surgically. Train: to form by instruction or discipline; to teach so as to make fit, qualified, or proficient; to make prepared (as by exercise) for a skill. Biofeedback training programs use electrical impulses to reduce Electro Dermal Response stress reactions. To treat is to attempt to cure by various means. In most jurisdictions, it is illegal for any unlicensed professionals to treat any medical or psychological condition, disease, or disorder.

TERMINOLOGY for PRACTITIONERS

The ED.X and NUCLEUS System utilizes five comfortable, conductive straps that envelop the wrists, ankles, and forehead. Once the Calibration process is complete, a three-minute Electro Dermal Response stress reaction test is performed with results displayed on the ED.X and NUCLEUS computer software interface.

Through the use of biofeedback training programs, Certified Biofeedback Technicians and Specialists can train their clients with stress reduction training programs for muscle re-education, relaxation training, pain management, and brainwave training to improve the quality of their lives.

The programs are “relaxation training for tri vector signature”. Trivector signatures appearing on the program should not be discussed with clients.

Certified Biofeedback Technicians and Certified Biofeedback Specialists and Unlicensed Practitioners are trained to be professional providers but are not qualified to diagnose, treat, cure or prevent any medical or psychological condition, disease, or disorder. By law, it is illegal in the US, Canada, and most other jurisdictions, it is legal for qualified professionals to practice Biofeedback as long as they do not claim to diagnose, treat, cure or prevent any medical or psychological disorder, disease, or condition. In addition, they may not use any terms that imply any of the aforementioned.

The following are guidelines for practicing with safe terminology:

DO NOT Assess	DO NOT Make Claims or Promises
DO NOT Counsel	DO NOT Prescribe
DO NOT Cure	DO NOT Recommend
DO NOT Detect	DO NOT Repair
DO NOT Diagnose	DO NOT Scan
DO NOT Evaluate	DO NOT Test
DO NOT Heal	DO NOT Treat
DO NOT Identify	DO NOT Use the word “Disease”

- The items mentioned above are not within the scope of biofeedback.
- We use biofeedback stress reduction training programs for muscle re-education, relaxation training, pain management, and brainwave training.
- We explain that health is an individual matter with many approaches to healing and no guarantees.
- Biofeedback is a practice that enables an individual to learn how to change some physiological activities for the purpose of improving health.
- Biofeedback instrumentation uses mild electrical impulses that measure skin temperature known as (EDR) Electro Dermal Response. Electro Dermal response has replaced the term (GSR) Galvanic Skin Response

Appendix B

Electrode Data Check Box Does Not Appear

The **Electrode Data Screen** box should appear at the **25%** indicator when the program is loading. (*This only takes place if you have the ED.X and NUCLEUS device connected to your computer*). Make sure the device is turned on and connected before you open the **ED.X and NUCLEUS** program from your desktop. You will see the circle around the power button light up, and then hear a musical tone as it initializes.

If you have a preloaded laptop that came with your purchase it is usually just a matter of having it connected to the designated USB port. A document will come with your computer explaining which port has been designated. The port is usually on the left side bottom or back. It may also be one on the back of your computer, generally the left side, depending on the Computer Model. The right side port is usually reserved for an external mouse.

BASIC TROUBLESHOOTING

MOST FREQUENT CAUSES OF FREEZE-UP:

Freezes occur much more frequently with a new practitioner than with a seasoned one. Why might that be? It's because they may tend to click indiscriminately at buttons, faster than the computer can handle. It boils down to a lack of knowledge about when to click, which buttons to click, how many clicks (single or double click), at what interval of clicking as well as what should happen after a button is clicked! Each click is essentially an instruction to the computer, and too many instructions coming at once can put the computer on "brain overload". If you double-click when all that is needed is a single click, that ends up being a double command and the computer tries to process both commands. If you add to that, another hasty click because it is not responding, what could result is a white screen and freeze-up.

The computer needs to take its time to process those instructions in the background, one by one until it is ready to proceed again.

If one waits **P-A-T-I-E-N-T-L-Y** and gives the computer time to catch up, it is then able to process further commands. Unfortunately, the freeze-up is often not just temporary and other measures must be taken.

Another common problem is that the touchpad on your computer can be very sensitive and can mistake a brush of the finger for a click. *A highly advisable solution would be to get an external USB or wireless mouse.* This can save you a lot of frustration and make it easier for you to navigate the software.

The experienced **ED.X and NUCLEUS** user has learned what can be called the Art or Zen of Clicking! Through hard-won experience, the seasoned practitioner knows all that is needed is a single click or a double click. The practitioner, in essence, learns to dance with the system. Getting to know the timing and sequencing of when and where to click becomes as important as learning the steps of the tango with a partner. As a newcomer, you will need to give yourself time to experiment – to discover – which buttons require a waiting time before the next button is clicked, which program needs to be clicked on before the other, etc.

TRAINING PROGRAMS PRONE TO LOCK UP

- Disease Probability and Dictionary
- Biofeedback Systemic Relaxation Training Programs
- Auto Trivector
- Auto Scalar
- Spinal Muscular Re-Education, Sarcode Feedback
- ED.X and NUCLEUS or AFE Training Programs
- And Other ED.X and NUCLEUS Training Programs

Sometimes a **WHITE SCREEN** appears in place of the panel image you were previously working on. Do not click anything! The program is attempting to process the information and is overloaded. Just be patient and wait a few

minutes. The panel image will return once the process that was engaged has been completed. If you are waiting and this is taking much longer than is necessary (a few minutes or the amount of time the program was set for). This can often be avoided if you are careful not to click a second time, too hurriedly, on a process that is already running or one that is not immediately responding. Please see *What to do if the Screen is Frozen*

DISEASE PROBABILITY AND DICTIONARY TRAINING

It is best to activate the Disease Probability and Dictionary training program before activating the Bodyviewer and Irid. Avoid using the Long Term Therapy choice if you want to be free to change to other training choices during your session.

Do not use the Minimize button at the top right of the panel. This goes for all panels including Iridology, Bodyviewer, Consciousness, Disease Dictionary, Irid, Bodyviewer, etc. (to toggle between all programs). Or you can press the Alt + Tab keys on your keyboard to easily rotate between programs.

When you have started the Disease Probability and Dictionary it is best to wait a few moments for it to begin the process before clicking on Consciousness (bottom task bar) to return to the **ED.X and NUCLEUS** Program. If you try to return to **ED.X and NUCLEUS** too soon, it can lock up. Just wait a few moments and try to click on Consciousness again. This should minimize the Disease Dictionary training program to the background while you continue to work in another panel.

BIOFEEDBACK SYSTEMIC RELAXATION TRAINING PROGRAM

Problems With Graphics or Multimedia

If the *Bodyviewer* is open and activated properly, video graphics will display in the middle of the *Bodyviewer* panel. Not all items have a video graphic, so in some cases, you may get an Unspecified Error message. In this situation, click Stop Animation in the top menu bar to remove the error message. To hack to see if the graphics are playing properly you can select any picture choice to the left of the panel, once you click on it, the video graphic should start to play.

Make sure *DivX* is installed on your computer

- Download the DIVX codec by clicking [HERE](#)
- Install the downloaded codec and you are good to go.

If the codec pack is set properly there should be no problems with any graphic or multi-media options. That means in Manfly.avi, Bodyviewer, Iridology, Body Scan and Face Therapy, or other ED.X and NUCLEUS Training Programs.

This codec is already included in your windows 10 installation, so there is no need to download it, the videos should run fine.

Biofeedback Systemic Relaxation

Sometimes it seems like the *Biofeedback Systemic Relaxation* training program gets frozen. The following processes may help alleviate problems with potential lock-ups in the Biofeedback Training Program:

- Activate the *Reduce Stress in Nerves* (drop-down menu) under the top menu bar of *Biofeedback Systemic Relaxation* (top menu bar). Do this before activating other training choices. This particular process seems to tune to the client quickly and you will see rectified numbers flickering and changing almost immediately after activating (in the top menu bar). This seems to lessen the chances of freezing in the rest of the Biofeedback Systemic Relaxation training programs. Rectification occurs when you have a value above 85

that stops changing. Then you can click the Click Here to Stop button (in the top menu bar) to stop the process.

- When you activate the next training protocol watch to see that rectified values are appearing and changing. If you see nothing within the first few seconds Click Here to Stop and try to re-activate the process. Sometimes the freeze-up occurs if it does not tune to the client quickly enough.
- If you are having trouble closing a Training Program and there is a smaller three-dimensional graphic box in the upper left-hand corner of your screen, you can try to close it first. Use the red X in the upper right-hand corner of this smaller graphic box and then attempt to Click Here to Stop and terminate the process.
- If the ED.X and NUCLEUS Working box is pulsing on and off during a particular Biofeedback Systemic Relaxation training program; there is a special technique to stop it. Click on the Click Here to Stop. If that doesn't work, you can click on the Click Here to Stop in the menu bar very quickly 20-30 times. An easy way to do this is to hold the mouse with one hand and click with the index finger of the other hand. One of these clicks will produce just the right timing to catch it on the color change (when the Click here to Stop button pulses darker).
- Another button to try: Stop Therapy Loop might also work. However, when you click it tells you to *Please Wait*, you still might have to click the *Close* button on the **ED.X and NUCLEUS is Working** panel. You might have to do this a few times. If the Close button gets a line around it, it means it is frozen.
- You can press the enter key (sometimes there is an open box in the background that you can't see, and it is waiting for a response).

Disease Probability And Dictionary

It is best to activate the Disease Probability and Dictionary training program before activating the Bodyviewer and Irid. Avoid using the Long Term Therapy choice if you want to be free to change to other training choices during your session.

Do not use the Minimize button at the top right of the panel. This goes for all panels including Iridology, Bodyviewer, Consciousness, Disease Dictionary, Irid, Bodyviewer, etc. (to toggle between all programs). Or you can press the Alt + Tab keys on your keyboard to easily rotate between programs.

When you have started the Disease Probability and Dictionary it is best to wait a few moments for it to begin the process before clicking on Consciousness (bottom task bar) to return to the **ED.X and NUCLEUS** Program. If you try to return to **ED.X and NUCLEUS** too soon, it can lock up. Just wait a few moments and try to click on Consciousness again. This should minimize the Disease Dictionary training program to the background while you continue to work in another panel.

Hints And Warnings

If you had to close **ED.X and NUCLEUS** using *the* Ctrl+Alt+Delete method when you restart a warning message will appear ED.X and NUCLEUS were shut down improperly do you want to restore the last client, **Yes** or **No**

Hint:

- If you are returning to the **ED.X and NUCLEUS** program after a lock up and you used Ctrl Alt Delete or another method to shut down, you will usually get a message ED.X and NUCLEUS were shut down improperly do you want to restore the last client?
- If you reload immediately with the same client and are trying to return to ED.X and NUCLEUS after a lock-up; click **Yes** to the Question.
- After Password, go directly to “**Test**”, from the Main Menu and resume. You would skip the Demographic and Calibration steps. This way the last client you were working with will be restored and appear as the **Client** on the **Main Matrix** screen.

Warning:

- If you clicked **Yes** but intend to load a different client you will end up with overlapping information. The last client’s info gets restored to the blue info box, but the new client is loaded as the **Testee**.
- Always click **No** to this message when your intention is to load a **new or different client** since the last session.
- If you intend to reload the **same client** and are trying to return to ED.X and NUCLEUS after a lock-up; click Yes to the Question.
- If, when you open the Password screen, there is a message asking you to enter your name or Country, **DO NOT DO IT**, as this could result in a loss or corrupt your client files. (The only time you should need to re-enter the data is after loading a new version or re-installing the software).
- If you are asked to insert your name or Country, just get out of the program using the Ctrl Alt Delete method to close the program.
- Shut down the Computer and turn off the ED.X and NUCLEUS devices.
- Wait two minutes and restart the computer.
- Restart. Reload the client as usual but do not retest. The test results will be restored. However, you will notice that the color values change and there will be restoration. However, you will notice that the color values change and there will appear to be more reds, purples, and yellows. This happens when it reverts to a default setting. Instead, you will have to pay more attention to the number value and placement of the items.

Setting The Com Port

If it is still not detected you may need to check and set the Com Port:

ATTENTION: ED.X and NUCLEUS ONLY work ON COM4!

Take note some of these steps require a right mouse click vs. the usual left mouse click. Keep the device turned off during this process.

- **Right** Click My Computer
- Click Properties
- Click Hardware Tab
- Click Device Manager
- Find Ports and Double Click to Open it may say Ports (Com & LPT)
- Click Find Prolific USB to Serial Com Port#
- (if it is Com4 this is correct, if not continue)
- Right Click on Prolific USB to Serial Port Com#
- Click Properties

- Click Advanced
- Click on the Com # to highlight it
- Use the up arrow on your keyboard and click up until you get Com4
- Click OK, and close
- Click OK again as you close out
- Close all windows and go back to the Desktop
- Turn the device on again and reopen the **ED.X and NUCLEUS** Program

Special ED.X and NUCLEUS Programs Or AFE

Several panels in the program have additional **ED.X and NUCLEUS Programs** (any button that has **ED.X and NUCLEUS or AFE** in the name). These are found on:

- **NLP** screen (Special **ED.X and NUCLEUS** Therapies)
- **Body and Face Scan**, [**ED.X and NUCLEUS** Aura Cleanse, **ED.X and NUCLEUS AFE** Treatment (Training) or Psychic Disturbances]
- **Anti Aging** (**ED.X and NUCLEUS** Interactive Anti-Aging pulse)
- **Check Current Rectifications** (Super **ED.X and NUCLEUS** Electrical Rectification)
- **Dark Field Analysis** (**ED.X and NUCLEUS** Pulse for Immune Stim), and **Blood Cleanse** (**ED.X and NUCLEUS** Pulse for Stem Cell Stimulation)

If you click on any of these **ED.X and NUCLEUS** training protocols, the **ED.X and NUCLEUS are Working** box comes up and does not go away. **You are not frozen!**

Using the AFE or Special ED.X and NUCLEUS Training Programs

- Once you have selected one of these Training Programs, you must pull the ED.X and NUCLEUS is Working box up by clicking and dragging on the top border of this panel. When the whole box is visible on your screen you will see a white area on the right bottom inside of the ED.X and NUCLEUS is a **Working** box. (The **Aura Cleanse** is an exception, this one only has a button to “**Stop Aura Cleanse**”, do not use the **Stop Therapy loop** button)
- Click the **Start** (name of the program) in the white area. The rectification will show to the right of these keys. This can take several minutes.
- Click on the Stop/End button (for the name of the program) in the white area.
- Click OK (on the pop-up message box)
- Click **Close** on the **ED.X and NUCLEUS is Working** box to remove it from the screen.
- **Click OK (if you get another pop-up message)**

HINT:

If you get into a situation where you are caught in a loop and the pop-up keeps repeating after clicking close, the following is a little trick you can try:

- Drag the pop-up message so that the OK is aligned with the Close button.

- Now place the cursor on the OK use one hand to do the mouse click and the other to quickly click Enter on your keyboard immediately after. This may take several attempts.

Training Programs Prone To Lock Up

- *Disease Probability and Dictionary*
- *Biofeedback Systemic Relaxation Training Programs*
- *Auto Trivector*
- *Auto Scalar*
- *Spinal Muscular Re-Education, Sarcode Feedback*
- *ED.X and NUCLEUS or AFE Training Programs*
- *And Other ED.X and NUCLEUS Training Programs*

Sometimes a **WHITE SCREEN** appears in place of the panel image you were previously working on. Do not click anything! The program is attempting to process the information and is overloaded. Just be patient and wait a few minutes. The panel image will return once the process that was engaged has been completed. If you are waiting and this is taking much longer than is necessary (a few minutes or the amount of time the program was set for).

This can often be avoided if you are careful not to click a second time, too hurriedly, on a process that is already running or one that is not immediately responding.

Please see *What to do if the Screen is Frozen*.

Tricks To Prevent And Deal With Freeze-Up

Be patient and give yourself (and the computer) **TIME**. The **TIMING** of your clicking is important.

- Completely shutting down both the **ED.X and NUCLEUS programs** and the computer between clients for at least two minutes clears the **RAM**.
- Consider running your computer maintenance functions.
- Use each screen's **Close** button as opposed to the red **X** at the top right of the window.
- When you are running the **ED.X and NUCLEUS** programs. Avoid **stacking** panels as much as possible. It is better to return to the **Main Menu** or **Main Matrix** screen before venturing off in a new direction. Bill Nelson has provided various links from panel to panel within the program, but using these links can add unnecessary complexity to the issue and will sometimes increase lockups. (Example: If you are in the Risk Profile training program and link to Spinal, yes you can check for vertebrae relationships and train them, but it may not be advisable to start going into other spinal training programs, or linking from there, to additional training program panels.)

Using Ctrl+Alt+Delete Tricks

- You must press the **Ctrl** Alt and Delete keys on your keyboard all at once and release.
- The Windows Task Manager box appears.
- Make sure you are in the Applications tab.

- The programs that you have opened will be listed: Consciousness, Disease, Bodyviewer, Manfly.avi, and/or Irid.
- In the Status column, check to see whether the programs are running or not responding. If the program is listed as running, that is good – you are more likely to be successful in getting out of the freeze-up.
- At the bottom of the box, click on New Task.
- The Create New Task box now appears. Close by clicking the red X on the top right-hand corner of this panel.
- Close by clicking the red X on the top right-hand corner of the Task Manager box.
- All the above steps may seem complicated at first but once you get the hang of it, there's really nothing to it. Doing this seems to trick the computer into coaxing a process to complete or for the OK message to appear on top, rather than behind the panel.
- Another variation of the above trick is to use the Switch To button (instead of New Task). If Consciousness is highlighted, change the selection by *clicking on the Bodyviewer, now click Switch To. Attempt this a few times* back and forth and it may unfreeze the computer.

What to do when computer screen is frozen

- Hold the **Alt** and click the **Tab** button on your keyboard. View the panels available by keeping one finger on **Alt** and then clicking and releasing the **Tab** you can move between the different panels that are open on your computer. Going through one or two rotations like this could un-stick you, but it only works for simple lock-ups.
- In some cases, a panel will freeze because a pop-up message box is waiting for an OK click, but it is stuck behind the panel. The most often occurs in The Spinal and Sarcode, and Risks Profile training programs and other programs linked to it (Ex: pain, nerves, muscles, sports, etc.) To resolve this simply put your cursor towards the bottom of your computer screen to bring up the taskbar so you can click on Consciousness. Clicking on this will often bring the message box into the foreground, so you can click **OK** and continue.
- Another option is to press the **Alt** and the **F4** keys together simultaneously. *Only hold them down briefly*, as a long hold can close you right out of the whole program. Attempt this a few times to see if it works.
- If you're running a training program with the **ED.X and NUCLEUS is Working** box on your screen you can try this option as a last resort:

Click and drag on the top border of the ED.X and NUCLEUS is Working box to see if you have the option to Stop Therapy Loop. When you click on it, you will get a message Please Wait. Be patient – it may take several minutes, or it may not work at all. It is worth trying if all other options fall.

Updates and chnages

Brain wave pattern therapy

Basic values are obtained by pressing the “TREAT” button

Electro-Physiological Oscillation

Device Diagnostics

Organ Electrical Patterns

ECG Fast Fourier Transform

Trivector EMG Patterns

Autonomic Nerval Systems

Brain Wave Pattern Therapy

Subspace Morphic Trans

000 Alpha, treats stress

000 Beta-Gamma, stabilizes

000 Theta, Consciousness Stim

000 Delta, For insomnia

000 Immature, slow thought

000 Normal

000 Cerebellar Dysfunc

000 Local Slow Waves

000 Amplitude Disorder

000 Epileptic Forms

000 General Asymmetry

000 Injury Profile

000 Bilateral Sync

000 Deep Brain Asymmetry

000 Personality Dysfunc

000 Attention Deficit

Treat clicked items

Stim Brain with Wave Form Clicked, 2 min.

Attention Deficit, brain Fatigue, LEAP

Average frequency:

0

Text on EEG_ECG

Music Therapy

Anti - Aging

START

Results

Notes

Fourier + Hysteresis Fractal Fuzzy Analysis NFP-HS

0%

Close

?

Subspace Active

Balancing begins with “Delta Waves” from the left column. The highest value in the right column will be automatically selected

Click on “Treat clicked items”

Electro-Physiological Oscillation

Device Diagnostics

Organ Electrical Patterns

ECG Fast Fourier Transform

Trivector EMG Patterns

Autonomic Nerval Systems

Brain Wave Pattern Therapy

Subspace Morphic Trans

95 ● Alpha, treats stress

61 ● Beta-Gamma, stabilizes

112 ● Theta, Consciousness Stim

113 ● Delta, For insomnia

113 ● Immature, slow thought

83 ● Normal

70 ● Cerebellar Dysfunc

100 ● Local Slow Waves

104 ● Amplitude Disorder

69 ● Epileptic Forms

68 ● General Asymmetry

69 ● Injury Profile

102 ● Bilateral Sync

69 ● Deep Brain Asymmetry

65 ● Personality Dysfunc

95 ● Attention Deficit

Treat clicked items

Stim Brain with Wave Form Clicked, 2 min.

Store Brain Wave Engram to Prayer

Attention Deficit, brain Fatigue, LEAP

Average frequency:

24569

Text on EEG_ECG

Music Therapy

Anti - Aging

START

Results

Notes

Fourier + Hysteresis Fractal Fuzzy Analysis NFP-HS

0%

Past Trauma Pattern Detected

Worry Pattern Detected

Obsessive Pattern Detected

Past Trauma Pattern Detected

Vitamin Deficiency Pattern Detected

Depression Pattern Detected

Mineral Deficiency Pattern Detected

Worry Pattern Detected

Sadness Pattern Detected

Make Report of Medical Reference

Close

?

Subspace Active

After the value selected in the right column drops by 10 or more below the level of the “NORMAL” value, we move to “Theta waves”

The procedure is the same, the next highest value in the right column is selected and treatment is performed until the selected value drops by 10 or more below the “NORMAL” value, by clicking on “Treat clicked items”

Also balance the cell waves (“Beta Gamma”, “Alpha”)

Electro-Physiological Oscillation

Device Diagnostics

Organ Electrical Patterns

ECG Fast Fourier Transform

Trivector EMG Patterns

Autonomic Nerval Systems

Brain Wave Pattern Therapy

Subspace Morphic Trans

95 ● Alpha, treats stress

61 ● Beta-Gamma, stabilizes

112 ● Theta, Consciousness Stim

113 ● Delta, For insomnia

113 ● Immature, slow thought

83 ● Normal

70 ● Cerebellar Dysfunc

100 ● Local Slow Waves

89 ● Amplitude Disorder

69 ● Epileptic Forms

68 ● General Asymmetry

69 ● Injury Profile

102 ● Bilateral Sync

69 ● Deep Brain Asymmetry

65 ● Personality Dysfunc

95 ● Attention Deficit

Treat clicked items

Stim Brain with Wave Form Clicked, 2 min.

Store Brain Wave Engram to Prayer

Attention Deficit, brain Fatigue, LEAP

Average frequency:

24569

Text on EEG_ECG

Music Therapy

Anti - Aging

START

Results

Notes

Fourier + Hysteresis Fractal Fuzzy Analysis NFP-HS

0%

Past Trauma Pattern Detected

Worry Pattern Detected

Obsessive Pattern Detected

Past Trauma Pattern Detected

Vitamin Deficiency Pattern Detected

Depression Pattern Detected

Mineral Deficiency Pattern Detected

Worry Pattern Detected

Sadness Pattern Detected

Make Report of Medical Reference

Subspace Active

Close ?


Careful! The “NORMAL” value in the right column cannot be selected because it is the reference value for the other values.

[Changes made with Update no.43 & 44](#)

From the latest update made this last week, Update no. 43&44, some basic changes have been made.



As soon as you click the close button on this panel where the device is found, you shall enter the password panel directly, without having to choose the password option on the main panel.

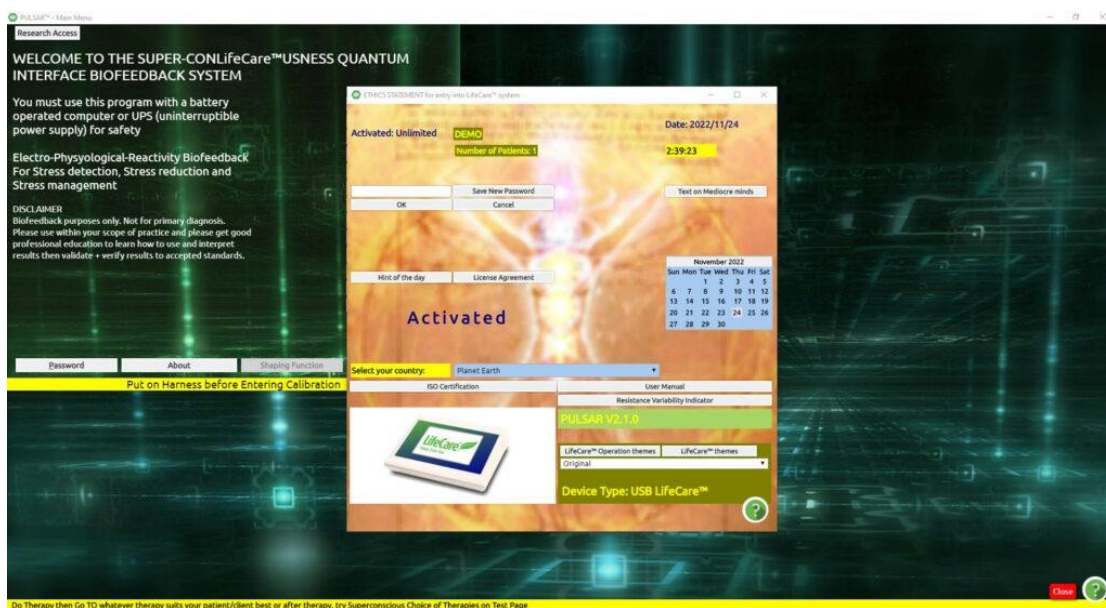

Do Therapy then Go TO whatever therapy suits your patient/client best or after therapy, try Superconscious Choice of Therapies on Test Page

As soon as you put in your password and close, the Demographics panel opens directly.

FUSART™ Universal Biofeedback System

Programs: Additional: Subacute Programs: Reporting: Auto Focus Zap: Retest: Order of Remedies: BPR Enhancements: Library: System Power Settings: Audio Player: Help

John Doe **V** **A** **W** **S** **U** **R** **E** **N** **D** **E** **R** **E** **M** **O**

regulatory dysfunction in organ system
 explore risk of degeneration
 explore risk of reactivity dysfunction
 reactivity: normal
 cellular vitality index: normal, below 3 chronic

Time remaining:
 Last visit: 62
 Next visit: 65
 Resistance speed index: 4
 Data integrity: 99
 Chase angle: 2
 Push response resistance: 1025
 Impedance: 1751

Body Viewer
 Degeneration Scan
 Disease Dictionary (Library)
 Biology (Library)
 Personal Evolution
 Organ Generator
 Quantum Prayer
 Supernatural Choice of Therapies
 Self-Defence

Expand Name **Prepare**

Data for John Doe **Soc Index: 126**
 Owner: DEMO
 Regulatory dysfunction in organ system
 Explore risk of degeneration
 Explore risk of reactivity dysfunction
 Reactivity: normal
 Return to Main Matrix
 GH Foods
 Guna Remedies
 HerbaPlus
 Hormones
 Horny Enzymes
 Imponderables
 Life Care™

Sys POW: Normal
 Percent Done
 0%
 09:00
 72%
 71%

Search for item:
 Search Above Tray
 Clear Search
 Close

Situation Hold:
 Search for item:
 Search Above Tray
 Clear Search
 Close

Harmonic Coupling of 2 items:
 Empty below

When Done with Data Gets Therapy, or if therapy done, then go to specialty screen

For now, let us proceed to test and we shall go through each change thereafter.

FUSART™ Universal Biofeedback System

Programs: Additional: Subacute Programs: Reporting: Auto Focus Zap: Retest: Order of Remedies: BPR Enhancements: Library: System Power Settings: Audio Player: Help

John Doe **V** **A** **W** **S** **U** **R** **E** **N** **D** **E** **R** **E** **M** **O**

regulatory dysfunction in organ system
 explore risk of degeneration
 explore risk of reactivity dysfunction
 reactivity: normal
 cellular vitality index: normal, below 3 chronic

Time remaining:
 Last visit: 62
 Next visit: 65
 Resistance speed index: 4
 Data integrity: 99
 Chase angle: 2
 Push response resistance: 1025
 Impedance: 1751

Body Viewer
 Degeneration Scan
 Disease Dictionary (Library)
 Biology (Library)
 Personal Evolution
 Organ Generator
 Quantum Prayer
 Supernatural Choice of Therapies
 Self-Defence

Expand Name **Prepare**

Data for John Doe **Soc Index: 126**
 Owner: DEMO
 Regulatory dysfunction in organ system
 Explore risk of degeneration
 Explore risk of reactivity dysfunction
 Reactivity: normal
 Return to Main Matrix
 GH Foods
 Guna Remedies
 HerbaPlus
 Hormones
 Horny Enzymes
 Imponderables
 Life Care™

Sys POW: Normal
 Percent Done
 0%
 09:00
 72%
 71%

Search for item:
 Search Above Tray
 Clear Search
 Close

Situation Hold:
 Search for item:
 Search Above Tray
 Clear Search
 Close


Harmonic Coupling of 2 items:
 Empty below

When Done with Data Gets Therapy, or if therapy done, then go to specialty screen

After Testing, you shall notice now a new order of the drop-down menus. Firstly, in red, you shall see the drop-down menu NEW.



These are the new programs that have been added so that you can find them and refer to them faster with one glance. Of course, this is their temporary place just in order to inform you of a new program added. If you select one of the programs in this NEW drop-down menu, let's say the first one for example you will notice a message in red coming up below.

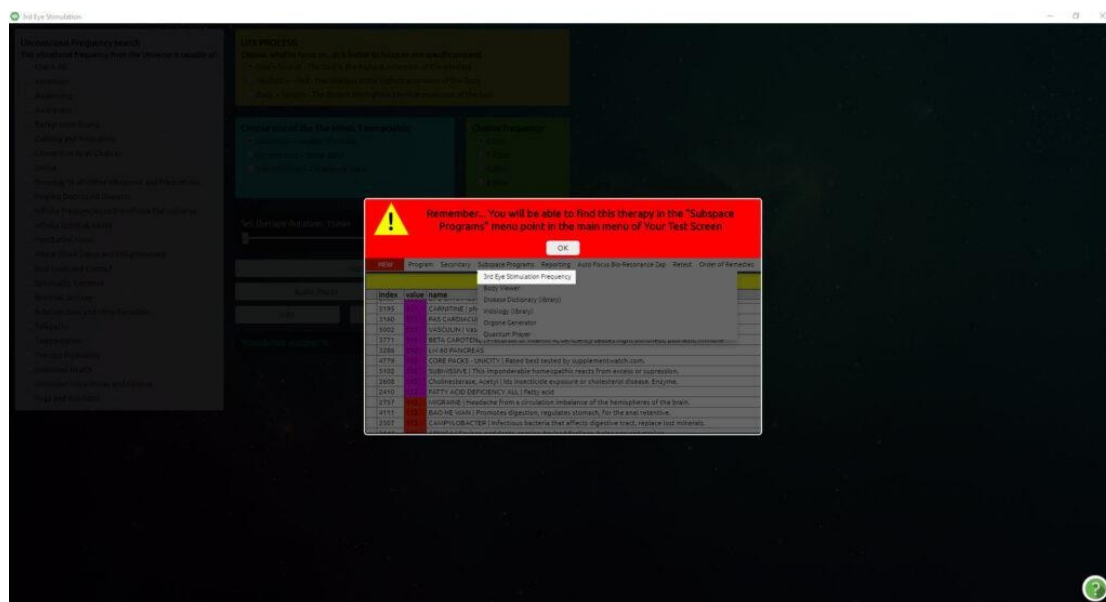

This informs us of where we will be able to locate it within the program. For instance, in this example, the 3rd Eye Stimulation Frequency, tells us that we will be able to find it under the dropdown menu Subspace programs. By clicking OK, we proceed in running this program. The same process is for each of the items under the NEW drop-down button. Remember: the NEW drop-down menu is only a temporary list of new applications added to our software. So please take note of the messages which appear, as, at some point later on in updates, you will find these items only in the specified areas.

The next drop-down menu which has changed is Library.

[illegible]

Here you can find library information, literature, things to read or information you may want to refer to, etc.

[illegible]

Enjoy your new update. This year so many new add-ons will appear.

With this update, several changes have been made.

AAS 344™ Universal Biofeedback System

New Programs Additional Subspace Programs Reporting Audio Video Display Order of Remedies EPR Influences Library System Power Settings Audio Player Help

V	A	R	H	O
79	61	32	90	90
Last visit:				
100	61	61	96	78

SOC index: 100
Vib rate: 32
Time remaining:
Pulse pressure: 60 Last visit: 63
Suction pressure: 74 Last visit: 75
Resonant frequency system: Relative accuracy:
#928 #94
Resistance speed index: Data integrity (%)
6 90
Phase angle: 4
Phase response reactance: 1142
Impedance: 1769
Close



Data for John Doe
Owner: DEMO
Regulatory dysfunction in spinal energy flow
Explore risk of acupuncture transfer dysfunction
Explore risk of energy management dysfunction
Reactivity:
Return to Hahn Matrix
Allerxes:
Amino Acids
Ayurvedic
Back Flowers:
Bacteria
Bio War
Blood

Soc Index: 100
SYS POW: Normal
Percent Done
0%
09:00
Reactivity / Rectify
0%
Shaping instn als
0%

Mold:
Total Mold Reaction Test Items in Mold Trap above Both Traps
Optimize Optimal Test items as Mold Trap above

Search file item:
Search Above Trap Clear Search Close

Some features are inactive.
Please run "Prepare" and "TEST" to activate them.



For more in-depth explanations of the readings please refer to our video VARHOPE Score Interpretation on our support site <https://support.quantum-medical.eu/> or Chapter 2 in our user manual.

If you wish to refer to the red VARHOPE panel at some other point during your therapy, simply refer to the drop-down menu Retest, and click on Activate red panel with vital organs.

The next new and exciting feature of our software is the Coronavirus Immunization panel. As you already know from our previous updates, all new programs added to our software are initially and temporarily found under the NEW red drop-down menu, for quick reference.

Universal Biofeedback System

Home Programs Additional Subspace Programs Reporting Auto Focus Zap Retreat Order of Remedies BPR Enhancements Library System Power Settings Audio Player Help

Coronavirus Immunization

Remember... You will be able to find this therapy in the "Additional" menu point in the main menu of your Test Screen

RETEST

Data for John Doe

Soc Index: 109

Owner: DEMO

SYS POW: Normal

Percent Done

Regulatory dysfunction in signal energy flow

Explore risk of oxygenation transfer dysfunction

Explore risk of energy management dysfunction

Reactivity: normal

Return to Main Matrix

59%

Shaping up as

57%

Test Item in Hold Tray above

Both Trays

Search Above Tray

Clear Search

Close

Its permanent place, however, is located under the Additional drop-down menu.

Universal Biofeedback System

Home Programs Additional Subspace Programs Reporting Auto Focus Zap Retreat Order of Remedies BPR Enhancements Library System Power Settings Audio Player Help

Coronavirus Immunization

RETEST

Data for John Doe

Soc Index: 109

Owner: DEMO

SYS POW: Normal

Percent Done

Regulatory dysfunction in signal energy flow

Explore risk of oxygenation transfer dysfunction

Explore risk of energy management dysfunction

Reactivity: normal

Return to Main Matrix

59%

Shaping up as

57%

Test Item in Hold Tray above

Both Trays

Search Above Tray

Clear Search

Close

Coronavirus Immunization

Super boosting and stimulating the immune system with the special New Immune Booster from Life Vision™.

Our immune system is essential for our survival. Without a properly working immune system, our bodies would be open to attack from bacteria, viruses, parasites, and more. With Coronavirus, SARS and so many other dangerous viruses causing a pandemic, it is up to our immune system that keeps us healthy and protected as we encounter a sea of pathogens.

Proper diet and supplementation is essential, however may not be enough. The NEW immune booster works especially to harmonize the body while boosting the immune system in alert stage in order to decrease your susceptibility to the virus or other pathogens.

Unconscious Search TREAT Clear Close

0:00 / 0:00 Rectification: -

RETEST

Data for John Doe Soc Index: 109

Owner: DE/IO System: Normal

Regulatory dysfunction in signal energy flow Percent Done: 0%

Explore risk of oxygenation transfer dysfunction

Explore risk of energy management dysfunction Time remaining: 0:00

Reactivity: normal Reactivity / Rectify

Return to Main Matrix 50%

Allersodes Amino Acids Apurvedic Bach Flowers Bacteria Bio War Blood

Search for items

Search Above Tray Clear Search Close

Treat Indic. Reaction Test item in Hold Tray above Both Trays

Hold GetFeedback Load Report Items Test item in Hold Tray above

Search for items

Search Above Tray Clear Search Close

For more on how to use the Coronavirus Immunization panel, please refer to our video on our support site <https://support.quantum-medical.eu> on the Coronavirus Immunization or Chapter 17 in our user manual. The next feature that changed with this update is the Feedback Homeopathic Activation Panel. This can be located under the Program dropdown menu.

Feedback Homeopathic Activation Panel

Super boosting and stimulating the immune system with the special New Immune Booster from Life Vision™.

Our immune system is essential for our survival. Without a properly working immune system, our bodies would be open to attack from bacteria, viruses, parasites, and more. With Coronavirus, SARS and so many other dangerous viruses causing a pandemic, it is up to our immune system that keeps us healthy and protected as we encounter a sea of pathogens.

Proper diet and supplementation is essential, however may not be enough. The NEW immune booster works especially to harmonize the body while boosting the immune system in alert stage in order to decrease your susceptibility to the virus or other pathogens.

Unconscious Search TREAT Clear Close

0:00 / 0:00 Rectification: -

RETEST

Data for John Doe Soc Index: 109

Owner: DE/IO System: Normal

Regulatory dysfunction in signal energy flow Percent Done: 0%

Explore risk of oxygenation transfer dysfunction

Explore risk of energy management dysfunction Time remaining: 0:00

Reactivity: normal Reactivity / Rectify

Return to Main Matrix 50%

Allersodes Amino Acids Apurvedic Bach Flowers Bacteria Bio War Blood

Search for items

Search Above Tray Clear Search Close

Treat Indic. Reaction Test item in Hold Tray above Both Trays

Hold GetFeedback Load Report Items Test item in Hold Tray above

Search for items

Search Above Tray Clear Search Close

All therapies and usage remain the same as previously, only the appearance of the panel has changed, being more user-friendly.

[illegible]

PULSAR™ Universal Biofeedback System

NEW Programs Additional Subspace Programs Reporting Auto Focus Zap Retest Order of Remedies EPR Enhancements Library Syst

DISCLAIMER:

Biofeedback Purposes Only.
Not for Primary Diagnosis.
Please use within your scope of practice and please get good professional education to learn how to use and interpret results then validate + verify results to accepted standards.

Biofeedback Systemic Relaxation

Universal BiofeedBack Therapy

Hyper-Reactant Reactivity Training

Muscular Re-Education for Athletes

Muscular Re-Education for Injured Muscles

Relaxation Trainig


Relaxation Trainig for Spiritual

EMG EEG ECG GSR Biofeedback Chart


Close

10949 Resveratrol Plus. Life Care® | testosterone&inc. benilan. glutathione&elev. anti acina. antibacterial. antioxidant. blood circ. bloo

Click Biofeedback Systemic Relaxation to enter the program.

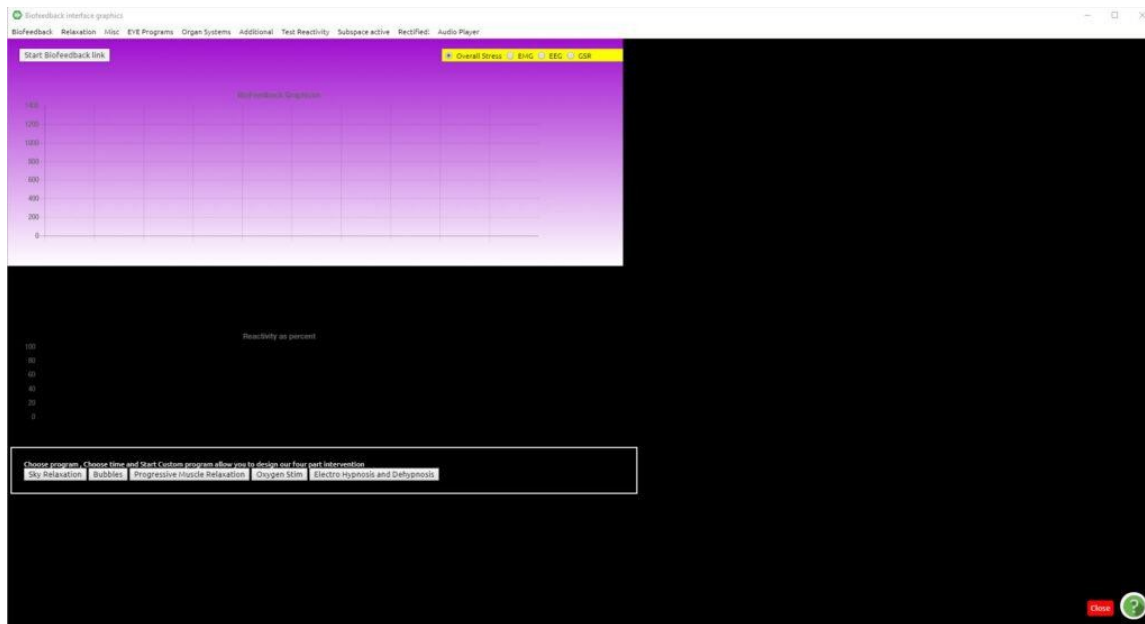

The feature which has changed from previously is the selection of unconscious choice under EACH drop-down menu. The feature now allows you to click the unconscious choice up to three times so that the client's unconscious chooses which therapy it needs.

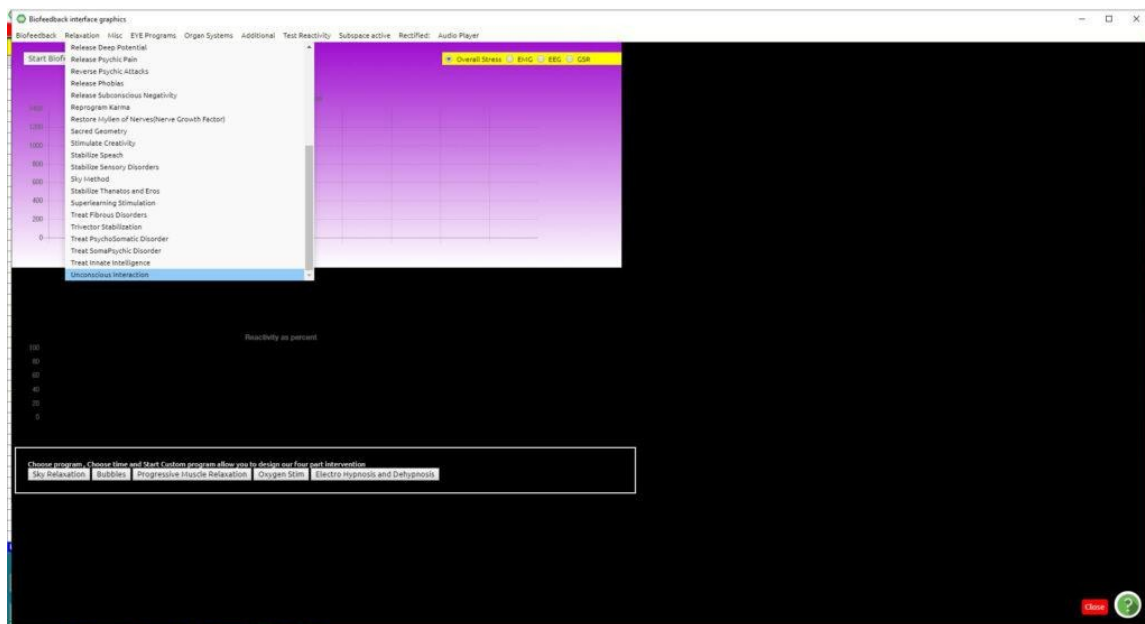


You shall notice a check beside each chosen therapy, while after three selections the “unconscious” is greyed out. The purpose is so that you can choose, as a therapist, what the client needs most. You can repeat this for each drop-down menu, choosing and getting an idea and information about what the client needs most.

For more on how to use the Biofeedback Systemic Relaxation Panel, please refer to our video on our support site <https://support.quantum-medical.eu/> or Chapter 14 in our user manual.

Changes made with Update no. 46


With this update, several changes have been made.

First off, you will notice a change in the VARHOPE panel. The same readings are still there, the only difference is the appearance, neater and more concise. If you wish to see the meanings of each reading, just hover your cursor over the number and the meaning shall appear.

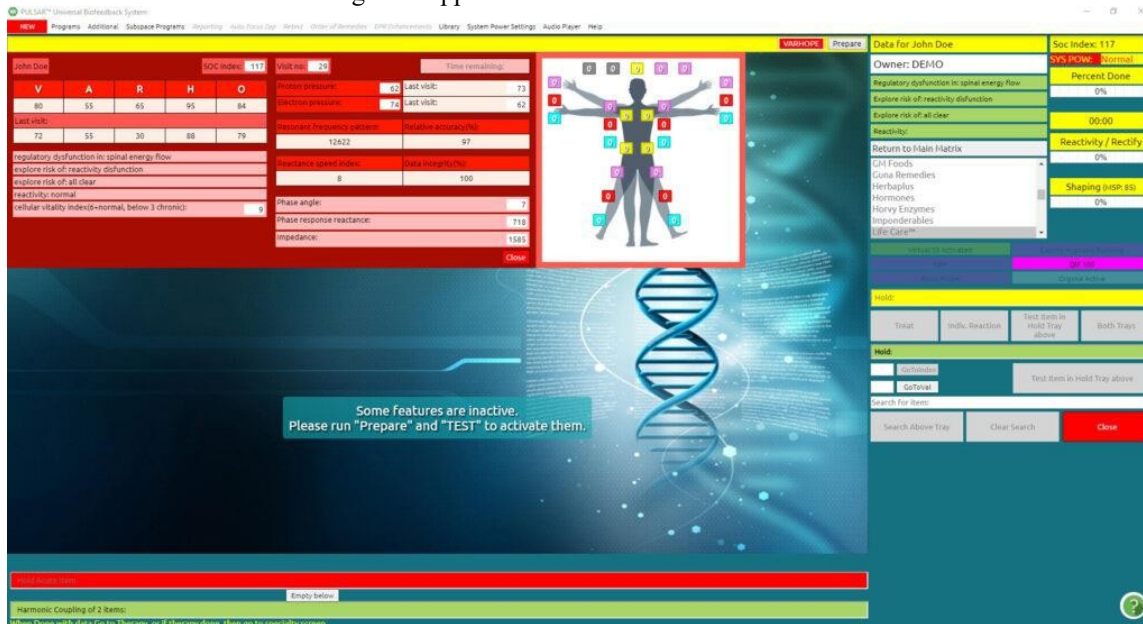

As below, I hovered my cursor over O (oxygenation) and you can see the meanings of the readings.



For more in-depth explanations of the readings please refer to our video VARHOPE Score Interpretation on our support site <https://support.quantum-medical.eu/>

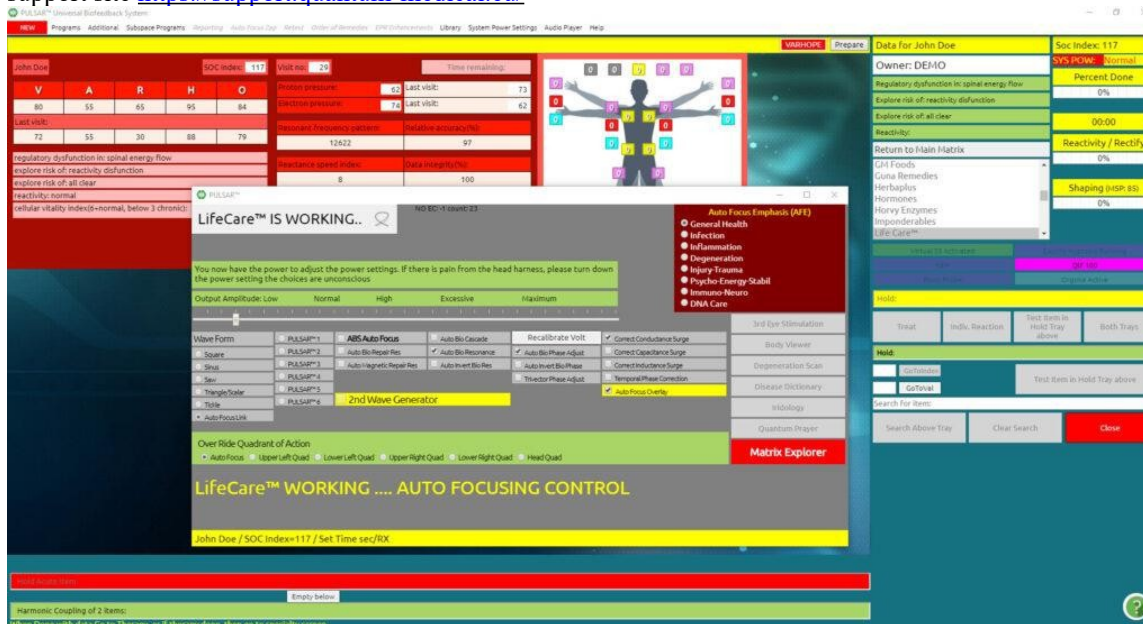

After testing, when the matrix appears, you will notice that the VARHOPE panel remains in view for your review and analysis. If you wish to close the panel, simply press the button close on the right bottom of the red panel.

If you wish to refer to the red VARHOPE panel at some other point during your therapy, simply press the VARHOPE button.

The next new and exciting feature of our software is the Coronavirus Immunization panel. As you already know from our previous updates, all new programs added to our software are initially and temporarily found under the NEW red drop-down menu, for quick reference.

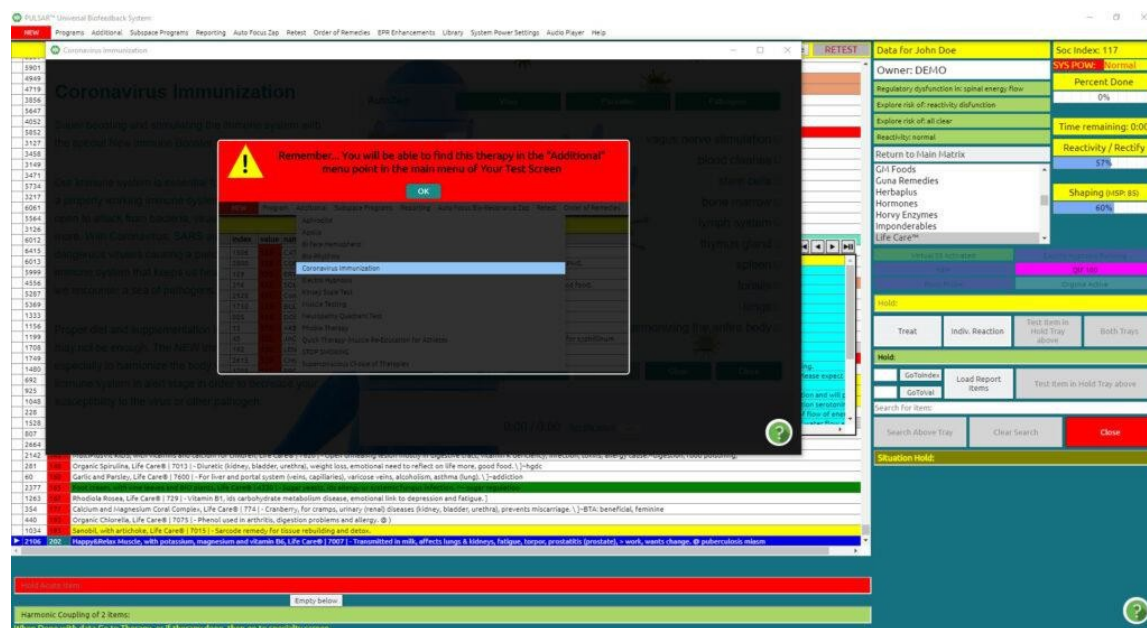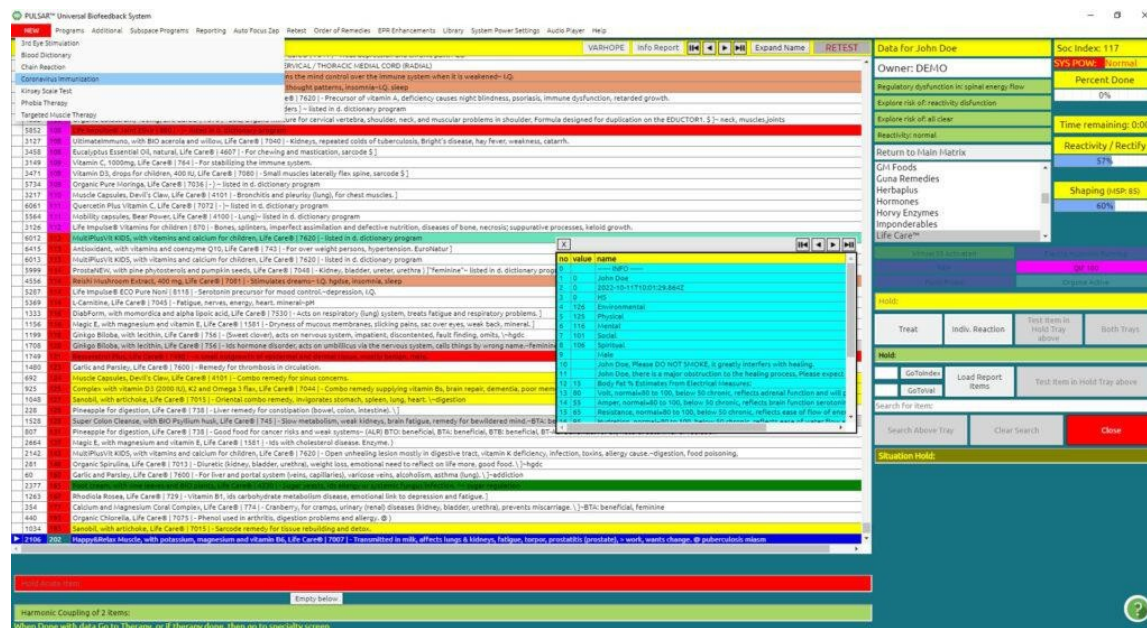

Its permanent place, however, is located under the Additional drop-down menu.

NEW Programs Additional: Subsystem Programs Reporting Auto Focus 2d Retest Order of Remedies EPB Enhancement: Library System Power Settings Audio Player Help

4407	Hagic Root, with BIO plasma, Life Care® 4450 - NERF E Stimul, CBNVAC, ThioCAMP, VEDNA, COMB (BIOPLAS)	Info Report	RETEST	Info Report
4408	Organic Juice, Aloe Vera with pulp, Life Care® 8004 - Strengthens the mind control over the immune system when it is weakened-Up			
4419	Reishi Mushroom Extract, 400 mg, Life Care® 7081 - Repetitive thought patterns, insomnia, sleep			
4420	MultiPhospho KIDS, with vitamins and calcium for children, Life Care® 7020 - Preschool of vitamin A deficiency causes night blindness, porosis, immune dysfunction, retarded growth.			
5647	Beauty Forte, with hyaluronic acid, Life Care® 7030 - Skin disorders - it is listed in a disconnection error.			
4052	Organic Catabolism, 400mg, Life Care® 17 - Coronavirus Immunization			
5013				
3127	Uromaximum, with BIO ascorbic acid and vitamin C, Life Care® 7040 -			
3458	Organic Essential Oil, natural, Life Care® 7040 -			
3459	Vitamin C, 1000mg, Life Care® 7041 -			
3471	Vitamin D3, drops for children, 400 IU, Life Care® 7042 -			
3758	Organic Pure Hydration, Life Care® 7034 -			
3759	Hydroxylation, 500 mg, Life Care® 7035 -			
6061	Quercetin Plus Vitamin C, Life Care® 7036 -			
6062	Ability, 1000mg, Life Care® 7037 -			
2428	Life Impact, 1000mg, Life Care® 7038 -			
6072	MultiPhospho KIDS, with vitamins and calcium, Life Care® 7020 -			
6415	Antioxidant, with vitamins and calcium, Life Care® 7020 -			
6073	Prostate, with vitamins and calcium, Life Care® 7020 -			
3999	Prostate, with vitamins and calcium, Life Care® 7020 -			
4124	Prostate, with vitamins and calcium, Life Care® 7020 -			
5207	Life Impact, ECO Pure, Life Care® 8116 - See			
3369	Life Impact, Life Care® 7041 -			
1333	Prostate, with vitamins and calcium, Life Care® 7020 -			
1156	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1159	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1709	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1709	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1480	Life Impact, with vitamins and calcium, Life Care® 7020 -			
460	Life Impact, with vitamins and calcium, Life Care® 7020 -			
923	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1048	Life Impact, with vitamins and calcium, Life Care® 7020 -			
2107	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1528	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1007	Life Impact, with vitamins and calcium, Life Care® 7020 -			
2464	Life Impact, with vitamins and calcium, Life Care® 7020 -			
2143	Life Impact, with vitamins and calcium, Life Care® 7020 -			
287	Life Impact, with vitamins and calcium, Life Care® 7020 -			
40	Life Impact, with vitamins and calcium, Life Care® 7020 -			
2377	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1381	Life Impact, with vitamins and calcium, Life Care® 7020 -			
354	Life Impact, with vitamins and calcium, Life Care® 7020 -			
440	Life Impact, with vitamins and calcium, Life Care® 7020 -			
1013	Life Impact, with vitamins and calcium, Life Care® 7020 -			

AutoZap: Virus Parasites Pathogens

ALPHA BETA GAMMA DELTA OMEGON

fractal auto scalar harmonizing the entire body

Unconscious Search TREAT CLEAR CLOSE

0:00 / 0:00 Rectification

Owner: DEMO

Regulatory dysfunction in spinal nerve flow

Explore risk of reactivity dysfunction

Time remaining: 0:00

Reactivity / Rectify

57%

Shaping (BIO) 80%

60%

blood cleanse

stem cells

bone marrow

lymph system

thymus gland

spleen

tonsils

lungs

Test items in hold tray above

Test items in hold tray above

Clear Search

Close

Super boosting and stimulating the immune system with the special New Immune Booster from LifeCare™.

Our immune system is essential for our survival. Without a properly working immune system, our bodies would be open to attack from bacteria, viruses, parasites, and more. With Coronavirus, SARS and so many other dangerous viruses causing a pandemic, it is up to our immune system that keeps us healthy and protected as we encounter a sea of pathogens.

Proper diet and supplementation is essential, however may not be enough. The NEW Immune booster works especially to harmonize the body while boosting the immune system in alert stage in order to decrease your susceptibility to the virus or other pathogen.

See index: 117

SYS POW: Normal

Percent Done

0%

PULSAR® Universal Biofeedback System

NEW Programs Additional Subspace Programs Reporting Auto Focus Zap Retest Order of Remedies ERP Enhancements Library System Power Settings Audio Player Help

VARIOPE Info Report Expand name RETEST

Data for John Doe Soc Index: 117
Owner: DEMO
Regulatory dysfunction in spinal energy flow
Explore risk of reactivity dysfunction
Explore risk of all clear
Reactivity: normal
Return to Main Matrix
GFI Foods
Guna Remedies
Herbapilus
Hormones
Horny Enzymes
Imponderables
Life Care™

Sys POW: Normal
Percent Done
0%
Time remaining: 0:00
Reactivity / Rectify
57%
Shaping (min 80)
60%

Hold

Treat Indiv. Reaction Test item in Hold Tray above Both Trays

Search For Item

Search Above Tray Clear Search Close

Situation Hold:

Empty below

Harmonic Coupling of 2 Items

When Done with data Go to Therapy, or if therapy done, then go to specialty screen

DISCLAIMER:
Biofeedback Purposes Only.
Not for Primary Diagnosis.
Please use within your scope of practice and please get good professional education to learn how to use and interpret results then validate + verify results to accepted standards.

Biofeedback Systemic Relaxation
Universal Biofeedback Therapy
Hyper-Reactant Reactivity Training
Muscular Re-Education for Athletes
Muscular Re-Education for Injured Muscles
Relaxation Training
Relaxation Training for Spiritual
EMG EEG ECG GSR Biofeedback Chart

Close

Empty below

Harmonic Coupling of 2 Items

When Done with data Go to Therapy, or if therapy done, then go to specialty screen

PULSAR® Universal Biofeedback System

NEW Programs Additional Subspace Programs Reporting Auto Focus Zap Retest Order of Remedies ERP Enhancements Library System Power Settings Audio Player Help

VARIOPE Info Report Expand name RETEST

Data for John Doe Soc Index: 117
Owner: DEMO
Regulatory dysfunction in spinal energy flow
Explore risk of reactivity dysfunction
Explore risk of all clear
Reactivity: normal
Return to Main Matrix
GFI Foods
Guna Remedies
Herbapilus
Hormones
Horny Enzymes
Imponderables
Life Care™

Sys POW: Normal
Percent Done
0%
Time remaining: 0:00
Reactivity / Rectify
57%
Shaping (min 80)
60%

Hold

Treat Indiv. Reaction Test item in Hold Tray above Both Trays

Search For Item

Search Above Tray Clear Search Close

Situation Hold:

Empty below

Harmonic Coupling of 2 Items

When Done with data Go to Therapy, or if therapy done, then go to specialty screen

DISCLAIMER:
Biofeedback Purposes Only.
Not for Primary Diagnosis.
Please use within your scope of practice and please get good professional education to learn how to use and interpret results then validate + verify results to accepted standards.

Biofeedback Systemic Relaxation
Universal Biofeedback Therapy
Hyper-Reactant Reactivity Training
Muscular Re-Education for Athletes
Muscular Re-Education for Injured Muscles
Relaxation Training
Relaxation Training for Spiritual
EMG EEG ECG GSR Biofeedback Chart


Close

Empty below


Harmonic Coupling of 2 Items

When Done with data Go to Therapy, or if therapy done, then go to specialty screen

Click Biofeedback Systemic Relaxation to enter the program.

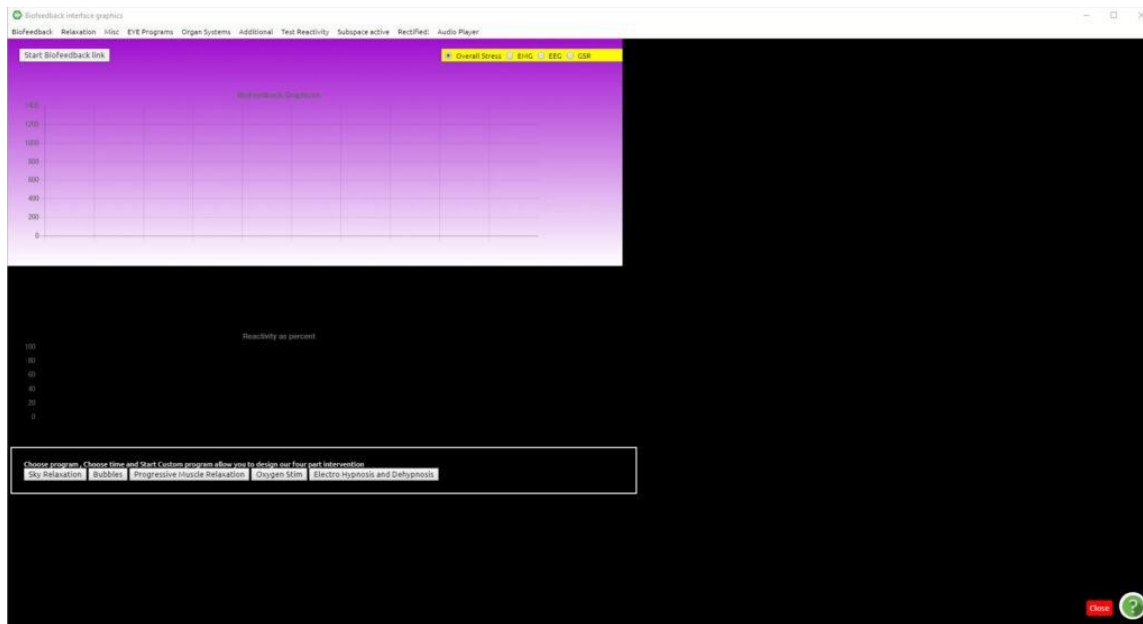

The feature which has changed from previously is the selection of unconscious choice under EACH drop-down menu. The feature now allows you to click the unconscious choice up to three times so that the client's unconscious chooses which therapy it needs.

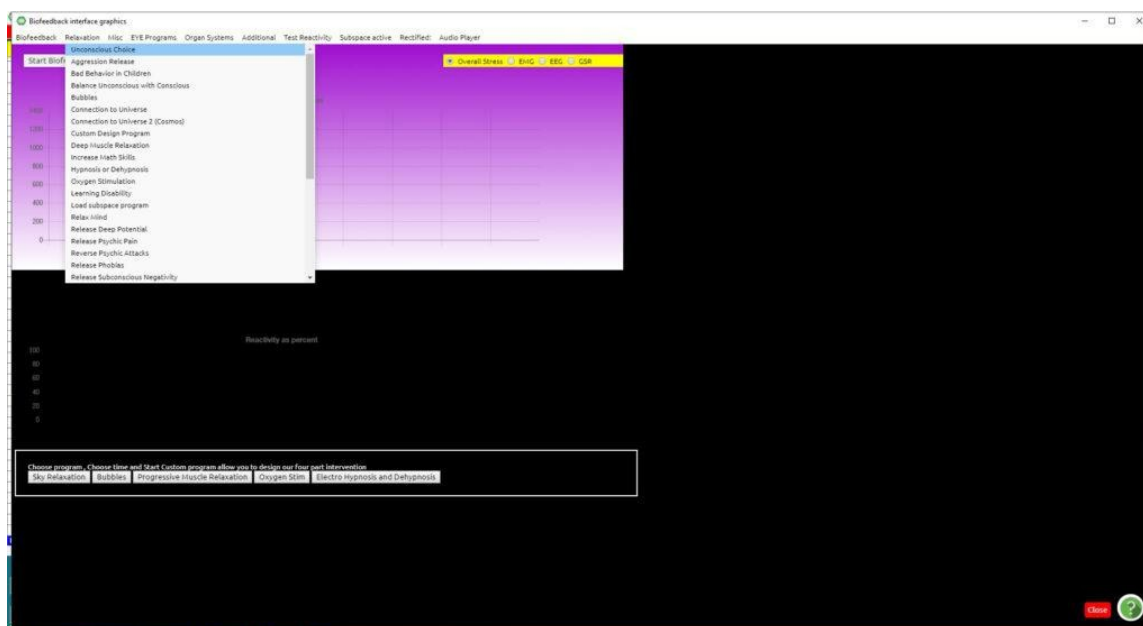

You shall notice a check beside each chosen therapy, while after three selections the “unconscious” is greyed out. The purpose is so that you can choose, as a therapist, what the client needs most. You can repeat this for each drop-down menu, choosing and getting an idea and information about what the client needs most.

For more on how to use the Biofeedback Systemic Relaxation Panel, please refer to our video on our support site <https://support.quantum-medical.eu/>

